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
Onglet 2 – Identité

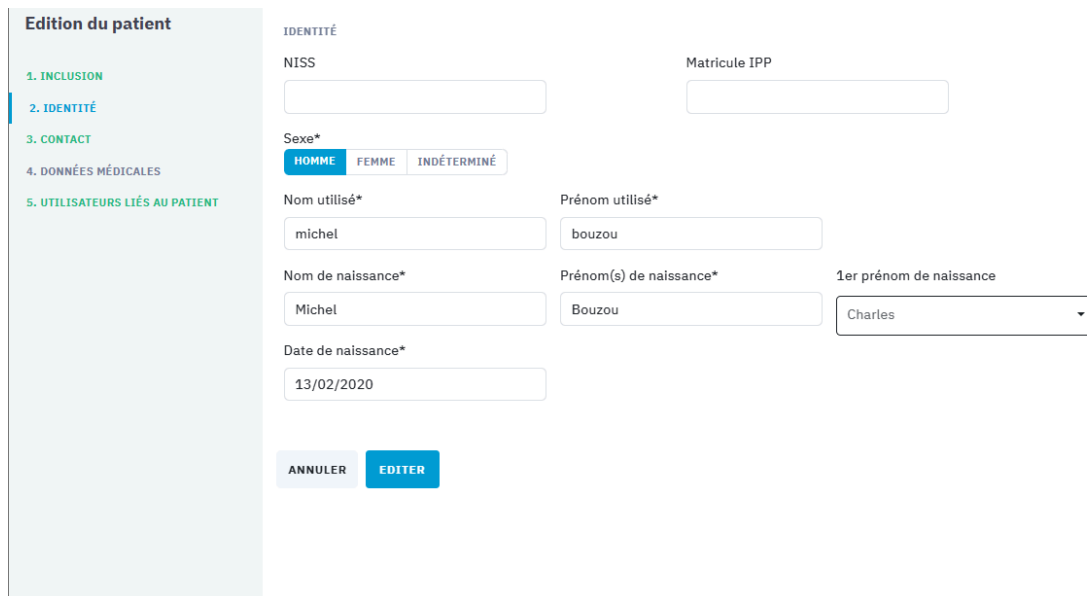A screenshot of a web application interface for editing a patient's information. On the left is a sidebar titled "Edition du patient" with a list of tabs: "1. INCLUSION", "2. IDENTITÉ" (which is highlighted with a blue bar), "3. CONTACT", "4. DONNÉES MÉDICALES", and "5. UTILISATEURS LIÉS AU PATIENT". The main area is titled "IDENTITÉ" and contains several form fields. At the top, there are two input fields: "NISS" and "Matricule IPP". Below these is a "Sexe*" section with three buttons: "HOMME" (highlighted in blue), "FEMME", and "INDÉTERMINÉ". Then, there are two input fields: "Nom utilisé*" (containing "michel") and "Prénom utilisé*" (containing "bouzou"). Below these are three input fields: "Nom de naissance*" (containing "Michel"), "Prénom(s) de naissance*" (containing "Bouzou"), and "1er prénom de naissance" (a dropdown menu showing "Charles"). At the bottom of the form fields is a "Date de naissance*" input field containing "13/02/2020". At the very bottom of the form area are two buttons: "ANNULER" and "EDITER" (highlighted in blue).

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The screenshot shows a web application dashboard. At the top, there's a header with a menu icon, the title 'Tableau De Bord', a search bar, and a button 'AJOUTER UN PATIENT'. Below the header, a sidebar on the left contains icons for home, notifications, and user profile. The main content area starts with a greeting 'Bonjour, Dr Test Soraya Rafa Newcard'. The primary section is 'Alertes médicales' with a link 'VOIR TOUT (9)'. It lists three alerts, each in a yellow box with a '1' in a circle, indicating a high priority. The alerts are for MARCHAND Ansberte, PAUL Émérance, and DUFOUR Alverède, all dated 05/03/2026, with the note 'Tension haute'. Below the alerts, there's a section 'Patients à renouveler' containing a table with two rows: Monique carpentier (+33311564439) and Esther nicard (+33302484779), each with a 'RENOUVELER' button. To the right of this table is a 'Statistiques patients' section showing two icons with counts 131 and 53.

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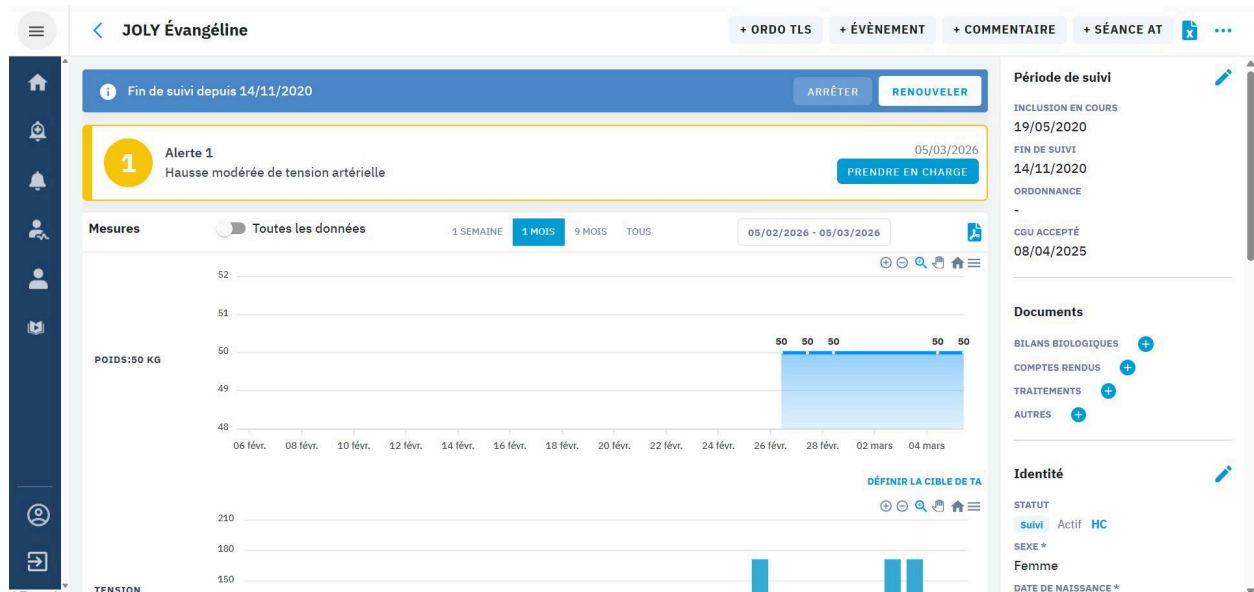
The screenshot shows the 'Prise en charge de l'alerte' (Alert management) form. It is titled 'Alerte 1 Déclenchée' (Alert 1 Triggered) and shows the patient 'Ansberte Marchand'. The form includes a section for 'PRENDRE EN CHARGE' (Take charge) with a dropdown for 'Action effectuée*' (Action performed*) and buttons for 'ANALYSE SUR LA PLATEFORME' (Analysis on the platform), 'APPEL DE PATIENT' (Patient call), and 'CONSULTATION' (Consultation). There's also a section for 'Résolution de l'alerte' (Alert resolution) with a dropdown for 'Retour Spontané à la normale' (Spontaneous return to normal). The form also includes a 'Commentaire' (Comment) field, a 'Pertinence de l'alerte' (Alert relevance) section with a star rating, and a 'HISTORIQUE' (History) section with a dropdown for 'Alerte en cours' (Alert in progress) and 'Alerte terminée' (Alert ended).

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USER GUIDE FOR THE NEWCARD PLATFORM

for Belgian users

Summary

1. Welcome	58
2. The basics of remote monitoring	59
3. First time logging in to the platform	64
4. Overview of the platform	67
a. Dashboard	67
b. Feature sidebar	68
c. Changes to personal data	73
d. Access the help section	76
5. Inclusion of a patient	78
6. Patient list	83
7. The patient record	84
a. Administrative information	84
b. Alert in progress	86
c. Access to remote monitoring data	86
d. Visualisation of medical and biological data	89
e. Stocks on the charts	91
f. Patient history	91
g. Generating a PDF export of remote monitoring data	94
h. Actions on the patient's profile	95
i. Sending text messages	96
8. Alert management	97
9. Notification management	99
10. Profile and user management	100
11. Contact us	115
	170

1. Welcome

Your practice/department has expressed an interest in using NEWCARD for remote patient monitoring.

This user guide has been designed to guide you, step by step, through the process of getting started with the platform.

Remote medical monitoring enables the early detection of signs that a chronic condition is worsening, allowing complications to be anticipated and care to be adapted. It also allows for closer monitoring of patients, thanks to the regular transmission of data from their homes. This system is not a substitute for emergency care.

At Newcard, we place particular importance on personalised support, offering tailored care for patients, carers and healthcare professionals.

Throughout your experience, you can count on the expertise and availability of our paramedical, administrative and sales teams.

About the guide

The guide forms an integral part of the “1 Minute pour mon Coeur” (1 Minute for my Heart) software (**basic IUD-ID 3770034348MPC6G**) and ‘1 Minute pour mes Reins” (1 Minute for my Kidneys) software (**basic IUD-ID 3770034348MPR7E**). It provides a detailed overview of the software, its configuration and how to use it. Its purpose is to help users understand how their remote monitoring works and to support them throughout the process.

If required, a printed version can be provided on request by emailing <https://newcard.io/contacts/> ou cardio@newcard.io.

There is no additional charge for this service, and the user manual will be delivered to you within seven calendar days of receipt of your request.

2. The basics of remote monitoring

Intended use of the solutions

- ‘1 Minute pour mon Cœur (1 Minute for my Heart) is a software-as-a-medical-device solution designed for the remote monitoring of patients with chronic heart failure through the collection, analysis, storage and visualisation of health data.’

Log in to the platform

You have access to the Newcard platform to register and monitor your patients remotely.

<https://app.newcard.io>

Your username is the email address you provided to us.

Note: eligibility criteria

Specific aspects in cardiology (heart failure):

Patients with chronic heart failure who meet the following eligibility criteria, as defined within the remote monitoring programme:

- **Hospitalisation** for heart failure;
- **OR Hospitalisation** complicated by heart failure
- **AND** classified as **NYHA class II, III or IV** following stabilisation of acute decompensation.

Contraindications: criteria for ineligibility

Patients in any of the following situations are considered ineligible for remote monitoring:

- physical or mental inability on the part of the patient or carer to use the digital medical telecare device and/or its data collection accessories, as determined by the doctor wishing to include the patient in the medical telecare programme ;
- the patient's refusal to provide the data required to monitor the actual use of the digital health platform and to obtain individual or national real-world usage results ;
- the patient's refusal to receive therapeutic support.

Description of the labelling / explanation of the symbols

	Medical device
	Unique device identification code
	Name and address of the authorised manufacturer
	Batch number
	CE marking
	See the user guide
	Trade reference
	Packaging and items not to be disposed of
	Fragile
	Keep away from heat

Precautions

6. **This device is for personal use only.**
7. **The kit is on loan to you from Newcard for the duration of your remote monitoring. At the end of your monitoring period, you are required to return the device as soon as possible and in the same condition as when you received it. We would also be grateful if you could keep the original boxes and packaging, as these will be useful when returning the equipment.**
8. **Privacy and security:** Login credentials and passwords are strictly personal and must never be shared with a third party. Users are advised to use a strong password, enable two-factor authentication (2FA), and ensure their devices are protected in order to guarantee the confidentiality and security of their data. Any suspicious or unauthorised use of the account must be reported promptly to the relevant support team.
9. **The Newcard service is not an emergency service;** it is not intended to replace the intervention of healthcare professionals. Recorded data or alerts are not read and interpreted instantly, but only during the days and times specified by the remote monitoring operator. If you have a medical query or your condition suddenly worsens, please contact a healthcare professional or dial 15.
10. **Preventive and routine maintenance operations:** Newcard may modify or adapt the NEWCARD app or the Terms and Conditions of Use. In this case, you may need to update the app. In accordance with the Terms and Conditions of Use, this update must be carried out as soon as possible.

Intended users

Newcard's '1 minute pour mon coeur' (1 minute for my Heart) platform is intended for use by:

- Healthcare professionals involved in the treatment and management of chronic heart failure.

The patient app is intended for use by:

- Eligible patients with chronic heart failure.

Target patient groups

- Patients, both men and women, with chronic heart failure.

Clinical benefits

The purpose of telemonitoring is to enable a medical professional to remotely interpret the data necessary for the medical follow-up of a patient and, if necessary, to make decisions relating to their care. Using the Newcard software offers the following clinical benefits:

- **Improved patient survival (reduced mortality)**
(among patients with heart failure followed up via remote monitoring)
- **Improved patient independence**
(as part of remote monitoring for heart failure and kidney failure)

Technical specifications and requirements

An internet connection is required to use the app.

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Browser	Compatible version
Google Chrome	≥ 126.0
Hardware	
Desktop computer	Screen size: minimum 17 inches with a resolution of 1920 x 1080 (Full HD)
Laptop	Screen size: minimum 15 inches with a resolution of 1920 x 1200 (Full HD+)
Tablets	Screen size: 10 inches with a resolution of 800 x 1280
Smartphones	• Screen size 6.5 inches with a resolution of 412 x 915: for Android • Screen size 6.5 inches with a resolution of 414 x 896: for iPhone.

Mobile operating system	Compatible version
Android	≥ 10
iOs	≥ 10
Hardware	
Tablets	Screen size: 10 inches with a resolution of 800 x 1280
Smartphones	• Screen size 6.5 inches with a resolution of 412 x 915: for Android • Screen size 6.5 inches with a resolution of 414 x 896: for iPhone.

In the event of a medical device vigilance incident

Any incident or potential serious incident that has resulted in, or could result in, the death or serious deterioration in the health of a patient, user or third party involving Newcard's services must be reported immediately to Newcard at the following email address: contact@newcard.io and to the competent authority **(in France, the French National Agency for Medicines and Health Products Safety – ANSM)**.

Residual risks

The software has no direct physical contact with the patient and poses no mechanical, electrical or biological risk. The residual risks identified are of a technical or software-related nature; although rare, they may, in certain situations, have an indirect impact on patient care. The main residual risks associated with the use of this software are described below.

5. Appropriate technical and organisational measures have been put in place to ensure the platform remains available. Despite these controls, there remains a residual risk that the platform may be temporarily unavailable until it is restored (in the event of an internet connection failure, a technical error, etc.).

6. The platform incorporates mechanisms designed to ensure the reliability of the alerts it generates. Nevertheless, a residual risk of alert errors—including, in particular, false positives, false negatives, duplicates or alerts failing to appear—cannot be entirely ruled out.
7. Despite the controls in place, errors may occur in the transmission of patient measurements (failure to transmit data, or data that is inaccurate or unreliable).
8. Newcard has implemented appropriate cybersecurity measures to ensure the security of all data transmitted. Nevertheless, a residual risk of loss, alteration or unauthorised access to this data cannot be entirely ruled out. In all cases, Newcard is committed to safeguarding users' privacy and complying with applicable data protection regulations, in particular the GDPR.

Patient inclusion

Patients are enrolled directly via the platform.

Once a patient has been enrolled, Newcard takes care of:

- contact the patient;
- organise training on how to use the equipment;
- monitor compliance (ensuring measurements are taken regularly);
- deal with any technical issues that may arise.

Handling alerts

Remote monitoring data and triggered alerts are directly accessible on the platform.

When an alert is triggered :

- you receive an email;
- a link takes you to the patient's file;
- you decide on the course of treatment to follow.

3. First time logging in to the platform

When you log in to the Newcard remote monitoring platform for the first time, simply enter your login details.

You will then be asked to accept the Terms and Conditions of Use (T&Cs).

By accepting these T&Cs, you acknowledge that you have read and understood how the service works and you agree to:

- keep your login details confidential;
- use the platform for its intended purpose;
- not interfere with the proper functioning of the service.

The Terms and Conditions are available at any time in the **‘Documents’** tab of your profile.

Two-factor authentication

Two-factor authentication (2FA) enhances security when accessing Newcard by adding an extra verification step on top of the password.


In light of the growing risks of phishing, data breaches and compromised login credentials, this measure helps to protect the health data processed on the platform.

It is required on first login and then every fortnight by default for all users.

Login and two-factor authentication steps


1. Log in

Log in using your email address and password.

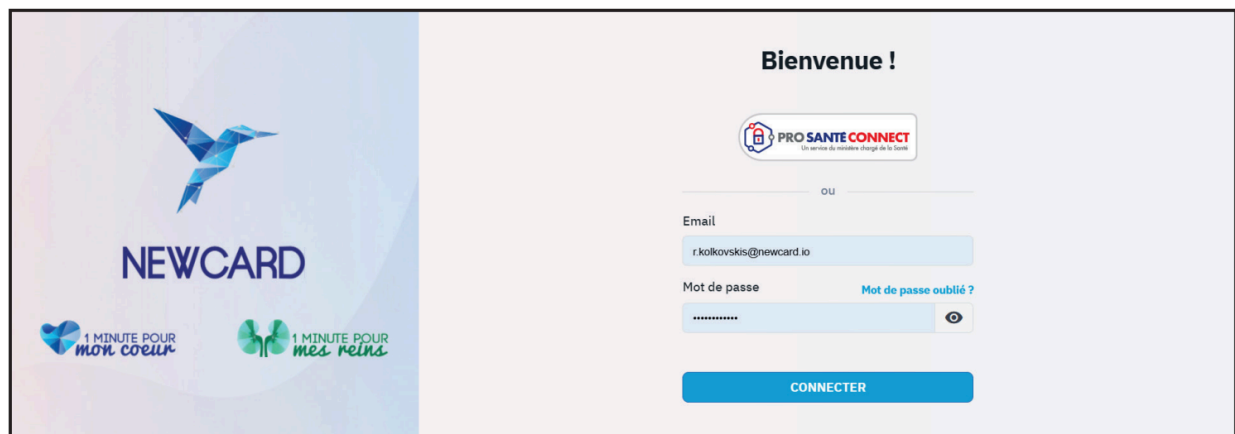The image shows the Newcard login interface. On the left, there is a blue and white graphic with a stylized bird logo and the text 'NEWCARD'. Below the logo, there are two small icons: a heart with the text '1 MINUTE POUR MON COEUR' and a kidney with the text '1 MINUTE POUR MES REINS'. On the right, there is a light blue background with the text 'Bienvenue !' at the top. Below this, there is a logo for 'PRO SANTE CONNECT' with the tagline 'Un service du ministère chargé de la Santé'. Underneath the logo, there is a horizontal line with the text 'OU' in the center. Below this, there are two input fields: 'Email' with the value 'r.kolkovskis@newcard.io' and 'Mot de passe' with a masked password '*****'. To the right of the password field, there is a link that says 'Mot de passe oublié ?'. At the bottom, there is a blue button with the text 'CONNECTER'.

2. Request for verification


A screen appears asking you to enter a verification code.

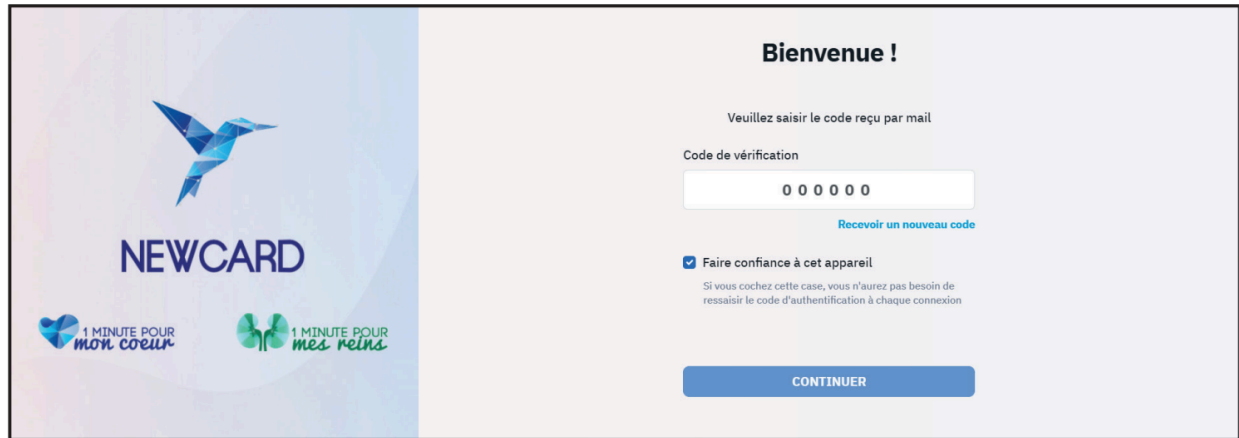The image shows a mobile application screen for Newcard. On the left, there is a blue bird logo and the text 'NEWCARD'. Below it are two logos: '1 MINUTE POUR mon cœur' and '1 MINUTE POUR mes reins'. On the right, the screen says 'Bienvenue !' and 'Veuillez saisir le code reçu par mail'. There is a text input field for 'Code de vérification' with the placeholder '0 0 0 0 0'. Below the field is a link 'Recevoir un nouveau code'. There is a checkbox labeled 'Faire confiance à cet appareil' which is checked. Below it is a note: 'Si vous cochez cette case, vous n'aurez pas besoin de ressaisir le code d'authentification à chaque connexion'. At the bottom is a blue button labeled 'CONTINUER'.


3. Receipt of the code

A 6-digit code will be sent to you by email.

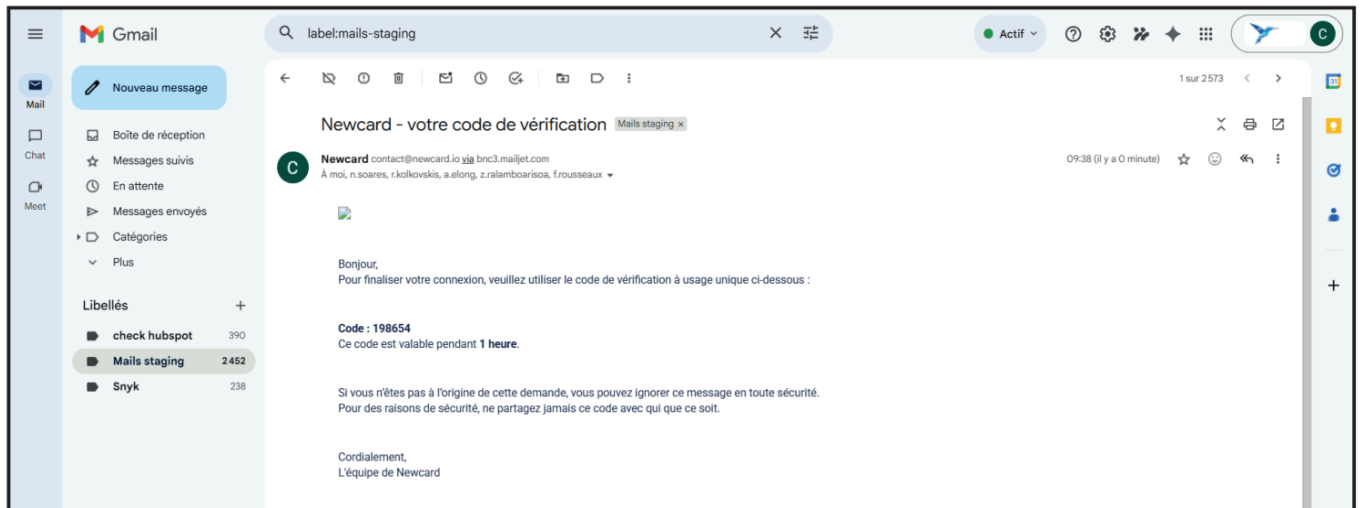4. Enter the code

Enter this code to confirm your access

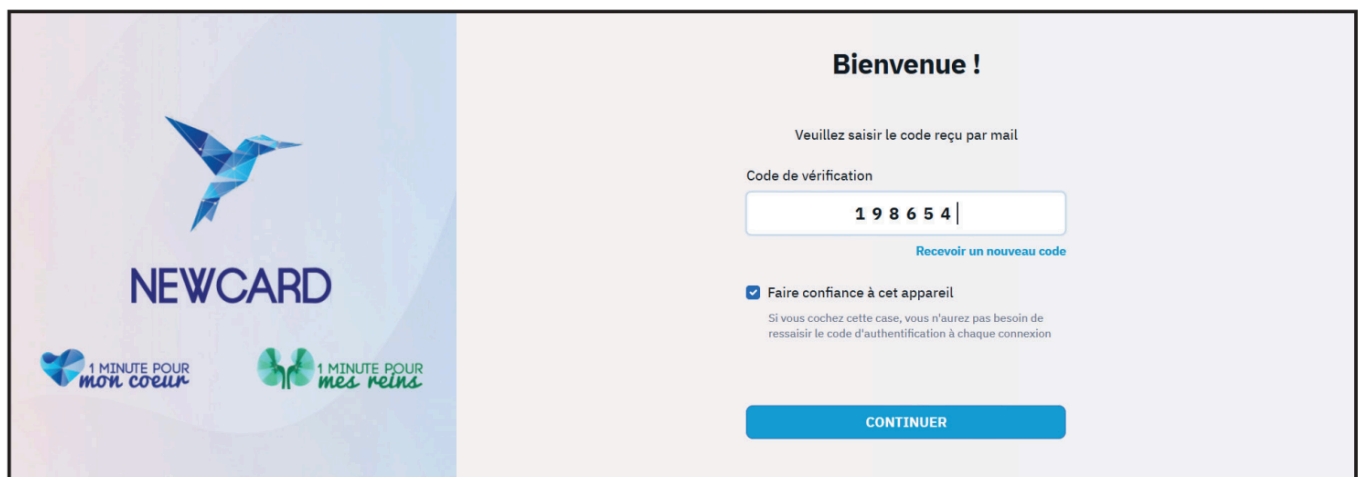The image shows the same mobile application screen as in step 2, but now the verification code '1 9 8 6 5 4' has been entered into the 'Code de vérification' field. The 'CONTINUER' button is still visible at the bottom.

5. Access to the platform

Your login has been verified.

You can now access the Newcard platform securely

Frequency of checks

By default, two-factor authentication is required every 15 days. You can adjust this frequency (from 1 to 30 days) in your profile, under the “Settings” section.

The screenshot displays the 'Profile' page of the Newcard platform, with the 'Configuration' tab selected. The page is organized into several sections: 'NOTIFICATIONS', 'Medical Alerts', 'Patient Hospitalization', 'FILTERS', and 'AUTHENTICATION'. The 'AUTHENTICATION' section is highlighted with a red box, showing a dropdown menu for 'Duration of validity of the two-factor authentication' set to '15 days'. A 'SAVE' button is located at the bottom of the page.

Notification Type	Main Email	Secondary Email	SMS
Medical Alerts	☐	☒	☐
Patient Hospitalization	☒	☐	☐

FILTERS

Only see patients where I am the prescriber ☐

AUTHENTICATION

Duration of validity of the two-factor authentication: 15 days ▾


SAVE

4. Overview of the platform

The platform is a digital tool that allows you to access your patients' data and view alerts. You will also find a list of healthcare professionals involved in the remote monitoring of your patients.

a. Dashboard

On the platform's login page, you will be taken to the remote monitoring platform's dashboard.

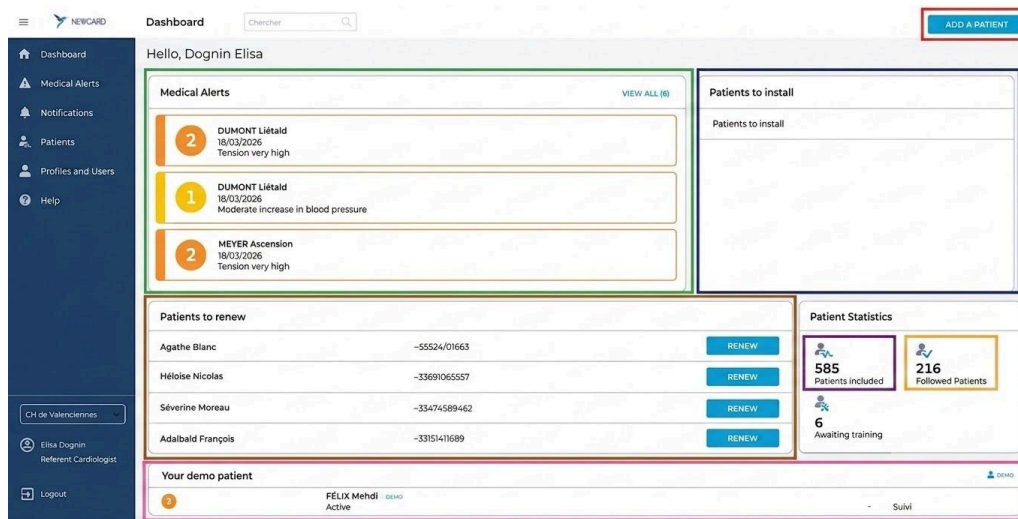This way, you can easily find :

The button to **add a patient**,

- ✓ Current **medical alerts**,
- ✓ The list of **patients to be admitted**
- ✓ The number of **patients enrolled**,
- ✓ The number of **patients monitored**,
- ✓ The list of your **'demo' patients**.

All these links are clickable and will take you directly to the relevant pages. These items can also be found in the categories on the sidebar of your platform.


Please note: some sections are conditional and are displayed or hidden depending on the permissions assigned to each user. As a result, the dashboard may vary from one user to another.

→ For further details on managing accounts and user rights, please refer to section 11 of this guide.

In this example, we outline a configuration with full access (a user with all permissions).

b. Feature sidebar

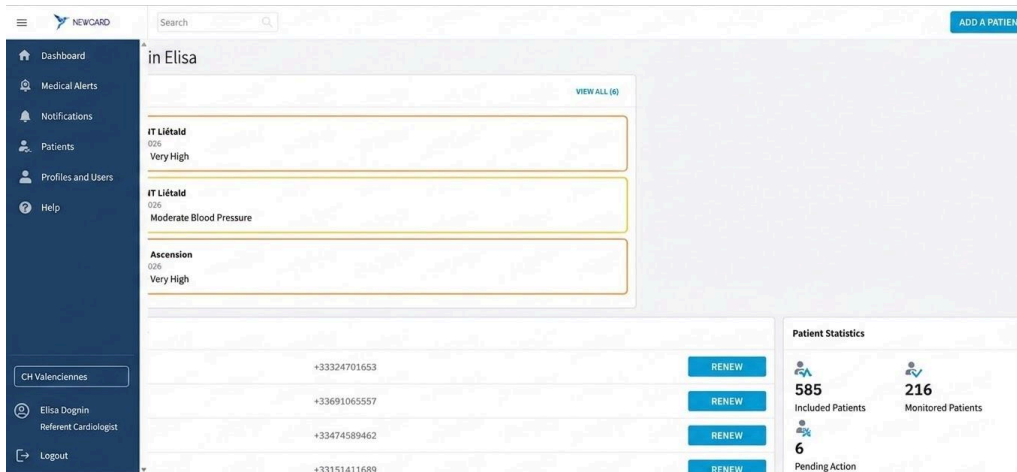The sidebar contains various categories of features.

Please note: Access to the tabs depends on the permissions assigned to each user (with the exception of the **Dashboard** and **Help**, which are accessible to everyone). As a result, the sidebar may vary from one user to another.

For example :

- Only users with the '**Profile and User Management**' permission can access the '**Profiles and Users**' tab. We recommend assigning this permission to the school's designated contact person.t.
- Only users with **billing** permissions can access the '**Billing**' tab.
- etc.

→ For further details on managing accounts and user rights, please refer to section 11 of this guide.

Overview of the different tabs

When you click on one of the categories in the sidebar, you are taken to the corresponding tab, which gives you access to a specific feature of the platform.

The screenshot shows the 'Dashboard' tab selected in the sidebar. The main content area displays 'Hello, Dognin Elisa' and a 'Medical Alerts' section with three alerts. Below this is a 'Patients to renew' table with four rows. To the right is a 'Patient Statistics' box showing 585 Patients included and 216 Monitored Patients. At the bottom, there is a 'Your demo patient' section for 'FÉLIX Mehdi'.

Medical Alerts

Alerts	Alerts	Alerts
2	DUMONT Liétald 18/03/2026 Very high blood pressure	
1	DUMONT Liétald 18/03/2026 Moderate blood pressure increase	
2	MEYER Ascension 18/03/2026 Tension très haute	

Patients to renew

Patients to renew	Patients to renew	Patients to renew
Agathe Blanc	+33324701653	RENEW
Héloïse Nicolas	+33699706557	RENEW
Séverine Moreau	+33474589442	RENEW
Adalbald François	+33151411689	RENEW

Patient Statistics

Patient Statistics	Patient Statistics	Patient Statistics
585 Patients included	216 Monitored Patients	
6 En attente de formation	Awaiting training	

Your demo patient

Your demo patient	Your demo patient	Your demo patient
FÉLIX Mehdi Active	DEMO	Monitored

6. The 'Medical Alerts' tab allows you to view all your current potential alerts and provides an overview of these alerts.

The screenshot shows the 'Medical Alerts' tab. It features a search bar and a table with columns: ALERT, FULL NAME, STATUS, LAST TAKEN RISE EN CHARGE, VIEWED BY SPÉCIALIST, TRIGGERED BY, and ACTIONS. The table lists six alerts, with the first three being 'Triggered' and the last two being 'In Progress'. The 'ACTIONS' column contains buttons for 'TAKE CHARGE' and 'MANAGEMENT IN PROGRESS'.

Medical Alerts

ALERT	FULL NAME	STATUS	LAST TAKEN RISE EN CHARGE	VIEWED BY SPÉCIALIST	TRIGGERED BY	ACTIONS
2	DUMONT Liétald	Triggered	18/03/2026	20/01/2026	Very High Blood Pressure	TAKE CHARGE
1	DUMONT Liétald	Triggered	18/03/2026	20/01/2026	Moderate Blood Pressure Rise	TAKE CHARGE
2	MEYER Ascension	Triggered	18/03/2026	27/01/2025	Moderate Blood Pressure Rise	TAKE CHARGE
1	MEYER Ascension	Triggered	18/03/2026	27/01/2025	Moderate Blood Pressure Rise	TAKE CHARGE
1	PIERRE Joëlle	In Progress	24/02/2026	09/02/2026	Very Gradual Weight Gain	MANAGEMENT IN PROGRESS
1	MICHEL Agathon	In Progress	05/02/2026	09/02/2026	Clinical Degradation	MANAGEMENT IN PROGRESS

Showing element 1 to 6 of 6

→ See section 9 for further information

7. The '**Notifications**' tab allows you to interact with Newcard and receive alerts about any new information concerning your patients (lab results, comments, documents, etc.).

ALERTE	NAME, SURNAME	TYPE	NOTIFICATION DATE	ACTIONS
	BARBIER Guenièvre	Comments	19/03/2026	PROCESS
	BARBIER Guenièvre	Bio Records	19/03/2026	PROCESS
	BARON Antonin	Documents	19/03/2026	PROCESS
	BARON Antonin	Documents	19/03/2026	PROCESS
	BARON Antonin	Documents	19/03/2026	PROCESS
	AUBRY Lothaire	Documents	19/03/2026	PROCESS
	AUBRY Lothaire	Documents	19/03/2026	PROCESS
	AUBRY Lothaire	Documents	19/03/2026	PROCESS

8. The '**Patients**' tab lists all the patients you are managing, whether they are awaiting training, currently being monitored, due for renewal, or deactivated. Using the search bar, you can search for a patient directly by name.

ALERTE	NAME, FIRST NAME	TELEPHONE	STATE	SPECIALIST	INCLUSION DATE	ACTION(S)
Deactivated	ADAM Agrippin	+33355739636	Followed	MERCIER Thomas	11/02/2025	
Deactivated	ADAM Archibald	+33348900875	Followed	MERCIER Thomas	22/09/2025	
Active	ADAM Olive	+33172153060	Followed	MERCIER Thomas	12/08/2025	RENEW
Active	ADAM Alexandre	+33640606022	Followed	PINCHON Emilie	21/08/2025	RENEW
Active	ADAM Danielle	+33642566865	Followed	DOGNIN Elisa	14/09/2025	RENEW
Deactivated	ANDRE Anthelme	+33322279053	Followed	MERCIER Thomas	08/09/2023	
Active	ANDRE Aurélienne	+33167761925	Followed	MERCIER Thomas	13/08/2025	RENEW
Active	ANDRE Judicaël	+33579767604	Followed	MERCIER Thomas	21/09/2025	RENEW
Deactivated	ARNAUD Avoys	+33574120764	Followed	MERCIER Thomas	12/11/2022	
Deactivated	ARNAUD Charlaine	+33197703878	Followed	MERCIER Thomas	17/11/2022	

This button  allows you to edit your patient's key details, namely their **registration details**, **personal details**, **contact information** and **medical records**.

You can also filter your patients by their medical condition using the drop-down menu at the top of your screen, next to the search bar.


Tous	Awaiting training	Monitored	Follow-up	Hospitalized	Absence	Back from hospitalization	Back from an absence	Installation in progress	Deactivated
585	6	2						0	363
ALERT	LAST NAME, FIRST NAME		ÉTATE	SPÉCIALIST	INCLUSION DATE				
	ADAM Agrippin Deactivated		Follow-up	MERCIER Thomas	11/02/2025				
	ADAM Archibald Deactivated	+33348900875	Follow-up	MERCIER Thomas	22/09/2025				

If you wish to view a patient's medical record, simply click on the relevant row.

- The '**Profiles and Users**' tab allows you to manage user accounts and permissions within your organisation (creation, modification, deactivation).

It enables you to organise access rights according to your needs, including for non-medical staff.

For a self-employed prescriber working alone, this tab is generally not required, except in cases of collaboration (secretarial support, billing, etc.).

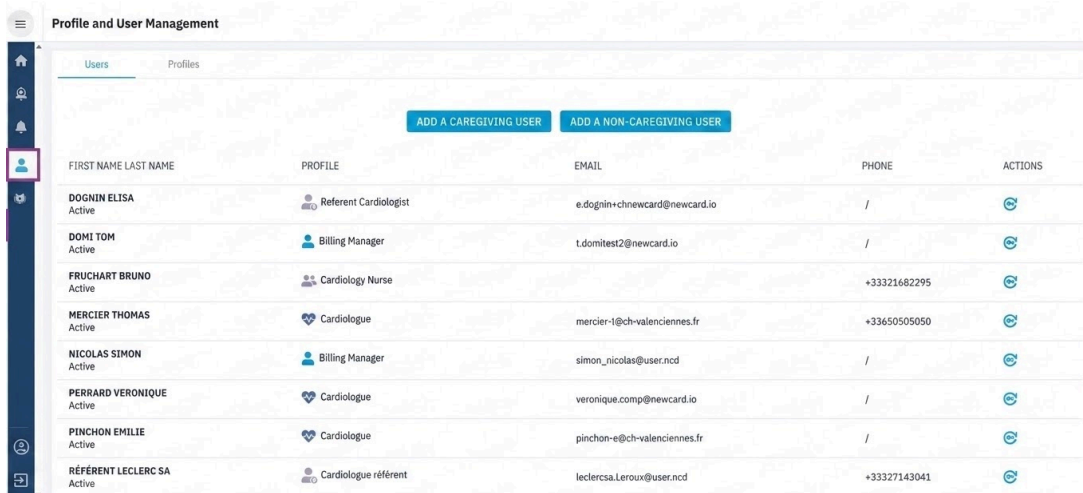FIRST NAME LAST NAME	PROFILE	EMAIL	PHONE	ACTIONS
DOGNIN ELISA Active	Referent Cardiologist	e.dognin+chnewcard@newcard.io	/	
DOMI TOM Active	Billing Manager	t.domitest2@newcard.io	/	
FRUCHART BRUNO Active	Cardiology Nurse		+33321682295	
MERCIER THOMAS Active	Cardiologue	mercier-1@ch-valenciennes.fr	+33650505050	
NICOLAS SIMON Active	Billing Manager	simon_nicolas@user.ncd	/	
PERRARD VERONIQUE Active	Cardiologue	veronique.comp@newcard.io	/	
PINCHON EMILIE Active	Cardiologue	pinchon-e@ch-valenciennes.fr	/	
RÉFÉRENT LECLERC SA Active	Cardiologue référent	leclercsa.Leroux@user.ncd	+33327143041	

→ See section 11 for further information.

10. The 'Help' tab gives you access to the Newcard help section, which is divided into two categories (cardiology and nephrology), where you will find the information you need to help you use the platform.

NEWCARD Dashboard

Good Morning, Dognin Elisa

Medical Alerts [VIEW ALL \(6\)](#)

- 2 DUMONT Liétald 18/03/2026 very high blood pressure
- 1 DUMONT Liétald 18/03/2026 Moderate rise in blood pressure
- 2 MEYER Ascension 18/03/2026 Very high blood pressure

Patients to Renew

Agathe Blanc	+33324701653	RENEW
Héloïse Nicolas	+33691065557	RENEW
Séverine Moreau	+33474589462	RENEW
Adalbold François	+33151411683	RENEW

Patient Statistics

- 585 Patients included
- 216 Patients followed
- 6 Pending training

CH de Valenciennes

Elisa Dognin Lead Cardiologist

Log out

Your demo patient

2	FÉLIX Mehdi DEMO Active	Follow-up
---	-------------------------	-----------

SERVICES • NEWS • LEARN MORE • CAREER • CONTACTS

MEDICAL PLATFORM

CARDIOLOGY HELP CENTER


NEPHROLOGY HELP CENTER

c. Changes to personal data

To access your personal data, click on your **name at the very bottom** of the left-hand sidebar.

This will take you to a tab where you can edit your key details, such as:

- Your username,
- Your password,
- Configuration details,
- Legal and regulatory documents such as the Terms of Service, legal notices, privacy policy and CE marking label,
- And the default units for reading biological data.

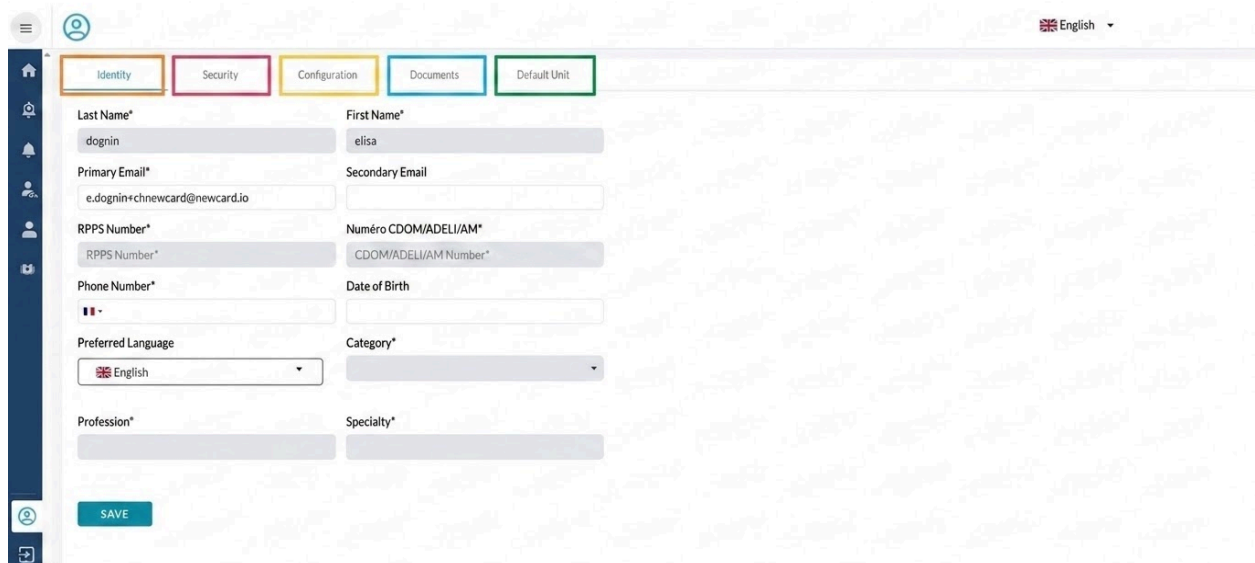The screenshot shows a web application interface for editing a user profile. On the left is a dark blue sidebar with icons for home, notifications, and user profile. The user profile icon at the bottom is highlighted. The main content area has a light blue header with five tabs: Identity (highlighted), Security, Configuration, Documents, and Default Unit. Below the tabs, the form is organized into two columns. The left column contains fields for Last Name* (dognin), Primary Email* (e.dognin@chnewcard.io), RPPS Number* (RPPS Number*), Phone Number* (with a French flag icon), Preferred Language (English), and Profession*. The right column contains fields for First Name* (elisa), Secondary Email, Numéro CDM/ADELI/AM* (CDOM/ADELI/AM Number*), Date of Birth, Category*, and Specialty*. A green SAVE button is located at the bottom left of the form area.

Any changes made to your profile must be **saved** using the button at the bottom of each page.

12. Identity

This tab allows you to edit your account details: **surname, first name, postal address and email address.**

The screenshot shows the 'Profile' page with the 'Identity' tab selected. The page has a sidebar on the left with icons for home, notifications, contacts, and a profile icon. The main content area has a header with 'Profile' and a language dropdown set to 'English'. Below the header are tabs for 'Identity', 'Security', 'Settings', 'Documents', and 'Default Unit'. The 'Identity' tab contains the following fields:

- Last Name*: dognin
- First Name*: elisa
- Primary Email*: e.dognin+chnewcard@newcard.io
- Secondary Email: (empty)
- RPPS Number*: (empty)
- CDOM/ADELI/AM Number*: (empty)
- Phone Number*: (empty)
- Date of Birth: (empty)
- Preferred Language: English (dropdown)
- Category*: (empty)
- Profession*: (empty)
- Speciality*: (empty)

A blue 'SAVE' button is located at the bottom left of the form.

13. Safety

This section allows you to change your password. You will then need to save these changes by clicking on [Enregistrer](#)

The screenshot shows the 'Profile' page with the 'Security' tab selected. The page has the same sidebar and header as the previous screenshot. The 'Security' tab contains the following fields:


- Current password: (masked with dots)
- New password: (empty)
- Confirmation: (empty)

A blue 'SAVE' button is located at the bottom left of the form.

14. Setting

This section allows you to specify the types of alerts you wish to receive, as well as the channel through which you wish to receive them.


For centres, it also allows you to filter the list of patients so that you can view only those for whom you are the prescriber.

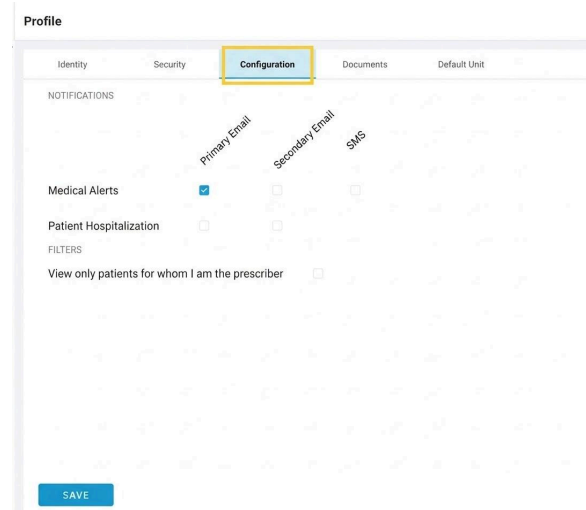The screenshot shows the 'Profile' page with the 'Configuration' tab selected. The page is divided into two main sections: 'NOTIFICATIONS' and 'FILTERS'. In the 'NOTIFICATIONS' section, there are three columns: 'Primary Email', 'Secondary Email', and 'SMS'. Under 'Primary Email', 'Medical Alerts' is checked with a blue square. Under 'Secondary Email', 'Patient Hospitalization' is checked with a blue square. In the 'FILTERS' section, there is a checkbox for 'View only patients for whom I am the prescriber' which is currently unchecked. A 'SAVE' button is located at the bottom left of the configuration area.

15. Documents

This section contains Newcard's regulatory documents: the Terms and Conditions, the Legal Notice, the regulatory label for CE marking, and the Privacy Policy.

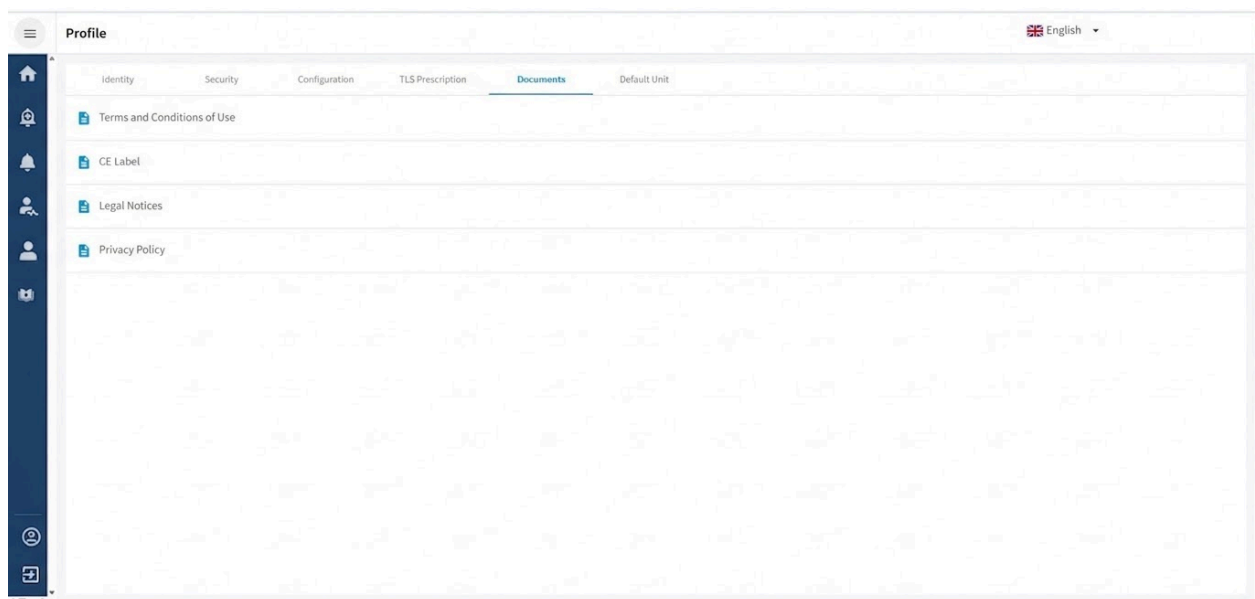The screenshot shows the 'Profile' page with the 'Documents' tab selected. The page displays a list of four documents: 'Terms and Conditions of Use', 'CE Label', 'Legal Notices', and 'Privacy Policy'. Each document is represented by a blue icon and a text label. The 'Documents' tab is highlighted in the top navigation bar. A sidebar with various icons is visible on the left side of the page.

16. Default unit


When entering and viewing your patient's biological data, you can select your preferred default units.

The screenshot shows the 'Profile' page with the 'Default Unit' tab selected. The page is divided into two columns of parameters, each with a dropdown menu for selecting units. The parameters and their current units are:

Parameter	Unit
Creatinine	mg/L
Potassium (K+)	mEq/L
Sodium (Na+)	mEq/L
Bicarbonates or Alkaline Reserve (HCO3-)	mEq/L
Hemoglobin	g/dL
Platelets or Thrombocytes	elements/mm3
Glucose	g/L
Total Cholesterol	g/L
LDLC	g/L
Triglycerides	g/L
Total Bilirubin	mg/L
Calcium	mg/L
Phosphorus	mg/L
24h Proteinuria	mg/24h
Spot Proteinuria	mg/g
Urea	g/L
Magnesium	mg/g
ALR: Albuminuria / Creatininuria Ratio	

d. Access the help section

The 'Help' tab gives you access to the Newcard help section, which is divided into two categories (cardiology and nephrology), where you will find the information you need to help you use the platform.

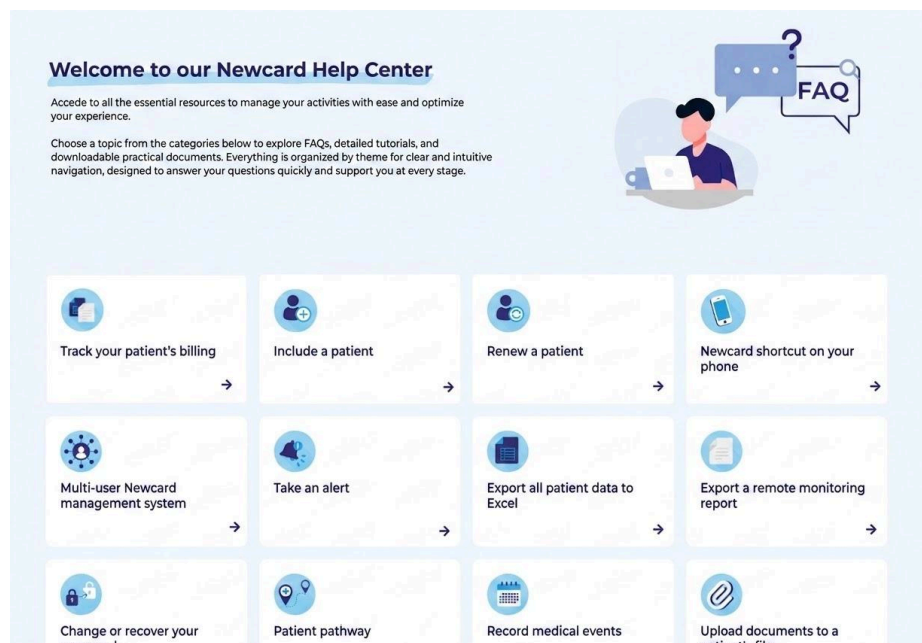Each tab contains an FAQ, a video tutorial and a downloadable document to guide you through the platform's various features.

Include a heart failure patient

Welcome to the page dedicated to including a heart failure patient!
Here you will find all the resources to support you at every stage, with answers to your questions, a tutorial, and documents.

FAQ

Your most frequent questions, simply answered.

What are the eligibility criteria for remote monitoring?


TUTO

A practical video for including a patient.

DOCUMENTS


Find documents to download for a summary of all the stages and additional tips for including a patient.

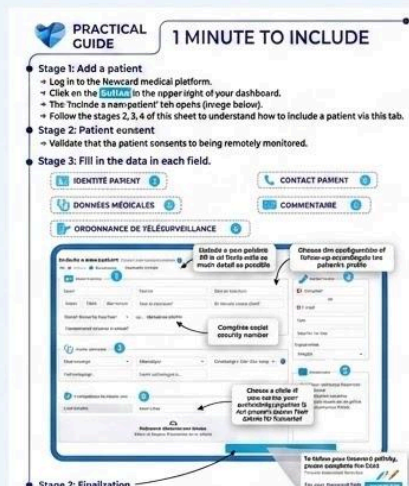5. Inclusion of a patient

The implementation of remote monitoring for a patient is subject to compliance with the eligibility criteria set out in the standard terms and conditions.

To add a patient to the platform, start by clicking the **ADD A PATIENT** button in the top right-hand corner of your dashboard.

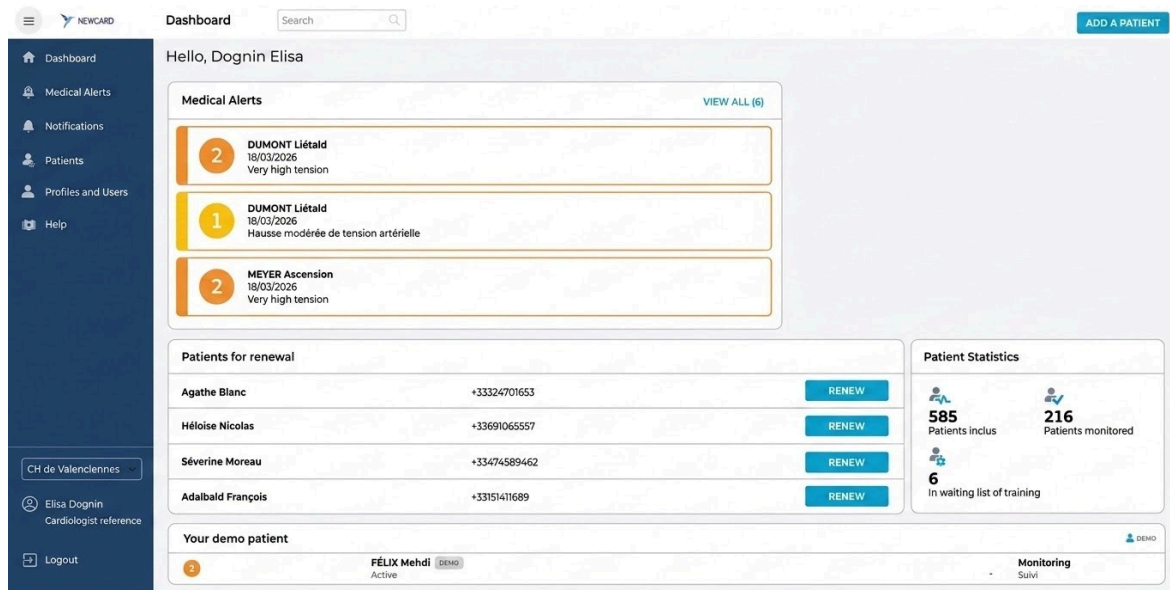Specific features of cardiology inclusion

A window will open on the right-hand side of your screen. There are five fields to complete, all of which are mandatory, in order to finalise the patient's registration.

Enroll a new patient within the framework of the INAMI agreement 

Kit ☐ Balance ☒ Tensiometer ☐ Mobile application

PATIENT IDENTITY

Last Name* First Name* Maiden name

MALE ☐ FEMALE ☐ UNDETERMINED ☐ Date of birth* NISS*

☐ Consent given by the patient*

PATIENT CONTACT

 Mobile*

or

 Landline*

Email

Address, CP, City

Spoken language

English

MEDICAL DATA

Current rhythm* Pacing*

Initial weight (kg)* Systolic function (%)

TELEMONITORING PRESCRIPTION

Barre et Royer FERNANDEZ Gaspard

 **Prescription required for billing.**
Drag and drop the prescription here or click.

COMMENT

Useful information for your patient's journey
Example: hospitalized patient, details about the patient's entourage, contact of the self-employed nurse

ENROLL

 **Please note:** A prescription is not required for patients covered by the INAMI.

Once you have entered this information, click to **INCLUDE** complete the addition.

A confirmation message will appear.

You have two options:

- **“Complete enrolment”**: the patient is immediately added to the Newcard platform and appears in the list of patients awaiting training. Once enrolled, NEWCARD will contact them and provide them with the necessary equipment.
- **“Complete details”**: you are redirected to the patient’s detailed profile

To help Newcard, complete the information.
Otherwise, Newcard will take over.

COMPLETE THE INFORMATION **FINISH THE INCLUSION**

Please complete the details

Patient Edition

1. INCLUSION ✓
2. IDENTITY
3. CONTACT
4. MEDICAL DATA
5. USERS LINKED TO PATIENT

TYPE OF INCLUSION

☒ Inclusion within the framework of the INAMI convention:

- Hospitalization for a flare-up of heart failure;
or
- Currently in NYHA class 2 or higher.

☒ Consent given by the patient

TELEMONITORING PRESCRIPTION *

Patient's telemonitoring prescriber*

FERNANDEZ Gaspard

Telemonitoring prescription to be dated 16/04/2026

Add your telemonitoring prescription
model here.
Drag and drop a file or click.

DOCUMENTS TO DOWNLOAD

☒ IC Prescription for inclusion ☒ Consent Form... ☒ car

KIT

☒ Balance ☒ Tensiometer Mobile application ▼

CANCEL EDIT

If you choose to add further details, a new window will open, allowing you to provide more detailed information about your patient's inclusion criteria. If you choose to add further details, you will see the same five tabs:

- Inclusion
- Identity
- Contact
- Medical data
- Related users

Once you have finished, click 'Edit' to save your changes. Your patient is now added to the list of all your patients and can be easily accessed there.

In the full registration process, you can download the consent form to give to the patient and have them sign it.

DOCUMENTS TO DOWNLOAD

IR Inclusion Prescription

Consent form...

Tab 2 – Identity

Patient Edit

- ✓ 1. INCLUSION
- 2. IDENTITY
- 3. CONTACT
- 4. MEDICAL DATA
- 5. USERS LINKED TO THE PATIENT

IDENTITY

NISS IPP Registration Number

Sex* ☒ MALE ☐ FEMALE ☐ UNDETERMINED

Used name* Used first name*

Birth de naissance* Birth first name(s)* 1st birth first name

Date of birth*

- Please provide: gender, surname, first name, birth name, date of birth, NISS.

Click on

NEXT STEP >

Tab 3 – Patient contact

Patient Edition

- 1. INCLUSION
- 2. IDENTITY
- 3. CONTACT
- 4. MEDICAL DATA
- 5. PATIENT-RELATED USERS

IDENTITY

Mobile* or Landline* Other phone

Address* Email

Spoken language

EMERGENCY CONTACT

Name Telephone

Email Relationship

- Mobile and landline numbers (if available)
- Email address and full postal address
- Contact person in case of absence (with link)

Click on

NEXT STEP >

Tabs 4 & 5 – Medical details

- Please complete the mandatory and optional medical details for the patient's follow-up.

Patient Edition

1. INCLUSION
2. IDENTITY
3. CONTACT
4. MEDICAL DATA
5. PATIENT-RELATED USERS

NYHA / BNP
NYHA
BNP (pg/mL)
Enhanced vigilance
or NT-pro BNP

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Current rhythm Renal insufficiency (clearance in ml/min)
SINUSAL FA TROUBLE SV NON < 30 30 - 60 > 60

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each tab.

Tab 7 – Linked users

- The centre and the referring doctor have already been entered.
- Add other professionals if necessary: registered nurse, another nephrologist, administrative staff, etc.

Patient Edition

1. INCLUSION
2. IDENTITY
3. CONTACT
4. MEDICAL DATA
5. PATIENT-RELATED USERS

PATIENT-RELATED USERS
Healthcare Institution*
Barre et Royer
Patient's Remote Monitoring Prescriber*
FERNANDEZ Gaspard

Healthcare or Administrative Staff
Search for a user

Poirier France Remote Monitoring Nurse in Cardiology X

CANCEL EDIT

Once you've finished, click on


INCLUDE

Your patient is now included in and can be viewed within the list of all your patients.

6. Patient list

The Patients list makes it easy to navigate between patient records. It allows you to quickly find a specific patient, or view a list of patients due for a follow-up appointment, for example.

The **'Patients'** tab will be visible if your profile permissions include **'View heart failure patients'**.

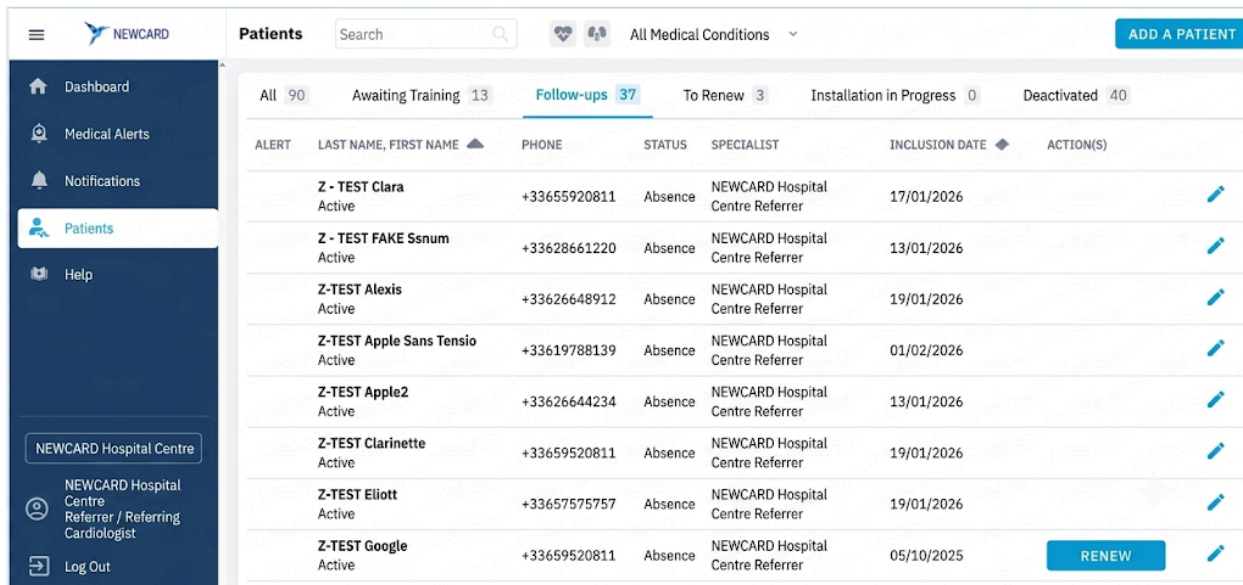ALERT	LAST NAME, FIRST NAME	PHONE	STATUS	SPECIALIST	INCLUSION DATE	ACTION(S)
	Z - TEST Clara	+33655920811	Absence	NEWCARD Hospital Centre Referrer	17/01/2026	
	Z - TEST FAKE Ssnum	+33628661220	Absence	NEWCARD Hospital Centre Referrer	13/01/2026	
	Z-TEST Alexis	+33626648912	Absence	NEWCARD Hospital Centre Referrer	19/01/2026	
	Z-TEST Apple Sans Tensio	+33619788139	Absence	NEWCARD Hospital Centre Referrer	01/02/2026	
	Z-TEST Apple2	+33626644234	Absence	NEWCARD Hospital Centre Referrer	13/01/2026	
	Z-TEST Clarinette	+33659520811	Absence	NEWCARD Hospital Centre Referrer	19/01/2026	
	Z-TEST Elliott	+33657575757	Absence	NEWCARD Hospital Centre Referrer	19/01/2026	
	Z-TEST Google	+33659520811	Absence	NEWCARD Hospital Centre Referrer	05/10/2025	RENEW

The list of patients being monitored is displayed by default:

Other filters are also available, such as: search by patient's surname or first name, filter by medical condition, or filter by stage for patients with kidney failure (stage IV, V or transplant recipients).

Almost all lists are sorted alphabetically by default, but it is entirely possible to sort them by date

of inclusion by clicking on the arrow to the right of the inclusion date

The 'All' patient list includes patients in all different statuses (Setup in progress, Awaiting training, Under monitoring, Deactivated).

The 'Awaiting training' patient list includes all patients who have just been added and have not yet taken any measurements.

The list of patients 'Under installation' includes all patients whose monitoring had been stopped and who are currently being re-equipped.

The list of patients 'Monitored' includes all patients who have taken measurements.

The list of patients 'Deactivated' includes all patients for whom monitoring has been stopped (exclusion from criteria, death, request to stop from the patient, etc.)

The list of patients 'To be renewed' includes all patients requiring renewal: the prescription's expiry date is within 30 days or less, or the prescription's expiry date has already passed. Here,

the list is sorted chronologically: the first entry corresponds to the renewal due earliest. A 'Renew' button is directly available in this list.

If you are authorised to add patients, the button

ADD PATIENT

is also available from this tab.

7. The patient record

The patient record is the key tool for monitoring your patients remotely. On this record, you will find the data regularly transmitted by the patient, the history of their alerts, periods of hospitalisation or absence, and the period of remote monitoring. It is also on the patient record that you can log the resolution of an alert and note any changes to the patient's medical details or treatment. The patient record is displayed differently depending on the patient's condition, particularly in the data access section. Furthermore, the actions available depend on the permissions associated with the user's profile (see 11. Profile and User Management).

a. Administrative information

You will find all the essential information for patient care in the column on the right-hand side of your screen.

Monitoring Period 

INCLUSION IN PROGRESS

15/12/2025

END OF MONITORING

12/06/2026

PRESCRIPTION

Prescription

TERMS & CONDITIONS ACCEPTED

10/12/2025

Documents

LAB REPORTS 

MEDICAL REPORTS 

TREATMENTS 

DIAGNOSTIC IMAGES

OTHERS 

Identity 

STATUS

Absent Active IV

SEX *

Female

You will find:

- The follow-up period, including the start date and the end date of follow-up,
- A copy of the current remote monitoring prescription,
- The date on which the patient accepted the Terms and Conditions,
- Any documents you have attached to your patient's file,
- The patient summary, including their contact details and medical condition,
- The healthcare team, with a summary of the healthcare professionals involved in their remote monitoring,
- The history of remote monitoring prescriptions.

This button  allows you to edit your patient's key details, such as:

- Their personal details,
- Their identity,
- Their contact details,

DATE OF BIRTH *
16/02/1999 - 27 years

INST. REF.
RA-38bf

NEW INST. REF.
b74-38bf

NIR National ID
299020427753884

CODE OF THE AFFILIATED ORGANIZATION
-

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OTHER RELATIONSHIP TYPE
Friend

Supervision

 **HEALTHCARE FACILITY**
PRACTICE OF DR REY FLORIAN

 **TLS PRESCRIBER**
REY Florian

 **INSTALLER**
Official Installer

TLS Prescriptions

Prescription 18/06/2025

Prescription 20/12/2024

Prescription 05/06/2024

Prescription 08/12/2023

b. Alert in progress

You will find details of the alerts at the top of the relevant patient's record. You can use the platform to indicate the actions you have taken to resolve the alert.

→For this, please refer to Section 9: Alert Management. You will then find this information in your patient's medical history.

2

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05/03/2026

PRENDRE EN CHARGE

c. Access to remote monitoring data

If you are a non-clinical user, you do not have access to the various charts containing medical data, but you do have access to this chart, which shows whether or not the patient has taken their measurements:


If you are a healthcare professional:

- Patients with heart failure have a weight graph, a blood pressure graph, a heart rate graph and a symptoms graph.
- Patients with kidney failure (stage IV, stage V and transplant recipients) have a serum creatinine graph, a blood pressure graph, a weight graph and their responses to the two questions asked of the user.

☐ All data

By using the function ☐ All data You can see the time the measurement was taken, as well as any measurements that have been excluded because they were deemed outliers.

On the graphs, for each measurement, you can see whether it was taken manually by the patient or entered manually by Newcard. If nothing is indicated, it means the measurement was taken using a connected device (blood pressure monitor or scales).

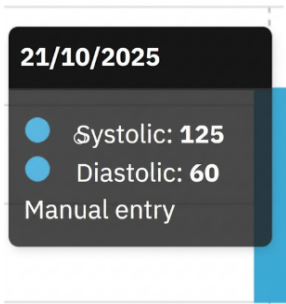On the patient record, you can view the following graphs simultaneously:

- Weight,
- Blood pressure,
- Heart rate
- Symptom monitoring:

Each time a reading is taken, the patient can choose “I feel the same

as usual” or “I feel worse than yesterday”.

- If the patient feels worse than yesterday, a list of symptoms appears:
 - I’m coughing more
 - I’m more out of breath
 - I’m using an extra pillow when I sleep
 - I’m more tired
 - My ankles are more swollen
 - Nothing to report
- If the patient feels as normal, the measurement is complete

So, if the indicator is blue, this means that the patient has replied ‘**I feel as usual**’ or ‘**Nothing to report**’, or has not selected that symptom from the list. If it is red, the patient has selected that particular symptom.

You can adjust the data display period by holding down the left mouse button on one of the graphs. Then, drag the cursor to change the period.

All the graphs are linked, so the period you select will be applied to each of them.

You can also select a predefined period at the top of the patient record, above the first graph

For the weight, blood pressure and heart rate graphs, you can view the date/time and measurement value in a tooltip.

You can adjust the data display period by holding down the left mouse button on one of the graphs. Then, drag the cursor to change the period.

All the graphs are linked, so the period you select will change on each of them.

You can also select a predefined period at the top of the patient record, above the first graph.

d. Visualisation of medical and biological data

Display of medical data

If the data has not been entered, it will not appear in this box.

Medical Data				
17/02/2026				
PRESCRIPTION DATE 17/02/2026	NYHA 2	NT-PROBNP 3460 pg/mL	HEIGHT 186 cm	EJECTION FRACTION 63 %
HYPERTENSION No	DIABETES No	CHRONIC ANEMIA No	RESPIRATORY IMPAIRMENT No	RENAL IMPAIRMENT Clearance 30 - 60
HISTORY OF SUPRAVENTRICULAR ARRHYTHMIA No	CURRENT RHYTHM Sinus	CARDIAC DEVICE Pacemaker	BETA-BLOCKER Yes	OTHER BRADYCARDIC AGENT No
ACE INHIBITORS / ARBS Yes	DIURETICS Yes	MRA (MINERALOCORTICOID RECEPTOR ANTAGONIST) No	SACUBITRIL No	SGLT2 INHIBITOR Yes
SHOW MORE				

Visualisation of biological data

Biochemical test results are displayed at the bottom of the patient record in a table organised by test date. Use the arrows to scroll back through the history of biochemical test results.

Biological Data			+	
ANALYSIS DATE*	<	27/01/2026	31/10/2025	>
TRANSMISSION DATE		27/01/2026	31/10/2025	
SOURCE		APICRYPT	INCLUSION	
CREATININE mg/l		22.7	22.96	
eGFR (CKD-EPI) ml/min/1.73 m ²			26	
eGFR (MDRD) ml/min/1.73 m ²			29	
24H PROTEINURIA mg/24h			200	
POTASSIUM (K+) meq/l		5.39		
SODIUM (NA+) meq/l		138		
SHOW MORE				

Only the data that has been entered is displayed.


Biological data may come from a number of sources:

- **Newcard:** data entered by a Newcard nurse.
- **Patient:** data entered into the mobile app by the patient
- **User:** data entered by a user at the healthcare facility
- **Enrolment:** data entered by the prescriber when the patient is enrolled.


To add biological data, click either the 'Add biological data' button in the action menu or the

button  here.

To edit them, you can use the pencil icon. Don't forget to click Save at the bottom of the page to save your changes.

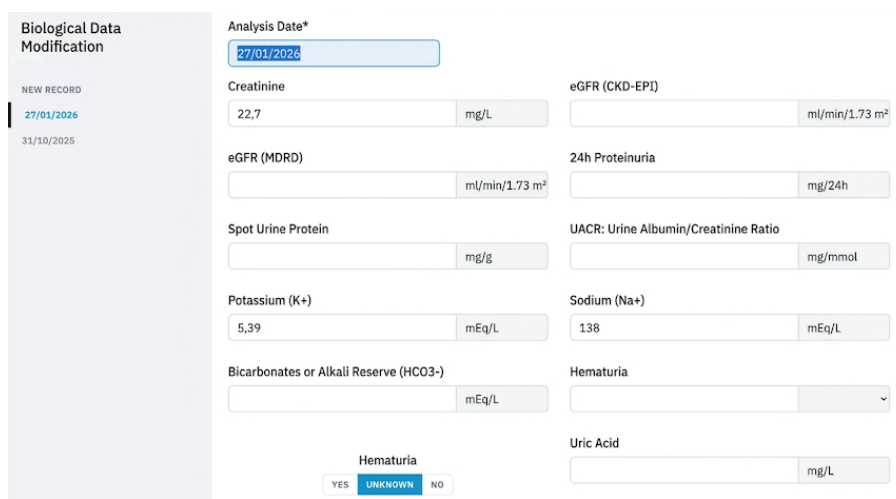You can change the default units in your profile:

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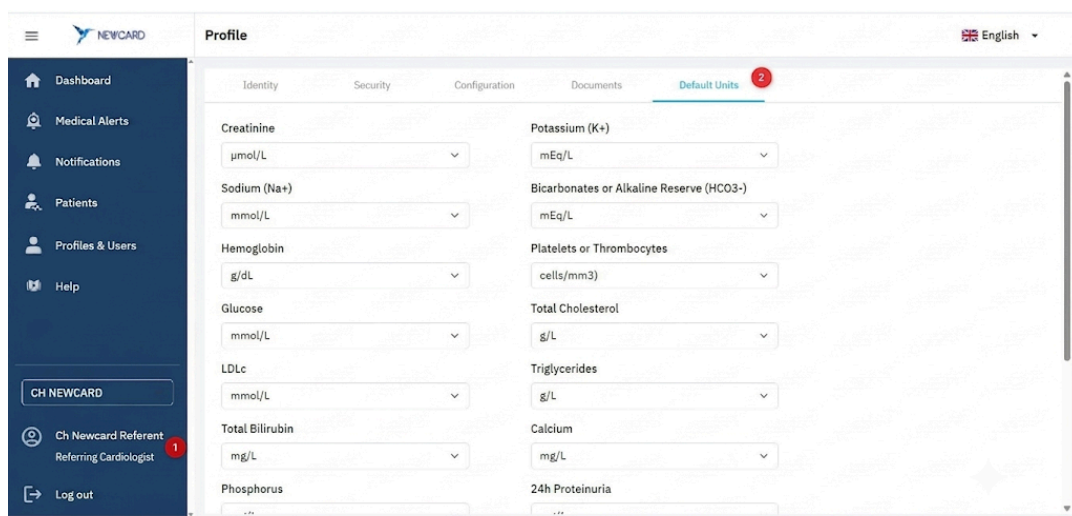f. Patient history

On the patient record, you will be able to add medical comments in the '**Medical Events**' tab. You can note instructions for our nurses, reports following consultations, changes to treatment, etc.

In turn, our team can also add comments for you, such as the postponement of a telephone consultation with a patient or the decision not to trigger a medical alert because the patient was contacted and showed no signs of distress, etc.

This tab also contains the history of medical alerts, technical support sessions, and technical information:

The medical history continues :

- Level 1 and 2 alerts that have been triggered,
- Notes and comments that have been added,
- Changes in the patient's medical status, including the associated dates: hospital admissions and discharges.

Medical events	Therapeutic support	Technical information	Identity information
Comment			
Rejected alerts history			
COMMENT 12/10/2023 11:15 MV left with patient regarding the receipt of blood work results			
13/04/2025 11:55	2 ALERT 2 Systolic blood pressure > than the threshold (3 consecutive measurements)	Completed Monitoring and absence of therapeutic changes	
PLATFORM ANALYSIS 14/04/2025 22:22			
29/10/2024 09:11	1 ALERT 1 Hemoglobin levels between 8 and 10 g/dL		
SHOW MORE			

The patient's history of therapeutic support is as follows :

- The dates and topics of the sessions,
- Comments left by the clinical psychologist responsible for the therapeutic support.
-

Medical Events	Therapeutic Support	Technical Information	Patient Information
Comment			
THERAPEUTIC SESSION 13/12/2022 19:31 Theme: Post-consultation re-evaluation - Conducted on: 13/12/2022 Discussion today with Mr. Coquant. Mr. Coquant is doing well. He saw Dr. Hoffmann last week; everything is fine, blood test results normal. No treatment modifications. Mr. Coquant has a cardio check-up scheduled for next week, and an annual check-up / transplant follow-up in February. He had the flu shot but has not had a new covid booster dose; he is considering it. Mr. Coquant does not plan to travel for the holidays. He will remain cautious regarding covid. Mr. Coquant will see Dr. Hoffmann in 6 months. Good treatment and therapeutic lifestyle self-care compliance. Tablet issues, so the kit was changed. Jessica IDEC Newcard			
Ca . SHOW MORE			
Theme: re-evaluation - Conducted on: 18/08/2022 Discussion with Mr. Coquant's wife.			

The technical history allows you to:

- To find out about any difficulties he may have encountered when using his equipment,
- To monitor the measures taken by Newcard to resolve the issue.

Medical Events	Therapeutic Support	Technical Information	Identity Information
Comment			
ABSENCE 22/09/2025 08:47	Patient hospitalization declared from the platform. From 18/09/2025 to 16/09/2025		
09/09/2025 14:07	T	ALERT T	
ABSENCE 18/08/2025 07:20	Patient hospitalization declared from the platform. From 21/08/2025 to 23/08/2025		
RENEWAL 23/06/2025 09:45	The patient has been renewed		
RENEWAL 18/05/2025 08:45	The patient has been renewed		
SHOW MORE			

You can add a comment to the patient's medical history or treatment history at any time. This provides additional information to the healthcare professionals involved in the remote monitoring process. Once you have written your comment, click on [REGISTER](#) to add it.

Medical Events

Therapeutic Support

Technical Information

Identity Information

Comment

☐ History of rejected alerts

3 Oct 2023
11:15

COMMENT: MV left for the patient regarding the reception of blood test results

12 Apr 2023
13:55

2

ALERT 2

ALERT 2: Systolic blood pressure > norm (3 consecutive measurements)

14 Apr 2023
22:22

Completed

PLATFORM ANALYSIS

PLATFORM ANALYSIS: Surveillance and absence of therapeutic modification

29 Nov 2022
09:11

1

ALERT 1

ALERT 1: Hemoglobin level between 8 and 10 g/dL

SHOW MORE


g. Generating a PDF export of remote monitoring data

You can generate a PDF export of the remote monitoring data by clicking this button in the

patient record  .

This window will open, and you can choose:


- the extraction period
- the sections to be included in the report

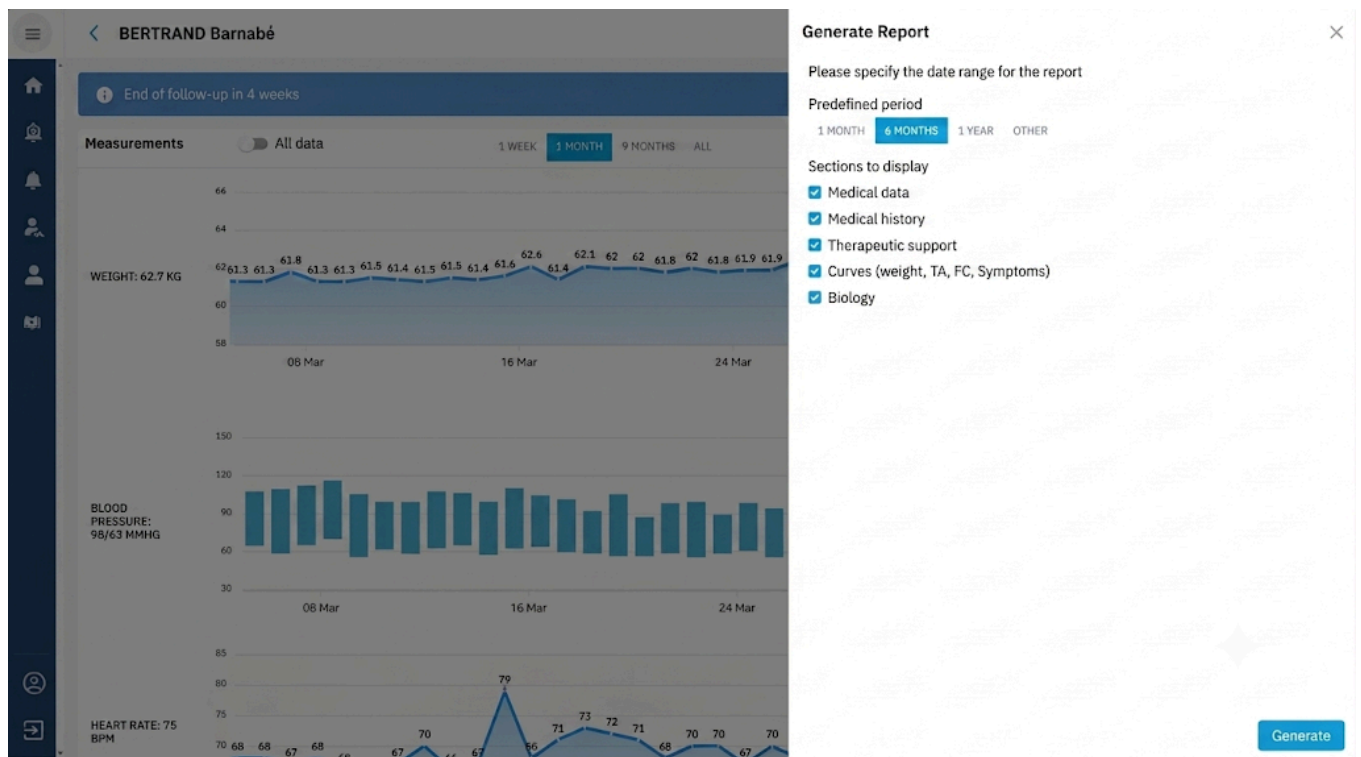By clicking the button **GENERATE**, The PDF export is downloading.

It contains the patient's full medical history: their alerts, clinical data, laboratory results, medical records, adherence, and, for patients with kidney failure, their quality of life.

h. Actions on the patient's profile

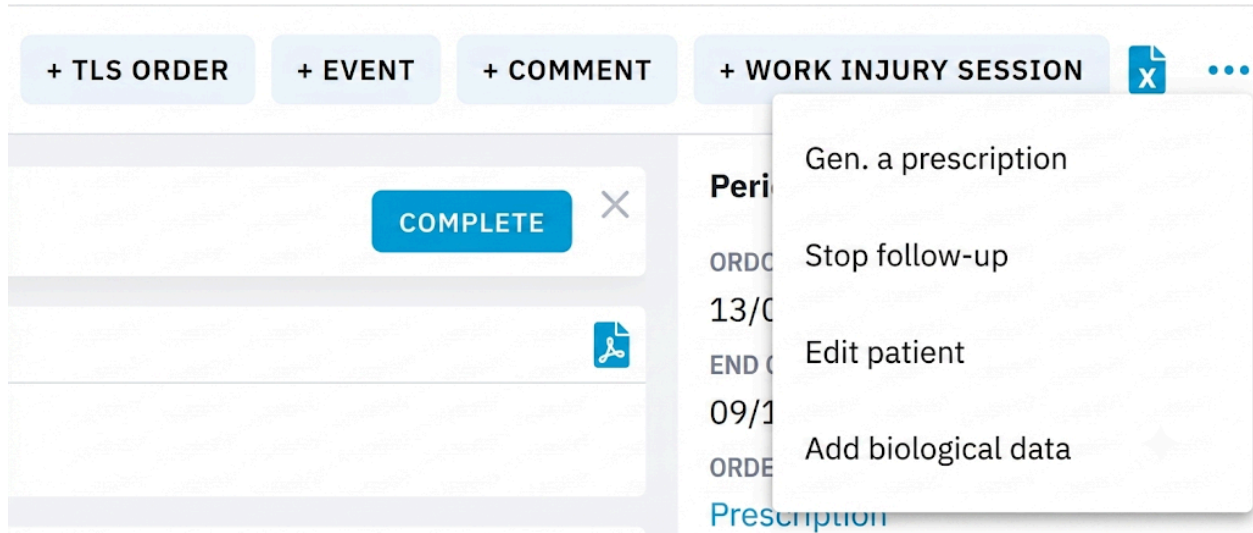At the top of the patient record, you will find several clickable action buttons that provide the following functions:

- The **'Event'** button allows you to indicate whether the patient has been admitted to hospital or is absent for another reason
- The **'Comment'** button allows you to go directly to the bottom of the patient record to add a comment directly to Medical Events, Therapeutic Support or Technical Information, depending on the profile's permissions
- The **'AT SESSION'** button allows you to add a new Therapeutic Support session.
- The X button allows you to export all your patient's data as an Excel document,
- The button  allows you to access a drop-down menu containing several options:
 - Generate a prescription for remote monitoring, medication or laboratory tests
 - Stop remote monitoring
 - Edit the patient's details to amend their administrative or medical information,
 - Add laboratory test results manually,


Each button is subject to a permission

→ See Chapter 11. Profile and User Management.

i. Sending text messages

You can send messages directly to patients with mobile phones via the platform. These are one-way text messages; patients cannot reply to you and your personal number will not be displayed. On the patient's profile page, in the right-hand column, click on the text message icon.

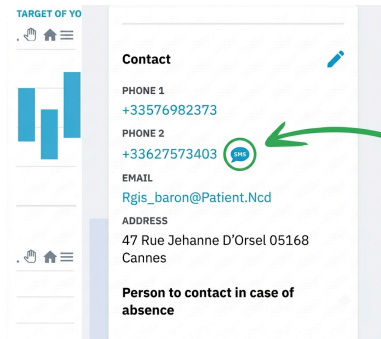A window will appear asking you to enter your message in the text box. The system will automatically include the patient's name and your name as the sender.

You can see the maximum number of characters allowed for this text message in the bottom right-hand corner of the text box.

Then click on

SEND

The patient receives a message from you without your phone number being disclosed.


This feature requires the “**Send SMS**” permission.

8. Alert management

If an alert is triggered, you will be notified via a non-personalised email and, if you wish, via text message as well. The email will contain a link that allows you to access the relevant patient's record directly.

When you log in to the platform, you are taken directly to the dashboard, which displays:

- The severity level of the alert (1 or 2 for patients with heart failure),
- The patient's name,
- The date the alert was triggered,
- The type of data that triggered the alert (weight, blood pressure, heart rate, symptoms or laboratory results, depending on the patient's condition),

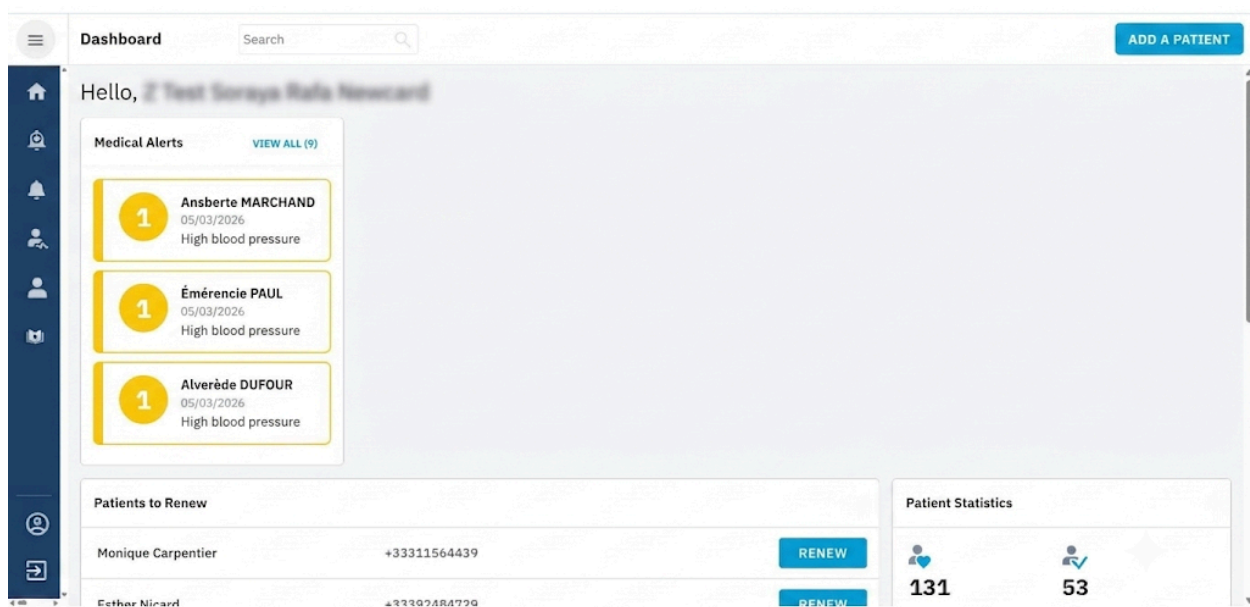
You can also view all alerts directly in the '**Medical Alerts**' tab on the left-hand sidebar.

There are three alert levels :

- **Triggered:** the alert has been sent but has not yet been acknowledged,
- **Viewed:** you have logged into the platform and opened the relevant patient record, but you have not acknowledged the alert,
- **In progress:** you have acknowledged an alert, but it has not yet been resolved.

To view a patient's data when an alert is active, click on their name.

The alert is displayed at the very top of their record. Below that, you will find the patient's data charts.

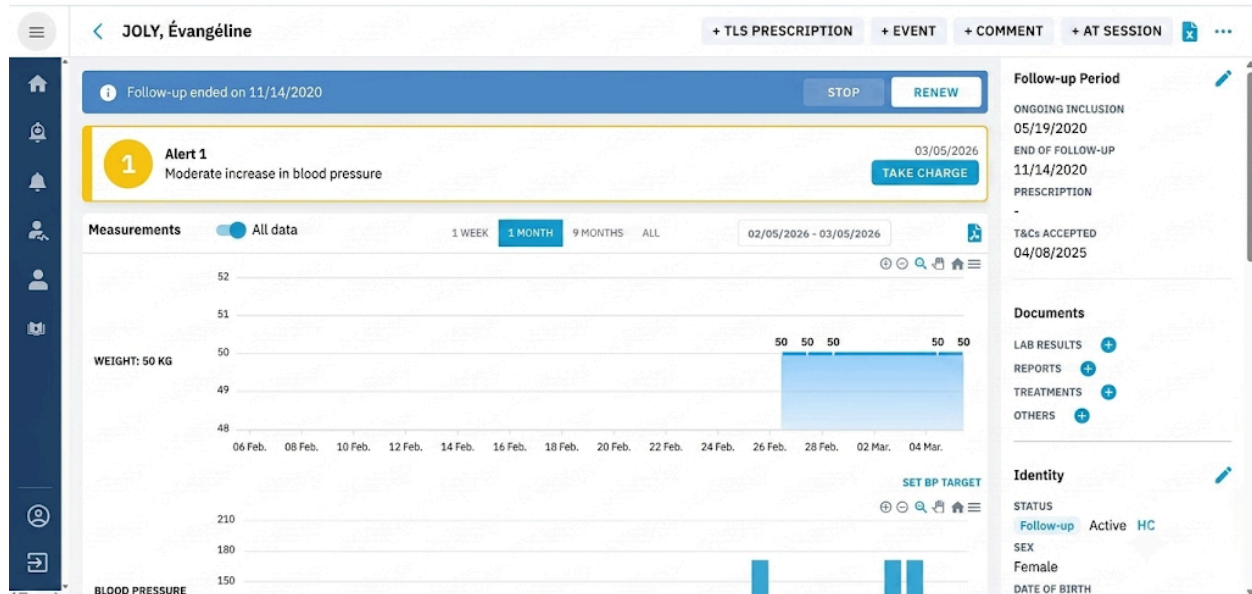Once you have reviewed the alert and the patient's details, you can confirm that you are taking care of the case on the platform by clicking on **TAKE CARE OF**. A window will open, allowing you to specify how this alert should be handled.

Recording patient care activities on the platform makes it easier and more efficient for our team to filter future alerts. On the platform, you must provide the following details:

- The action taken:
 - Analysis on the platform,
 - Patient call,
 - Consultation.
- The type of resolution of the alert, depending on the patient's condition:
 - Spontaneous return to normal,
 - Change in treatment
 - Hospitalisation
 - Patient cannot be contacted
 - etc.
- A comment on how this alert was handled
- An assessment of the alert's relevance (1 star: not very relevant; 5 stars: very relevant).

If the alert has not been resolved, click on '**Active alert**': you can specify the action to be taken later. However, Until an alert has been confirmed, no other alert can be triggered.

If the alert has been resolved, click '**Alert closed**'.

Once you have closed the alert, it will appear in the patient's medical history, as well as in the "Medical Alerts" tab, under the "Closed" sub-tab.

User access to viewing and managing alerts is subject to the '**View medical alerts**' and '**Manage medical alerts**' permissions.

These permissions are only available to healthcare professional profiles (see section 11. Profile and user management).

If a user only has access to "View medical alerts", they can see the alert banner and the alert history, but they do not have access to the "Take charge" button, which allows them to manage alerts.

9. Notification management

In the Notifications tab, you will find three types of information:

- adding a comment about a patient;
- adding a blood test result;
- adding a document.

The aim is to notify you as soon as the Newcard team enters new information about a patient, and vice versa.

Please note: once a notification is marked as ‘resolved’, it disappears.

Therefore, if several users are associated with the same organisation, the notification will no longer be visible to others once it has been resolved.

This process also applies to notifications relating to the addition of documents and laboratory results to a patient’s record.

ALERT	NAME, FIRST NAME	TYPE	NOTIFICATION DATE	ACTIONS
	ANDRE Judicaël	Biological Report	07/04/2026	PROCESS
	ARNAUD Constance	Document	07/04/2026	PROCESS
	ADAM alexandre	Comments	07/04/2026	PROCESS

Showing item 1 to 0 of 0

10. Profile and user management

User and profile management allows each user within your organisation to have permissions tailored to their role. In particular, this enables you to create non-medical accounts securely. Finally, you can manage these access rights independently. With multi-user management, you can:

- *Restrict access to patient data to authorised professionals,*
- *Provide the right information to the right people*
- *Tailor the system to your needs*
- *Deliver a consistent experience for all users*

This means that each institution can tailor its profiles and permissions to suit its own organisational structure.

USER GUIDE: SETTING UP USERS AND THEIR PERMISSIONS

Go to the new **Profiles and Users** tab to set up your teams.

Profiles and Users

Profile and User Management

Users Profiles

1 ADD CAREGIVER USER 2 ADD NON-CAREGIVER USER

LAST NAME - FIRST NAME	PROFILE	EMAIL	TELEPHONE	ACTIONS
BENESTEAL ALBERT ALST	Disabled	g.benesteal@newcard.ic	/	
BENETTON GEORGES ALST	Disabled	g.beneton@newcard.ic	/	
DELANGUE CARINE ALST	Disabled	c.delange12@newcard.ic	/	
DOONEN ELISA ALST	Reference User	d.degrete@mwac.din	+EXTERUC41	
DONS TOM ALST	Billing Manager	tubanfoat1@gmail.com	/	
DUFORT PORQUE ALST	Billing Manager	lmckstaly@gmail.com	/	
PRETTY OLIVER ALST	Cardiologist	labert@gmail.com	/	

Regenerate a password

For users who have lost their password, simply click on the **'Regenerate password'** option. A new automatic email will be sent.

Profile and User Management

Users Profiles

1 ADD CAREGIVER USER 2 ADD NON-CAREGIVER USER

LAST NAME - FIRST NAME	PROFILE	EMAIL	TELEPHONE	ACTIONS
BENESTEAL ALBERT ALST	Disabled	g.benesteal@newcard.ic	/	
BENETTON GEORGES ALST	Disabled	g.beneton@newcard.ic	/	
DELANGUE CARINE ALST	Disabled	c.delange12@newcard.ic	/	
DOONEN ELISA ALST	Reference User	a.degrete@newcard.ic	+2XTSHOC41	
DONS TOM ALST	Billing Manager	tulanologist@macol.io	/	
DUFORT PORQUE ALST	Billing Manager	lunskater@gmail.com	/	
PRETTY OLIVER ALST	Cardiologist	labert@gmail.com	/	

1

ADD A HEALTHCARE USER

- 1 Complete the requested information
(last name, first name, email, phone...)

Add A User

Identity

☒ Foreign user

Affiliated Facility
Lucas EURL

Unique PS Number

Last Name*
Doe

First Name*
John

Primary Email*
j.doe@mail.com

Secondary Email
j.doe@mail.com

Phone
FR

CDOM/ADEL/AM Number

Language
FR French

Date of birth
08/11/1995

REGISTER

2 REGISTER

- The user automatically receives their password by email, at the address indicated during user creation.

Attention: If none of the existing profiles meet your needs, you can create a new one (see page 6).

Profile permissions apply to all associated users:

- to assign specific permissions to a user, a new dedicated profile must be created!

2

ADD A NON-CLINICAL USER

- 1 Choose a profile from the non-clinical profiles (e.g.: *administrative, billing*)
- 2 Complete the requested information

The screenshot shows the 'Add A User' interface. On the left is a sidebar with navigation links: Dashboard, Medical Alerts, Notifications, Patients, Billing, Profiles and Users, and Help. Below these are user-specific links: SIMBAIRCO, Elise Bagon (Role: Nursing Cardiologist), and Logout. The main form area is titled 'Add A User' and has a sub-tab 'Identity'. It contains several input fields: 'Affiliated Facility' (pre-filled with 'CH Newcard'), 'Last Name*' (pre-filled with 'Doe'), 'Primary Email*', 'Phone' (with a dropdown menu showing 'French'), and 'Date of birth'. A 'Profile' dropdown menu is open, displaying a list of non-clinical profiles: 'Billing Manager', 'Disabled', 'Billing Manager' (highlighted in blue), 'Billing Profile 7', 'Secretariat', and 'TLS Secretary'. At the bottom of the form is a blue 'REGISTER' button, which is pointed to by a red arrow.

3 REGISTER

- The user automatically receives their password by email, at the address indicated during the user creation.

Attention : If none of the existing profiles meet your needs, you can create a new one (see page 6).

The permissions of a profile apply to all associated users:

- to assign specific permissions to a user, a new dedicated profile must be created!

Deactivate a user

Assign the **“Disabled”** profile to remove all access rights.

Referent Newcard

HEALTHCARE WORKER

Active

Identity

Configuration

Documents

Default Unit

Search the RPPS directory * To update the user's identity, please go through the RPPS directory.

Last Name*

référent

First Name*

Newcard

Affiliated Facility

Profile

Cardiologist Referent

Cardiologist

Cardiology Specialized Nurse

Cardiology Telesurveillance Nurse

Disabled

Cardiologist Referent

Dietitian

Primary Email*

ref.cardio@newcard.io

Phone


 +33327143041

CDOM/ADELI/AM Number

RPPS Number

Edit a profile (caregiver or non-caregiver)

- 1 Select the profile linked to the user
- 2 Check/uncheck permissions according to their role

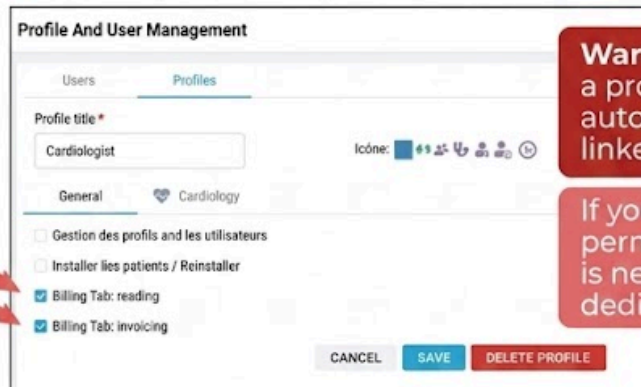Warning: any modification to a profile's permissions will automatically apply to all users linked to this profile.


If you want to assign specific permissions to a single user, it is necessary to create a new dedicated profile.


- 3 Save changes

 **Tip:** review permissions regularly to remain compliant with the department's needs.

Create a new profile (caregiver or non-caregiver)

- 1 Give a title to the profile
- 2 Choose a representative icon
- 3 Check the profile permissions

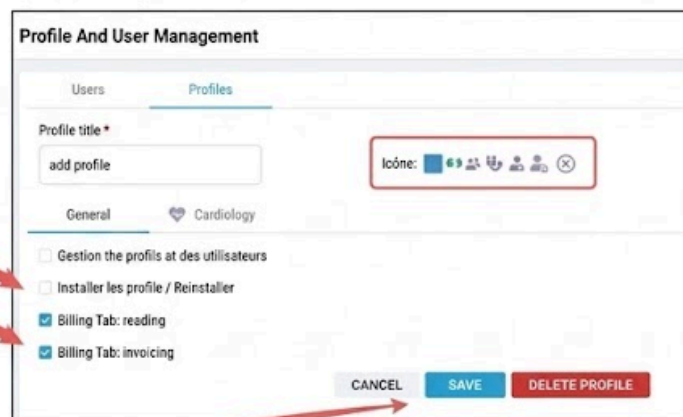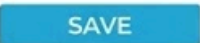- 4 

→ Useful for adapting the platform to new roles or services.


Configure notifications and alerts

Each user can:

- Choose the types of alerts they receive
- Define the patients concerned and the reception channel

Tip: Configure alerts upon patient inclusion so you don't miss anything.

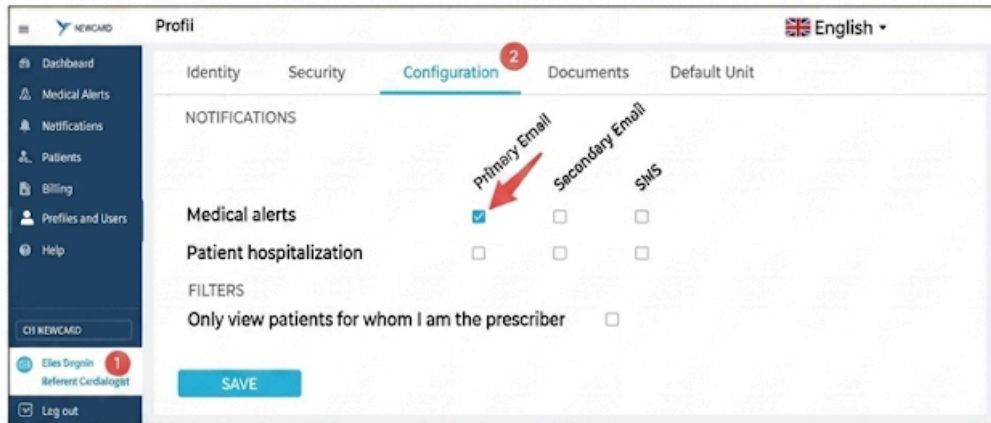NEWCARD Profil English

Identity Security **Configuration** Documents Default Unit

NOTIFICATIONS

Medical alerts ☒ Primary Email ☐ Secondary Email ☐ SMS

Patient hospitalization ☐ ☐ ☐

FILTERS


Only view patients for whom I am the prescriber ☐

SAVE

Filter patient view

In centers with multiple prescribers, all patients may be visible.

This filter allows you to display only the patients for whom you are the prescriber, for clearer navigation.

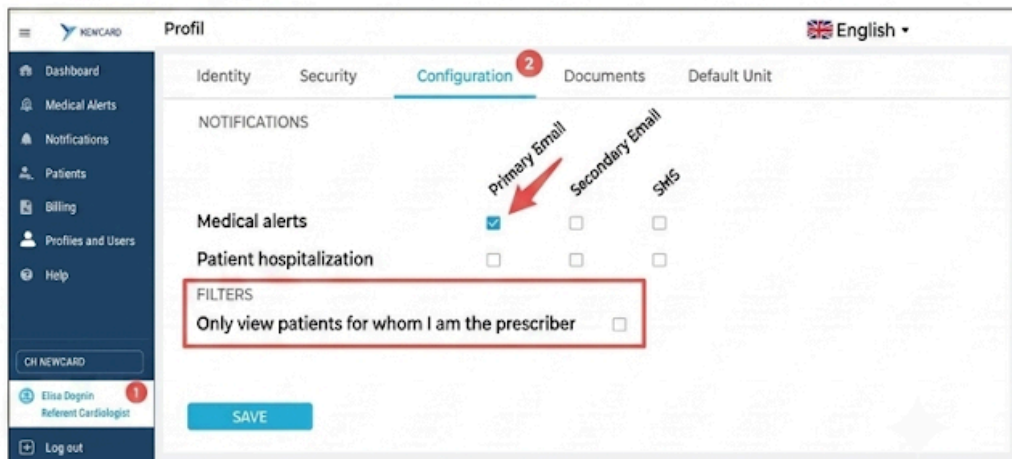NEWCARD Profil English

Identity Security **Configuration** Documents Default Unit

NOTIFICATIONS

Medical alerts ☒ Primary Email ☐ Secondary Email ☐ SMS

Patient hospitalization ☐ ☐ ☐

FILTERS

Only view patients for whom I am the prescriber ☐

SAVE

Only users who have the **‘Profile and User Management’** permission in their profile can access the **‘Profiles and Users’** tab. We recommend that the organisation’s designated contact person be granted this permission.

For a sole self-employed prescriber, this tab is not necessary; they will have all rights by default and will not need to configure other people’s rights.

List of rights and associated features:

Permission Types	Permission Name	PS Law	PNS Law	Definition
General	Profile and user management	x	x	The "Users and profiles" tab is active. You can create users and profiles and modify them within your facility.
General	Install / Reinstall patients	x	x	Ability to install the patient with this account using the email and password in the patient app, display of the list of patients to be installed, and ability to click on the "Reinstall" button in the patient file (provided there is access to the patient file with the activate cardio visualization or nephro visualization permissions).

Cardiology	Activate HF patients visualization	x	x	The "Patients" tab is visible displaying only the cardiology patients of the facility. View-only access to the patient file.
Cardiology	Can prescribe remote monitoring	x		The user account name appears in the drop-down list of remote monitoring prescribers during patient inclusion and editing. It is possible to use the pre-signed prescription under this prescriber's name. If the user can prescribe RM (Remote Monitoring), they automatically have all other rights related to the pathology.
Cardiology	Include patient	x		The "Add a patient" button is available on the dashboard and on the Patients list.

Cardiology	Re-include patient	x		The list of patients to be renewed is available on the dashboard. The "Renew" action is available on the dashboard and in the patient lists; the "To be renewed" list is visible, the renewal banner is visible, the RENEW button is available in the patient lists, and the "Renew" action is visible in the context menu (3 little dots).
Cardiology	Stop patient	x		Access to the "Stop follow-up" button in the patient file's context menu.
Cardiology	Declare an event (Hospitalized, absent,...)	x		Access to the "+ Event" button, which allows entering an absence, a hospitalization, or a return from hospitalization.

Cardiology	Generate and upload a prescription (remote monitoring, biology, therapeutics ...)	x		Access to the "Generate a prescription" button in the context menu + ability to upload a prescription in the semesters screen and during patient inclusion/editing.
Cardiology	View medical alerts	x		View the Medical alerts tab, see the list of medical alerts in this list, on the patient file, or on the dashboard, but without the ability to manage them.
Cardiology	Handle medical alerts	x		View the Medical alerts tab, see the list of medical alerts in this list, on the patient file, or on the dashboard; the "handle" button is available.
Cardiology	Send SMS	x	x	Access to the small "SMS" bubble, allowing an SMS to be sent from the platform to the mobile phone.

Cardiology	Consult documents and prescriptions	x	x	Grants read access to documents (document section) and prescriptions (prescriptions section, access to the semesters tab) + the Documents space is visible in the notification area (on the dashboard and in the Notifications tab).
Cardiology	Import a document to the patient file	x	x	Ability to add a document, open the document viewer, and upload a document (document section) and prescriptions (prescriptions section, access to the semesters tab) + the Documents space is visible in the notification area (on the dashboard and in the Notifications tab).

Cardiology	Modify patient administrative information	x	x	Access to the "Edit patient" button in the patient file's context menu, access to the pencil button to edit the patient in the Identity, Contact, and Supervision sections. For administrative editing, only the Identity, Contact, and Patient-linked users tabs are available for editing.
Cardiology	Modify patient medical information	x		Access to the "Edit patient" button in the patient file's context menu, access to the pencil button to edit the patient in the follow-up period and Medical data sections. For medical editing, only the Inclusion and Medical data tabs are available for editing. Access to the button to Add biological data.

Cardiology	View medical events and add comments	x		Access to Comments in the Notifications area on the dashboard and in the Notifications tab. Access to the patient file's medical events with the ability to add a comment directly in the section or with the "+Comment" button.
Cardiology	View technical information and add comments	x	x	Access to Technical information in the patient file with the ability to add a comment directly in the section or with the "+Comment" button.
Cardiology	Visualization of therapeutic support	x		Access to Therapeutic Support sessions. You can read the sessions but you cannot add a TS session via the "+TS Session" button (only reading the TS section is available). However, you can write a comment.

Cardiology	Add therapeutic support sessions	x		Access to Therapeutic Support sessions. You can read the sessions and add a comment directly in the section or with the "+Comment" button. You can also add a TS session via the "+TS Session" button.
------------	----------------------------------	---	--	--

11. Contact us

If you have any questions, please contact NEWCARD, a simplified joint-stock company, on **+33 3 20 14 71 14**

The registered office is located at 236 Avenue Clément Ader, 59118 Wambrechies, France.

MPR-NTU-PLTF-PS-BEL-00

This user guide (MPR-NTU-PLTF-PS-BEL-FR-00), created on 15 April 2026, may be subject to change depending on platform updates. Should this be the case, we will be sure to inform you of any changes.

The ‘1 Minute pour mon Coeur’ (1 Minute for my Heart) solution is a Class I medical device intended for the remote monitoring of heart failure. Manufactured by Newcard, this regulated healthcare product bears the CE mark in accordance with the relevant regulations. Newcard also provides patients with equipment to complement the remote monitoring software, in order to facilitate data collection and analysis. This device is not an emergency system and must be used in accordance with the instructions in the user manual. It is 100% covered by the French National Health Service.

Head office: Wambrechies - 821 766 961 - Registered with the Lille Métropole B Trade and Companies Register.

Date of incorporation: 18 May 2026 - Image: Newcard - Do not dispose of on public land - MPR-NTU-PLTF-PS-BEL-00

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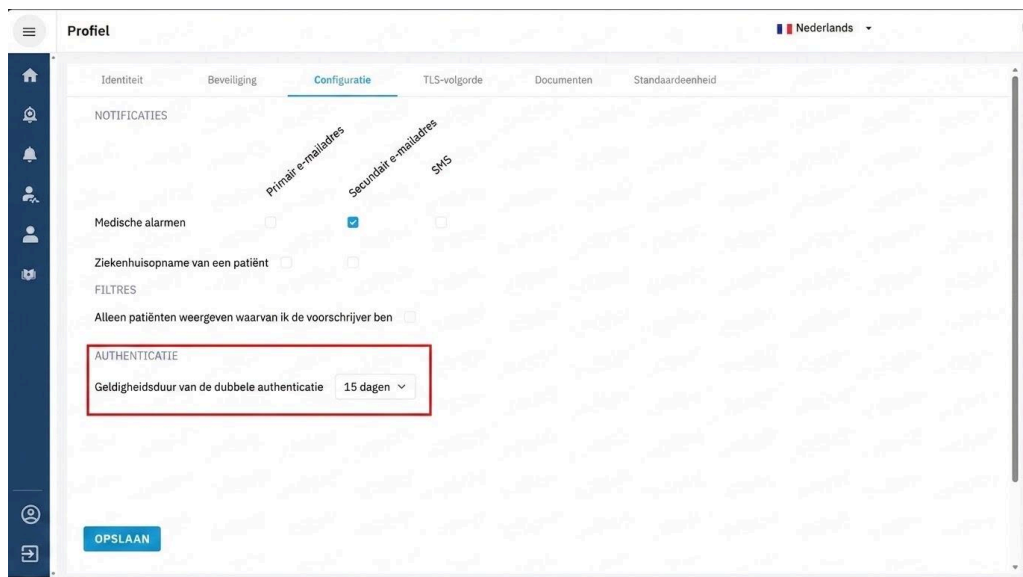
The image shows a screenshot of the 'Profiel' (Profile) settings page in the Newcard platform. The page has a dark blue sidebar with icons for home, notifications, settings, and other functions. The main content area is titled 'Profiel' and has a language dropdown set to 'Nederlands'. There are tabs for 'Identiteit', 'Beveiliging', 'Configuratie', 'TLS-volgorde', 'Documenten', and 'Standaardeenheid'. The 'Configuratie' tab is selected. Under 'NOTIFICATIES', there are options for 'Medische alarmen', 'Ziekenhuisopname van een patiënt', and 'FILTRES'. The 'FILTRES' section has a checkbox 'Alleen patiënten weergeven waarvan ik de voorschrijver ben'. The 'AUTHENTICATIE' section is highlighted with a red box and contains a dropdown menu for 'Geldigheidsduur van de dubbele authenticatie' set to '15 dagen'. A blue button labeled 'OPSLAAN' is at the bottom.

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The screenshot shows a user profile page with a sidebar and a main content area. The sidebar contains icons for home, notifications, and profile. The main content area is titled 'Profil' and has a language dropdown set to 'Nederlands'. Below the title are five tabs: 'Identiteit' (highlighted in orange), 'Beveiliging' (pink), 'Configuratie' (yellow), 'Documenten' (blue), and 'Standaardeenheid' (green). The 'Identiteit' tab is active, showing a form with fields for 'Naam*' (Last name, value: 'dognin'), 'Voornaam*' (First name, value: 'elisa'), 'Hoofd-e-mailadres*' (Primary email, value: 'e.dognin@newcard.io'), 'Secundair e-mailadres' (Secondary email), 'Numéro RPPS*', 'CDOM/ADEL/AM-nummer*', 'Telefoon*' (Phone), 'Geboortedatum' (Date of birth), 'Voorkeurstaal' (Preferred language, dropdown set to 'Nederlands'), 'Categorie*' (Category), 'Beroep*' (Profession), and 'Specialisatie*' (Specialization). A blue 'OPSLAAN' (Save) button is at the bottom left of the form.

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
The screenshot shows the 'Profiel' (Profile) page with the 'Identiteit' (Identity) tab selected. The page contains various input fields for personal information, including Name, First Name, Main Email Address, Secondary Email Address, RPPS Number, CDOM/ADELIA Number, Telephone, Date of Birth, Preferred Language, Category, Profession, and Specialization. A blue 'OPSLAAN' (Save) button is visible at the bottom left.

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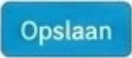The screenshot shows the 'Profiel' (Profile) page with the 'Beveiliging' (Security) tab selected. The page contains input fields for 'Huidig wachtwoord' (Current password), 'Nieuw wachtwoord' (New password), and 'Bevestiging' (Confirmation). Each field has a toggle icon to show or hide the password. A blue 'OPSLAAN' (Save) button is visible at the bottom left.

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The screenshot shows the 'Profiel' (Profile) page in the Newcard system. The 'Standaard eenheid' (Standard unit) tab is selected and highlighted with a green box. The page displays two columns of dropdown menus for selecting units for various biological parameters. The parameters and their current units are as follows:

Parameter	Unit
Creatinine	mg/L
Natrium (Na ⁺)	meq/L
Hemoglobine	g/dL
Glucose	g/L
LDLc	g/L
Totaal bilirubine	mg/L
Fosfor	mg/L
Protéinurie sur échantillon	mg/g
Magnesium	
Kalium (K ⁺)	mEq/L
Bicarbonaat of alkalireserve (HCO ₃ ⁻)	mEq/L
Bloedplaatjes of trombocyten	elementen/mm ³
Totaal cholesterol	g/L
Triglycériden	g/L
Calcium	mg/L
Eiwit over 24u	mg/24h
Ureum	g/L
UAC: Urine-Albumine-Creatinine ratio	

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Go to the new '**Profiles and Users**' tab to configure your teams.

Profiles and Users

LAST NAME / FIRST NAME	PEOPLE	EMAIL	TELEPHONE	ACTIONS
DRISTEAU ALBERT	Desired	albert@dristeau.com	/	[icon]
DRISTEAU SERGE	Desired	serge@dristeau.com	/	[icon]
PELLAGO CAROTTE	Desired	carotte@pellago.com	/	[icon]
BOONEN ELSA	Cardiology external	elsa.boonen@cardio.com	- CARDIOLOGY	[icon]
BOONEN TOM	Responsible to Sofia elise	tom.boonen@sofia.com	/	[icon]
SUPOT MONIQUE	Responsible to Sofia elise	monique@supot.com	/	[icon]
BLAIS OLIVIER	Cardiology	olivier@blais.com	/	[icon]

Regenerate a password

For users who have lost their password, simply click on the option '**Regenerate the password**'. A new automated email will be sent.

LAST NAME / FIRST NAME	PEOPLE	EMAIL	TELEPHONE	ACTIONS
DRISTEAU ALBERT	Desired	albert@dristeau.com	/	[icon]
DRISTEAU ALBERT	Desired	albert@dristeau.com	/	[icon]
DRISTEAU ALBERT	Desired	albert@dristeau.com	/	[icon]
BOONEN ELSA	Cardiology external	elsa.boonen@cardio.com	- CARDIOLOGY	[icon]
BOONEN TOM	Responsible to Sofia elise	tom.boonen@sofia.com	/	[icon]
SUPOT MONIQUE	Responsible to Sofia elise	monique@supot.com	/	[icon]
BLAIS OLIVIER	Cardiology	olivier@blais.com	/	[icon]

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