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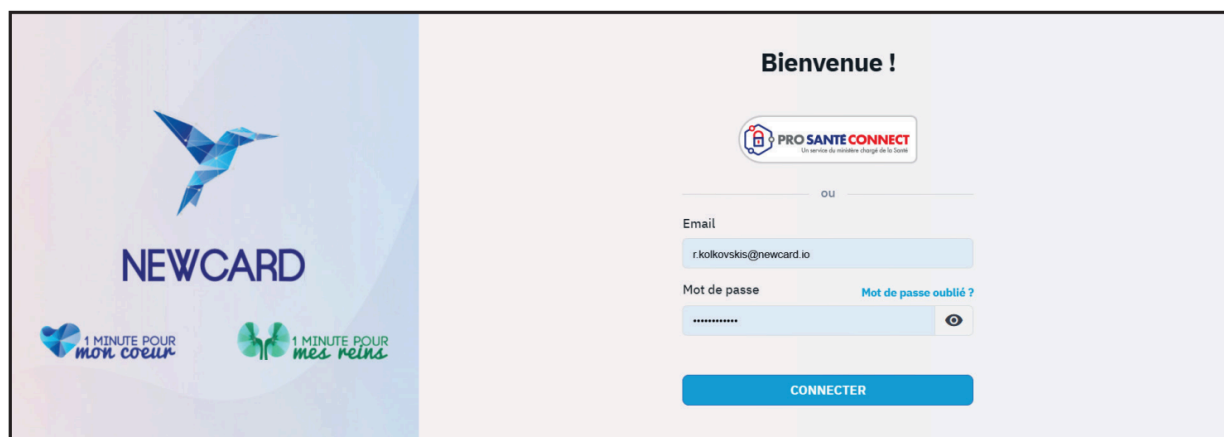
The login screen for Newcard features a blue and white background with a stylized bird logo. On the left, there are two circular icons: one with a heart and the text '1 MINUTE POUR mon cœur', and another with kidneys and the text '1 MINUTE POUR mes reins'. The main heading is 'Bienvenue !' followed by the 'PRO SANTE CONNECT' logo. Below this, there is a login form with fields for 'Email' (containing 'r.kolkovskis@newcard.io') and 'Mot de passe' (masked with dots). A 'Mot de passe oublié ?' link is next to the password field. A blue 'CONNECTER' button is at the bottom.

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The screenshot displays the NEWCARD application interface. On the left is a dark blue sidebar with a menu containing: Tableau De Bord, Alertes Médicales, Notifications, Patients, Facturation, Profils Et Utilisateurs, and Aide. Below the menu, there is a section for 'CH Valenciennes' and a user profile for 'Elisa Dognin, Cardiologue référent', with a 'Déconnexion' button at the bottom. The main content area has a search bar at the top right with the text 'Chercher' and a magnifying glass icon. Below the search bar, there is a header for 'in Elisa' and a button 'VOIR TOUT (6)'. The main content area is divided into three sections: 'IT Liétald' (026) with 'très haute' status, 'IT Liétald' (026) with 'modérée de tension artérielle' status, and 'Ascension' (026) with 'très haute' status. At the bottom, there is a table with four rows of patient data, each with a 'RENOUVELER' button. The table data is as follows:

Phone Number	Action
+33324701653	RENOUVELER
+33691065557	RENOUVELER
+33474589462	RENOUVELER
+33151411689	RENOUVELER

On the right side of the main content area, there is a 'Statistiques patients' section with two columns of data:

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The screenshot shows the NEWCARD dashboard interface. The sidebar on the left contains navigation links: Tableau de bord, Alertes Médicales (highlighted), Notifications, Patients, Facturation, Profils et Utilisateurs, and Aide. The main content area displays the following sections:

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AJOUTER UN PATIENTS

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The screenshot shows the NEWCARD dashboard. At the top right, there is a button labeled "AJOUTER UN PATIENTS" highlighted with a red box. The dashboard includes a sidebar with navigation options like "Tableau de bord", "Alertes Médicales", "Notifications", "Patients", "Facturation", "Profil et Utilisateurs", and "Aide". The main content area displays "Bonjour, Dognin Elisa" and "Alertes médicales" with three entries for patients DUMONT Liétald and MEYER Ascension. Below this is a "Patients à renouveler" table with four rows of patient information and "RENOUVELER" buttons. On the right, there is a "Statistiques patients" section showing metrics like "585 Patients inclus", "216 Patients suivis", "6 En attente de formation", and "94766€ En attente de facturation". At the bottom, there is a "Votre patient démo" section for FÉLIX Mehdi.

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The screenshot shows the 'Tableau De Bord' (Dashboard) of a medical platform. At the top, there's a search bar and a button 'AJOUTER UN PATIENT'. The main content area displays a greeting 'Bonjour, Z Test Soraya Rafa Newcard'. Below this, there's a section for 'Alertes médicales' with a link 'VOIR TOUT (9)'. Three medical alerts are listed, each with a severity level '1', a patient name, a date, and a description: 'MARCHAND Ansberte' (05/03/2026, Tension haute), 'PAUL Émérance' (05/03/2026, Tension haute), and 'DUFOUR Alverède' (05/03/2026, Tension haute). Below the alerts, there's a section for 'Patients à renouveler' with a table listing two patients: 'Monique carpentier' (+33311564439) and 'Eithar nicard' (+33309484729), each with a 'RENOUVELER' button. To the right, there's a 'Statistiques patients' section showing two icons with counts 131 and 53.

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The screenshot shows the 'Réseaux' (Network) configuration page of the jFSE application. The interface includes a sidebar with navigation links: 'Lecteurs', 'Réseaux' (selected), 'Maintenance', 'Cartographie', and 'Librairies'. The main content area is titled 'Réseaux' and contains three sections: 'Serveur JFSE', 'Serveurs locaux', and 'Proxy'. The 'Serveur JFSE' section has input fields for 'Adresse IP / Hostname' (filled with 'lecteur-newcard.jfse.fr'), 'Port' (filled with '443'), 'ID Communication' (filled with 'de8674b4-e1eb-4209-9b2c-8f5e9cadaf75'), and 'Mode de compatibilité' (filled with '2'). Below these are three toggle switches: 'Chiffrement des données échangées' (off), 'Paramétrage de l'adresse IP automatique du serveur FSE via internet' (off), and 'Affichage des animations lors d'une action sur le lecteur' (on). The 'Serveurs locaux' section has three toggle switches: 'Serveur web local' (on, port 8083), 'Service XML-RPC local' (off, port 8080), and 'Serveur web local TOMCAT' (on, port 8895). The 'Proxy' section has input fields for 'Host' (filled with 'Host') and 'Port' (filled with '0'). A button 'Enregistrer les modifications' is located in the top right corner. The sidebar also displays version information (Daemon: 3.28, Serveur: 3.7.0-DV-04-SNAPSHOT, Socle: 1.40.14) and environment details (OS: Windows 10 64 bits, Java: OpenJDK Runtime Environment 11.0.2).

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1. Welcome message

Your practice/department has expressed an interest in the remote monitoring of patients using NEWCARD.

This user guide has been designed to guide you, step by step, through the process of getting started with the platform.

Remote medical monitoring enables the early detection of signs that a chronic condition is worsening, allowing complications to be anticipated and care to be adapted. It also allows for closer monitoring of patients, thanks to the regular transmission of data from their homes. This system is not a substitute for emergency care.

At Newcard, we place particular emphasis on personalised support, offering tailored care for patients, carers and healthcare professionals.

Throughout your experience, you can count on the expertise and availability of our paramedical, administrative and sales teams.

About the guide

This guide forms an integral part of the '1 Minute for My Heart' software (base IUD-ID 3770034348MPC6G) and '1 Minute for My Kidneys' software (base IUD-ID 3770034348MPR7E). It provides a detailed overview of the software, its settings and how to use it. Its purpose is to help users understand how their remote monitoring works and to guide them through the process.

If needed, a printed version can be provided on request at the following address:
<https://newcard.io/contacts/> or cardio@newcard.io.

There are no additional costs involved and the user manual will be delivered to you at the latest within seven calendar days of receipt of the request.

2. The basics of remote telemonitoring

Destination of solutions

- *'1 Minute pour mon coeur (1 Minute for my heart) application is a software-as-a-medical-device solution designed for the remote monitoring of patients with chronic heart failure through the collection, analysis, storage and visualisation of health data.'*

- *'1 Minute Pour mes Reins (1 Minute for my kidneys) application is a software-as-a-medical-device solution designed for the remote monitoring of patients with chronic kidney disease through the collection, analysis, storage and visualisation of health data. The patient application allows the collection of your data for processing by the Newcard algorithm.'*

Logging in to the platform

You can access the Newcard platform to register and monitor your patients via remote monitoring.

<https://app.newcard.io>

Your username is the email address you provided to us.

You can also log in via **Pro Santé Connect**.

Indication: eligibility criteria

Telemonitoring is a service covered in full by the National Health Service (Assurance Maladie).

Specifics in cardiology (*heart failure*):

Patients with chronic heart failure, regardless of the type and etiology of the disease, who meet one of the following two conditions:

- Hospitalisation within the last 12 months due to an episode of heart failure; ;
- OR currently in NYHA class 2 or higher with elevated natriuretic peptide levels (BNP > 100 pg/mL or NT-proBNP > 100 pg/mL).

Specific in nephrology (*chronic kidney disease*):

Patients suffering from chronic kidney failure and eligible for remote monitoring:

- **Stage IV:** GFR between 15 and 30 ml/min/1.73 m²
- **Stage V:** GFR < 15 ml/min/1.73 m²
- **Kidney transplant patients:** transplant must have been performed more than 3 months ago

Contraindication: non-eligibility criteria (common to both pathologies)

Patients in any of the following situations are considered ineligible for remote monitoring:

- Physical or mental impossibility for the patient or caregiver to use the digital medical telemonitoring MD and/or its collection accessories according to the doctor's judgment wishing to include the patient in the medical telemonitoring project;
- Refusal by the patient to provide data necessary for monitoring the actual use of the DMN and for obtaining individualized or national real-life results;
- Refusal of the patient to have therapeutic support.

Specific in nephrology :

- Patients with stage IV chronic kidney disease who are not yet on dialysis,

Description of labelling/ explanation of symbols

	Medical device
	Unique identification code of the device
	Name and address of legal manufacturer
	Batch number
	CE mark
	view the user guide
	Product code
	Packaging and items not to be thrown away
	Fragile
	Keep away from the heat

Precautions

1. **This device is for personal use only.**
2. **The kit is a loan granted to you by Newcard for the duration of your monitoring.** At the end of your follow-up, you are required to make the device available to you as soon as possible and in the state in which you received it. We also thank you for keeping the original boxes and cardboard, which will be useful for returning the equipment.
3. **Confidentiality and security:** Login credentials and passwords are strictly personal and must never be shared with a third party. The user is asked to use a secure password, enable two-factor authentication (2FA), and ensure the protection of their equipment in order to guarantee data confidentiality and security. Any suspicious or unauthorized use of the account must be reported promptly to the relevant support team.
4. **The Newcard service is not an emergency service;** it is not intended to be a substitute for interventions by healthcare professionals. Recorded data or alerts are not read and interpreted instantly, but only during the days and hours set by the monitoring operator. If you have a medical issue or your condition suddenly worsens, please contact a healthcare professional or call 15.
5. **Preventive and regular maintenance operations:** Newcard may modify or adapt the NEWCARD application or the General Terms of Use. In this case, an update from you may be necessary. In accordance with the TOS, this update must take place as soon as possible.

Intended Users

The 1 Minute pour mon Coeur (1 Minute for my heart) platform in Newcard is intended to be used by:

- Health professionals involved in the treatment and follow-up of chronic heart failure.

The 1 Minute pour mes Reins (1 Minute for my kidneys) platform from Newcard is intended to be used by:

- Health professionals involved in the treatment and follow-up related to chronic kidney failure.

The patient app is intended to be used by:

- Eligible chronic heart failure patients.
- Eligible chronic kidney failure patients.

Target patient populations

- Male and female patients with chronic heart failure.

- Patients, men and women, patients suffering from chronic kidney failure.

Clinical benefits

The purpose of telemonitoring is to enable a medical professional to remotely interpret the data necessary for the medical follow-up of a patient and, if necessary, to make decisions relating to their care. Using the Newcard software offers the following clinical benefits:

- **Improved patient survival (reduced mortality)**
(among patients with heart failure followed up via remote monitoring)
- **Improved patient independence**
(as part of remote monitoring for heart failure and kidney failure)

Operating principle and mode of action

An internet connection is required to use the software.

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Browser	Compatible Version
Google Chrome	≥ 126.0
Hardware	
Desktop	Minimum 17-inch screen with 1920 x 1080 resolution (Full HD)
Laptop	Minimum 15-inch screen with 1920 x 1200 resolution
Tablets	10-inch screen size with 800 x 1280 resolution
Smartphones	• Screen size 6.5 with 412 x 915 resolution: for Android • Screen size 6.5 with resolution 414 x 896: for iPhone.

Mobile operating system	Compatible version
Android	≥ 10
iOs	≥ 10
Hardware	
Tablets	10-inch screen size with 800 x 1280 resolution
Smartphones	• Screen size 6.5 with 412 x 915 resolution: for Android • Screen size 6.5 with resolution 414 x 896: for iPhone.

In case of material vigilance

Any incident or risk of a serious incident that has led to or may lead to the death or serious deterioration in the health status of a patient, user or third party involving Newcard's services

must be reported without delay to Newcard at the following email address: contact@newcard.io and to the competent authority (in France, to the National Agency for the Safety of Medicines and Health Products - ANSM).

Residual risks

The software has no direct physical contact with the patient and poses no mechanical, electrical or biological risk. The residual risks identified are of a technical or software nature; although rare, they may, in certain situations, have an indirect impact on patient care. The main residual risks associated with the use of this software are described below.

1. Appropriate technical and organisational measures have been put in place to ensure the platform remains available. Despite these controls, there remains a residual risk that the platform may be temporarily unavailable until it is restored (in the event of an internet connection failure, a technical error, etc.).
2. The platform incorporates mechanisms designed to ensure the reliability of the alerts it generates. Nevertheless, a residual risk of alert errors—including, in particular, false positives, false negatives, duplicates or alerts not being displayed—cannot be entirely ruled out.
3. Despite the controls in place, errors may occur in the transmission of patient measurements (failure to transmit data, or data that is inaccurate or unreliable).
4. Newcard has implemented appropriate cybersecurity measures to ensure the security of all data transmitted. However, a residual risk of loss, alteration or unauthorised access to this data cannot be entirely ruled out. In any event, Newcard is committed to safeguarding users' privacy and complying with applicable data protection regulations, in particular the GDPR.

Patient inclusion

Patients are included directly through the platform.

Once a patient has been included, Newcard takes care of:

- contacting the patient;
- organising training on how to use the equipment;
- monitoring their compliance (ensuring measurements are taken regularly);
- managing any technical issues that may arise.

Alert management

Remote monitoring data and triggered alerts are directly accessible on the platform.

When an alert is triggered:

- you receive an email;
- a link takes you to the patient's record;
- you decide on the appropriate course of action.

Billing

Basic fee:

- **€28 per patient per month – TVB code** (excluding the surcharge for heart failure*)

***Surcharge on the basic fee for remote monitoring of heart failure:**

The adjustments provided for on the basis of the criteria set out in points 1 and 2 of Article L.162-54 and in Section II of Article R. 162-95 of the Social Security Code, applicable to the rates for Level 2 provider packages, resulting in the following increased rates:

- **Level 1: €56/month – code TVB2**

- Patients aged 80 or over with at least one of the comorbidities listed in point 2 of section III of this article;
- Patients under the age of 80 with at least two of the comorbidities listed in point 2 of section III of this article;
- Patients who have been hospitalized within the last 30 days for cardiac decompensation.

- **Level 2: €70/month – code TVB2.5**

- Patients aged 80 or over with at least two of the comorbidities listed in point 2 of section III of this article;
- Patients under the age of 80 with at least three of the comorbidities listed in point 2 of section III of this article;
- Patients with specific conditions or clinical histories:
 - Those who have suffered cardiogenic shock within the last 6 months;
 - Those on the waiting list for heart transplantation;
 - Cardiorenal syndrome and awaiting dialysis;
 - Refractory heart failure with a life expectancy of less than one year;
 - Restrictive cardiomyopathy, discharged from hospital with heart failure corresponding to at least NYHA Class II;
 - Cardiac amyloidosis, discharged from hospital with heart failure corresponding to at least NYHA Class II.

List of comorbidities:

- Cancer undergoing systemic treatment or radiotherapy
- Chronic kidney disease (only if stage IV or V → eGFR < 45 ml/min/1.73 m²)

- Severe anaemia or severe iron deficiency
- Severe malnutrition associated with sideroblastic anaemia or severe protein-energy malnutrition

3. First-time login to the platform

When you log in to the Newcard remote monitoring platform for the first time, simply enter your login details.


You will then be asked to accept the **Terms and Conditions of Use** (T&Cs).

By accepting these T&Cs, you acknowledge that you have read and understood how the service works and you agree to:

- keep your login details confidential;
- use the platform for its intended purpose;
- not interfere with the proper functioning of the service.

The T&Cs are available at any time in the '**Documents**' tab of your profile.

You can also log in via the '**Pro Santé Connect**' portal.

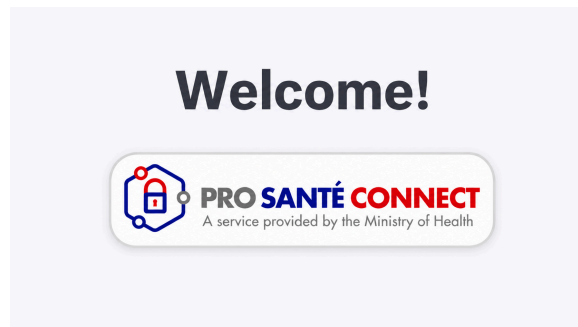Two-factor authentication

Two-factor authentication (2FA) enhances security when accessing Newcard by adding an additional verification step to the password.


In light of the growing risks of phishing, data breaches and compromised login credentials, this measure helps to protect the health data processed on the platform.

It is required upon first login and then every 15 days by default for all users.

Login and two-factor authentication steps

1. Login


Log in with your email address and password.

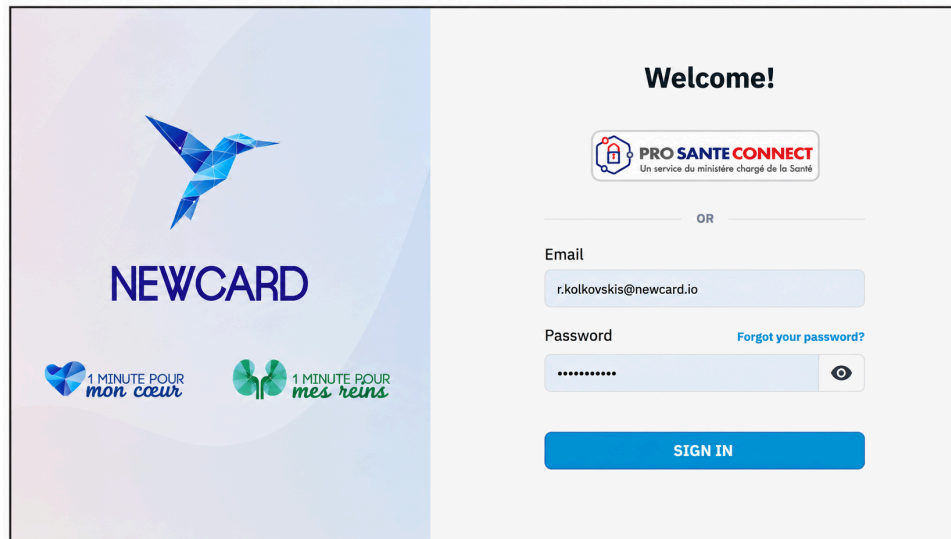The login screen is split into two panels. The left panel has a light blue background with a stylized blue bird logo at the top, the word 'NEWCARD' in large blue letters, and two logos at the bottom: '1 MINUTE POUR mon cœur' with a heart icon and '1 MINUTE POUR mes reins' with a kidney icon. The right panel has a light gray background. At the top, it says 'Welcome!' in bold. Below that is the 'PRO SANTE CONNECT' logo with the tagline 'Un service du ministère chargé de la Santé'. A horizontal line with 'OR' in the center separates this from the login fields. There are two input fields: 'Email' with the text 'r.kolkovskis@newcard.io' and 'Password' with a masked password '*****'. A link 'Forgot your password?' is next to the password field. At the bottom of the right panel is a blue button labeled 'SIGN IN'.

2. Verification request


A screen will appear asking you to enter a verification code.

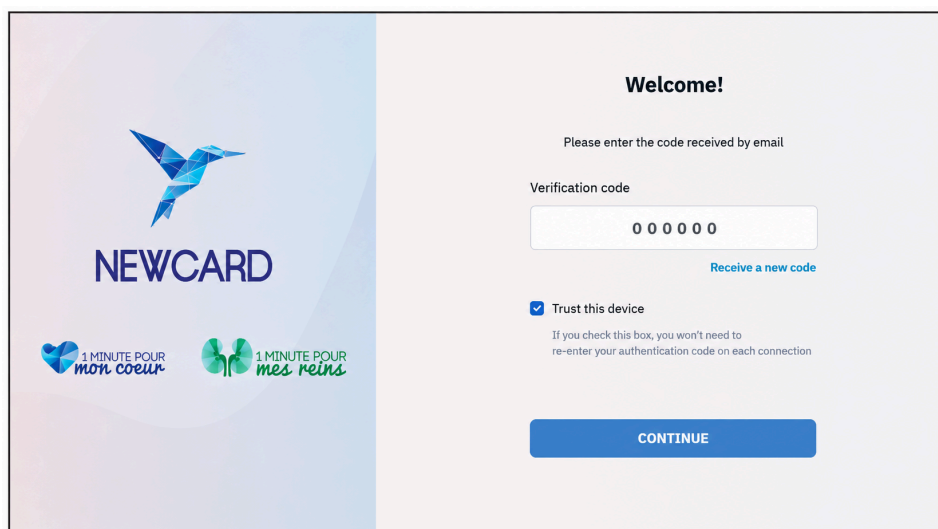The verification screen is split into two panels. The left panel is identical to the login screen. The right panel has a light gray background. At the top, it says 'Welcome!' in bold. Below that, it says 'Please enter the code received by email'. There is a 'Verification code' label above a text input field containing '0 0 0 0 0'. A link 'Receive a new code' is below the input field. Below that is a checkbox labeled 'Trust this device' which is checked. Underneath the checkbox, it says 'If you check this box, you won't need to re-enter your authentication code on each connection'. At the bottom of the right panel is a blue button labeled 'CONTINUE'.


3. Receiving the code

A 6-digit code will be sent to you by email.

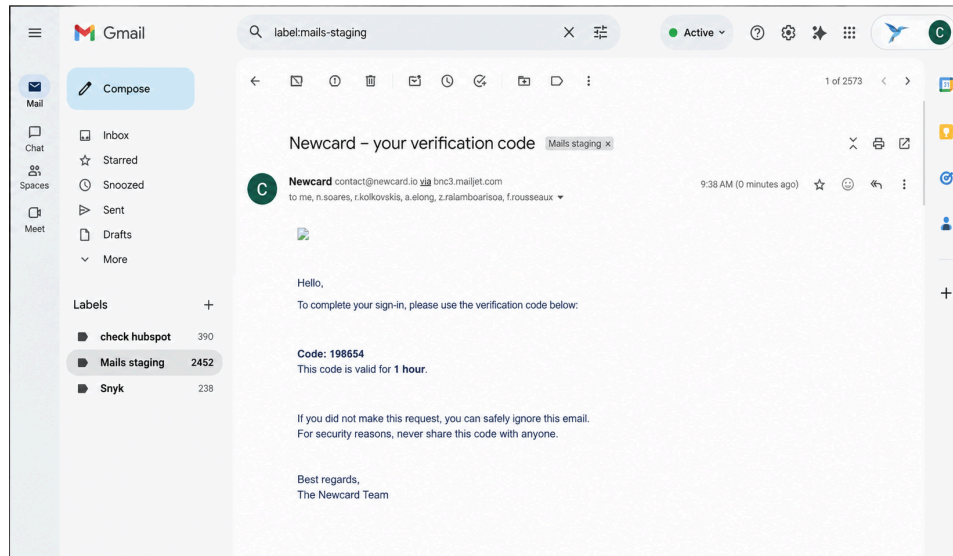4. Enter the code

Enter this code to confirm your access.

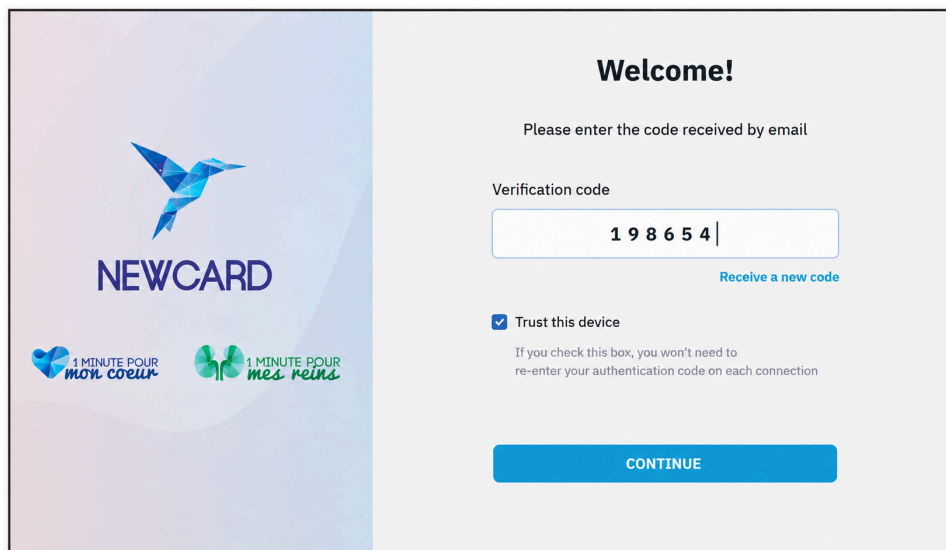5. Access to the platform

Your login has been verified.
You can now access the Newcard platform securely.

Verification frequency

By default, two-factor authentication is required every 15 days. You can adjust this frequency (from 1 to 30 days) in your profile, under the “Settings” section.

The screenshot shows a web interface for a user's profile. At the top, there's a 'Profile' header with a language dropdown set to 'English'. Below this is a navigation bar with tabs: Identity, Security, Configuration (active), TLS order, Documents, and Default unit. A left sidebar contains icons for home, notifications, profile, and other settings. The main content area is divided into two sections. The 'NOTIFICATIONS' section has a table with columns for 'Medical alerts', 'Patient hospitalization', and three notification methods: 'Primary email', 'Secondary email', and 'SMS'. In the 'Medical alerts' row, 'Secondary email' is selected with a blue checkmark. In the 'Patient hospitalization' row, no method is selected. Below this is the 'AUTHENTICATION' section, which contains a label 'Validity period for two-factor authentication' and a dropdown menu currently showing '15 days'. A red rectangle highlights this dropdown. At the bottom left of the configuration area is a blue 'SAVE' button.

	Primary email	Secondary email	SMS
Medical alerts	☐	☒	☐
Patient hospitalization	☐	☐	☐

AUTHENTICATION

Validity period for two-factor authentication: 15 days ▼


SAVE

4. General overview of the platform

The platform is a digital tool that allows you to access your patients' data and view alerts. You will also find a list of healthcare professionals involved in the remote monitoring of your patients.

a. Dashboard

On the platform's login page, you will be taken to the remote monitoring platform's dashboard.

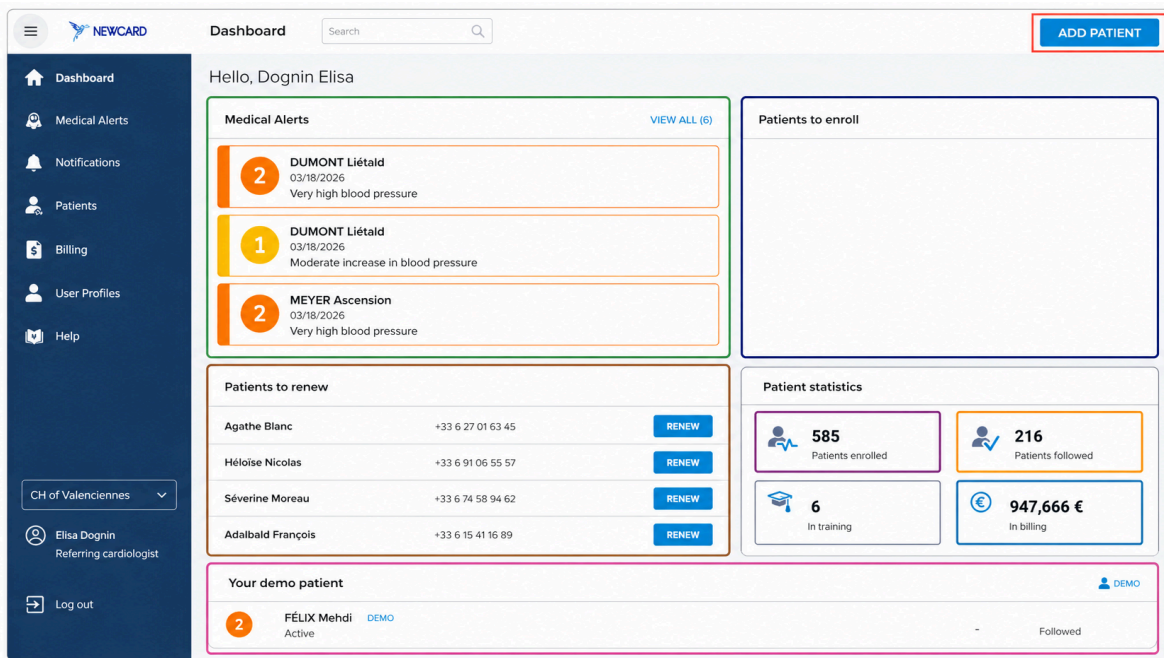This allows you to access the following directly:

- ✓ The button to **add a patient**,
- ✓ Current **medical alerts**,
- ✓ The list of **patients to be registered**,
- ✓ The number of **patients included**,
- ✓ The number of **patients being monitored**,
- ✓ The list of **patients due for renewal**,
- ✓ The total amount pending **invoicing** for your remote monitoring services,
- ✓ The list of your **'demo' patients**.

All these items are clickable and will take you directly to the relevant pages. These items can also be found in the categories on the sidebar of your platform.


Please note: some sections are conditional and will appear or be hidden depending on the permissions assigned to each user. As a result, the dashboard may vary from one user to another.

→ For further details on managing accounts and user rights, please refer to section 11 of this guide.

In this example, we outline a configuration with full access (a user with all permissions).

b. Feature sidebar

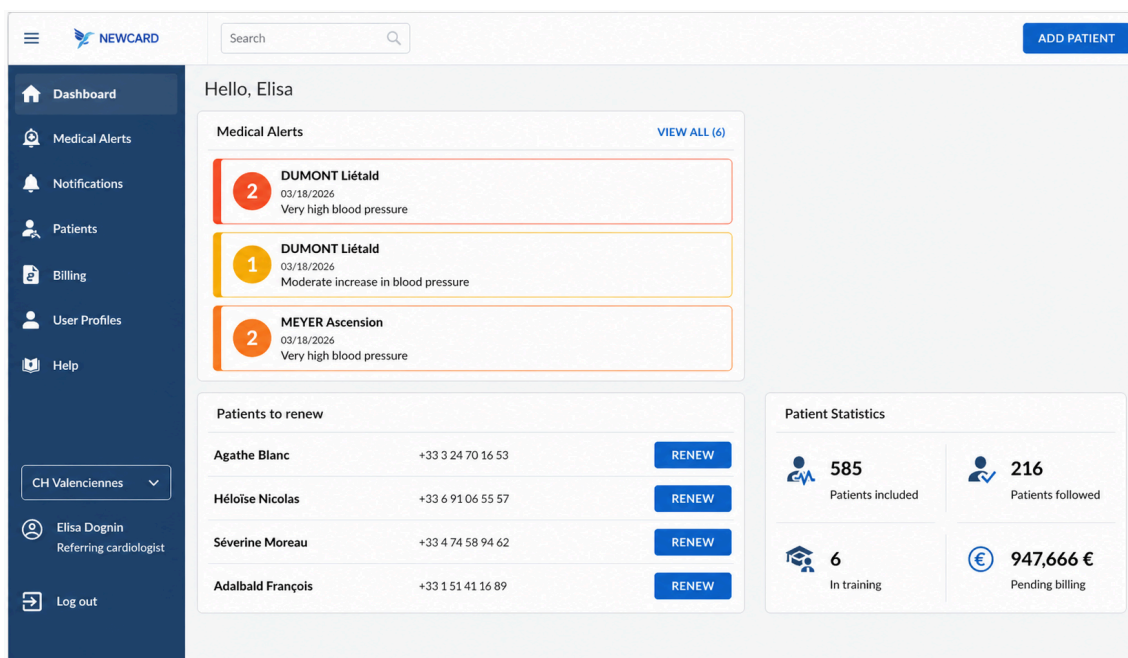The sidebar contains various categories of features.

Please note: Access to the tabs depends on the permissions assigned to each user (with the exception of the **Dashboard** and **Help**, which are accessible to everyone). As a result, the sidebar may vary from one user to another.

For example:

- Only users with the **'Profile and User Management'** permission can access the **'Profiles and Users'** tab. We recommend assigning this permission to the school's administrator.

- Only users with the **'Billing'** permission can access the **'Billing'** tab.
- etc.

→ For further details on managing accounts and user rights, please refer to section 11 of this guide.

Overview of the different tabs

When you click on one of the categories in the sidebar, you are taken to the corresponding tab, which gives you access to a specific feature of the platform.

The screenshot displays the NEWCARD Dashboard for user Elisa Dognin. The interface includes a sidebar with navigation options: Dashboard, Medical Alerts (highlighted), Notifications, Patients, Billing, User Profiles, and Help. The main content area is divided into several sections:

- Medical Alerts:** A list of alerts with a 'VIEW ALL (6)' link. The alerts are:
 - 2** DUMONT Liétald (03/18/2026): Very high blood pressure
 - 1** DUMONT Liétald (03/18/2026): Moderate increase in blood pressure
 - 2** MEYER Ascension (03/18/2026): Very high blood pressure
- Patients to renew:** A table listing patients with their phone numbers and 'RENEW' buttons.

Patient Name	Phone Number	Action
Agathe Blanc	+33 3 24 70 16 53	RENEW
Héloïse Nicolas	+33 6 91 06 55 57	RENEW
Séverine Moreau	+33 4 74 58 94 62	RENEW
Adalbald François	+33 1 51 41 16 89	RENEW
- Patient Statistics:** A summary of patient data.

Category	Count
Patients included	585
Patients followed	216
In training	6
Pending billing	€ 947,666
- Your demo patient:** A section for a demo patient, FÉLIX Mehdi, who is 'Active' and 'Followed'.

7. The **'Medical Alerts'** tab allows you to view all your current potential alerts and provides an overview of these alerts.

NEWCARD

Medical Alerts

Search

All 6

To take action 4

In progress 2

Terminated 1526

ALERT	PATIENT NAME	STATUS	LAST UPDATE	REVIEWED BY SPECIALIST	TRIGGERED BY	ACTIONS
2	DUMONT Liétald	Triggered	03/18/2026	01/20/2026	Very high blood pressure	TAKE ACTION
1	DUMONT Liétald	Triggered	03/18/2026	01/20/2026	Moderate increase in blood pressure	TAKE ACTION
2	MEYER Ascension	Triggered	03/18/2026	01/27/2025	Very high blood pressure	TAKE ACTION
1	MEYER Ascension	Triggered	03/18/2026	01/27/2025	Moderate increase in blood pressure	TAKE ACTION
1	PIERRE Joëlle	In progress	02/24/2026	02/09/2026	Very progressive weight gain	MANAGE IN PROGRESS
1	MICHEL Agathon	In progress	02/05/2026	02/09/2026	Clinical deterioration	MANAGE IN PROGRESS

Showing 1 to 6 of 6 entries

1

? → Please refer to section 9 for further information

- The **'Notifications'** tab allows you to interact with Newcard and receive alerts about any new information concerning your patients (lab results, comments, documents, etc.).

NEWCARD

Dashboard

Medical Alerts

Notifications

Patients

Billing

User Profiles

Help

Ella Duglin
Referring cardiologist

Log out

Notifications

Search

All 12

Comments 1

Lab Results 1

Documents 10

ALERT	PATIENT NAME	TYPE	NOTIFICATION DATE	ACTIONS
	BARBIER Geneviève	Comments	03/19/2026	VIEW
	BARBIER Geneviève	Lab Results	03/19/2026	VIEW
	SARON Antonin	Documents	03/19/2026	VIEW
	SARON Antonin	Documents	03/19/2026	VIEW
	SARON Antonin	Documents	03/19/2026	VIEW
	AUBRY Lethaïre	Documents	03/19/2026	VIEW
	AUBRY Lethaïre	Documents	03/19/2026	VIEW
	AUBRY Lethaïre	Documents	03/19/2026	VIEW

8. The **'Patients'** tab lists all the patients you are managing, whether they are awaiting training, currently being monitored, due for renewal or deactivated. Using the **search bar**, you can search for a patient directly by name.

Patients All medical conditions

[Add a patient](#)

All 58 Awaiting training 0 Follow-ups 23 To be renewed 6 Installation in progress 0 Disabled 35

ALERT	NAME, FIRST NAME ▲	PHONE	STATE	SPECIALIST	DATE D'INCLUSION ▲	ACTION(S)
	AUBRY Daphné Disabled	+33659520811	Follow up	ARNAUD Jean	20/05/2021	
	BARRE Fanny Active	+33659520811	Follow up	ARNAUD Jean	08/11/2025	RENEW
	BERGER Faustine Disabled	+33659520811	Follow up	ARNAUD Jean	19/11/2024	
	BERNARD Mathurin Active	+33659520811	Follow up	ARNAUD Jean	10/11/2025	RENEW
	BERNARD Véronique Disabled	+33659520811	Follow up	ARNAUD Jean	28/07/2022	
	BERTRAND Acanthe Disabled	+33659520811	Follow up	ARNAUD Jean	11/04/2023	
	BLANC Gislebert Active	+33659520811	Follow up	ARNAUD Jean	14/04/2026	
	BONNET Primerose Disabled	+33659520811	Follow up	ARNAUD Jean	11/06/2024	
	BOURGEOIS Camillien Active	+33659520811	Follow up	ARNAUD Jean	07/11/2025	RENEW

This button allows you to edit your **patient's key details**, namely their **registration details**, **personal details**, **contact information** and **medical records**.

You can also filter your patients by their medical condition using the drop-down menu at the top of your screen, next to the search bar.

Patients All medical conditions

[Add a patient](#)

All 58 Awaiting training 0 Follow-ups 23 Installation in progress 0 Disabled 35

ALERT	NAME, FIRST NAME ▲	PHONE	STATE	SPECIALIST	DATE D'INCLUSION ▲	ACTION(S)
	AUBRY Daphné Disabled	+33659520811	Follow up	ARNAUD Jean	20/05/2021	
	BARRE Fanny Active	+33659520811	Follow up	ARNAUD Jean	08/11/2025	RENEW

If you wish to view a patient's medical record, simply click on the relevant row.

9. The **'Billing'** tab allows you to track the months of remote monitoring that you can bill for; if you have set up automated billing, you can bill and submit claims directly to the National Insurance (sécurité sociale).

BILL	NAME-FIRST NAME	DATE OF BIRTH	NIR	ORGANIZATION CODE	PRESCRIBING PHYSICIAN
☐	BARRE FANNY	13/09/1967	267097425308730	01 444 5018	ARNAUD Jean
BILL	PRESCRIPTION DATE	INVOICE DATE	OBSERVANCE	STATUS	CODE - AMOUNT
☐	08/11/2025	05/02/2026	100%	Active	TVB 2 - 56€
☐	08/11/2025	07/03/2026	100%	Active	TVB 2 - 56€
☐	08/11/2025	06/04/2026	100%	Active	TVB 2 - 56€
☐	BERNARD MATHURIN	19/04/1990	190042652989019	/	ARNAUD Jean
BILL	PRESCRIPTION DATE	INVOICE DATE	OBSERVANCE	STATUS	CODE - AMOUNT
☐	10/11/2025	07/02/2026	100%	Active	TVB - 28€

→ Please refer to section 12 for further information

10. The **'Profiles and Users'** tab allows you to manage user accounts and permissions within your organisation (create, edit, deactivate).

It allows you to organise access rights according to your needs, including for non-medical users.

For a self-employed prescriber working alone, this tab is not usually necessary, except in cases where there is collaboration (secretarial support, billing, etc.).

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→ Please refer to section 11 for further information

11. The 'Help' tab gives you access to the Newcard help section, which is divided into two categories (cardiology and nephrology), where you will find the information you need to help you use the platform.

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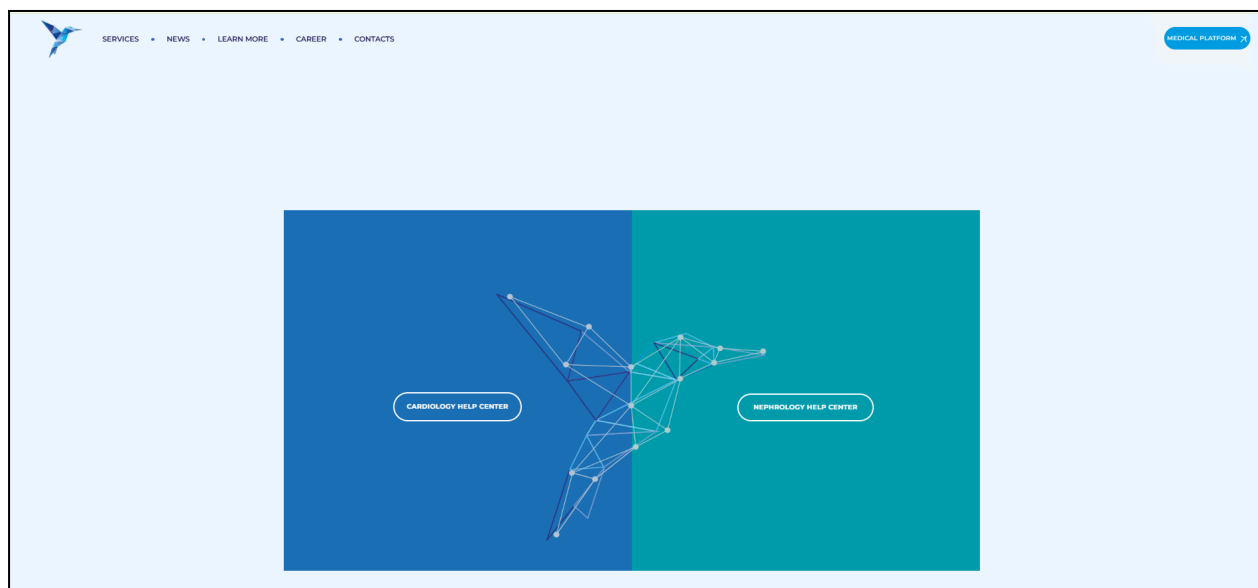
d. Updating personal information

To access your personal data, click on your name at the very bottom of the left-hand sidebar.

This will take you to a tab where you can edit your key details, such as:

- Your identity,
- Your password,
- Configuration details,
- Legal and regulatory documents such as the Terms of Service, legal notices, privacy policy and CE marking label,
- And the default units for reading biological data

The screenshot shows the 'Profile' page with the 'Identity' tab selected. The sidebar on the left contains navigation icons. The main content area has a header with a language dropdown (English) and a tab bar. The 'Identity' tab is active, displaying a form with two columns of input fields. The left column contains: 'Nom*' (Arnaud), 'Email principal*' (empty), 'N° RPPS*' (810000072222), 'Phone*' (with a French flag icon), 'Preferred language' (Anglais), 'Profession*' (Sage-Femme), and a 'REGISTER' button. The right column contains: 'First name*' (Jean), 'Secondary email' (empty), 'No CDM/Adeli/AM*' (empty), 'Date of birth' (empty), 'category*' (Civil), and 'Spécialité*' (empty).

Any changes made to your profile must be **saved** using the button at the bottom of each page.

6. Identity

This tab allows you to edit your general account details: **surname, first name, postal address and email address.**

This is a duplicate of the screenshot above, showing the 'Profile' page with the 'Identity' tab selected. The sidebar on the left contains navigation icons. The main content area has a header with a language dropdown (English) and a tab bar. The 'Identity' tab is active, displaying a form with two columns of input fields. The left column contains: 'Nom*' (Arnaud), 'Email principal*' (empty), 'N° RPPS*' (810000072222), 'Phone*' (with a French flag icon), 'Preferred language' (Anglais), 'Profession*' (Sage-Femme), and a 'REGISTER' button. The right column contains: 'First name*' (Jean), 'Secondary email' (empty), 'No CDM/Adeli/AM*' (empty), 'Date of birth' (empty), 'category*' (Civil), and 'Spécialité*' (empty).

7. Security

This section allows you to change your password. You will then need to save these changes by clicking on

REGISTER

The screenshot shows the 'Profile' page with the 'Security' tab selected. The page has a dark blue sidebar with icons for home, notifications, alerts, user profile, currency, and a folder. The main content area has tabs for Identity, Security (active), Configuration, Ordonnance de TLS, Documents, and Unité par défaut. Under the Security tab, there are three password input fields: 'Current password', 'New password', and 'Confirmation'. Each field has a toggle icon to show or hide the password. A blue 'REGISTER' button is located at the bottom left of the form.

8. Configuration


This section allows you to specify the types of alerts you wish to receive, as well as the channel through which you wish to receive them.

For clinics, it also allows you to filter the list of patients so that you can view only those for whom you are the prescribing doctor.

The screenshot shows the 'Profile' page with the 'Configuration' tab selected. The page has the same dark blue sidebar as the previous screenshot. The main content area has tabs for Identity, Security, Configuration (active), TLS Order, Documents, and Default unit. Under the Configuration tab, there are three sections: 'NOTIFICATIONS', 'FILTERS', and 'AUTHENTICATION'. The 'NOTIFICATIONS' section has three columns: 'Email principal', 'Secondary email', and 'SMS'. There are three rows of notifications: 'Medical alerts', 'Hospitalization of a patient', and 'View only patients for whom I am the prescriber'. Each row has checkboxes for each column. The 'Medical alerts' row has all three checkboxes checked. The 'Hospitalization of a patient' row has all three checkboxes unchecked. The 'View only patients for whom I am the prescriber' row has the 'Email principal' checkbox checked and the 'Secondary email' and 'SMS' checkboxes unchecked. The 'AUTHENTICATION' section has a label 'Double authentication validity duration' and a dropdown menu set to '15 days'.


9. Documents

This section contains Newcard's regulatory documents: the Terms and Conditions, the Legal Notice, the regulatory label for CE marking, and the Privacy Policy.

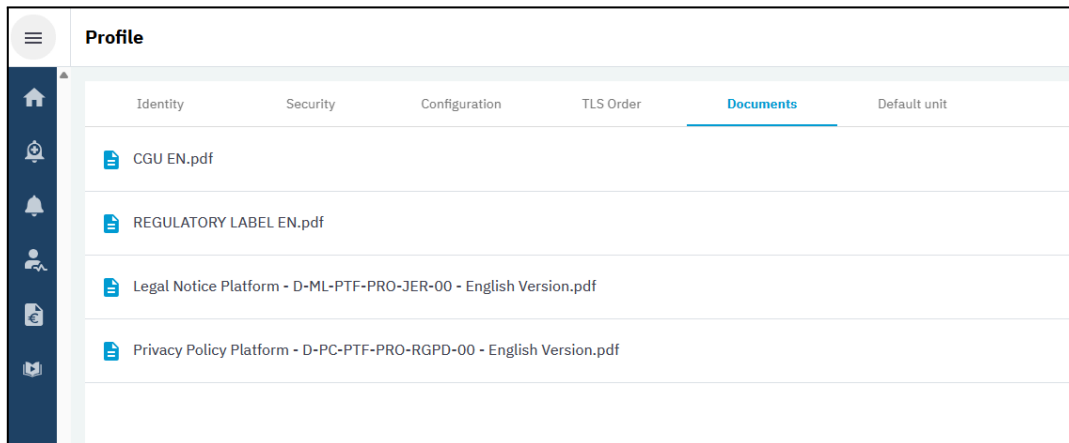10. Default units

When entering and viewing your patient's clinical data, you can select your preferred default units.

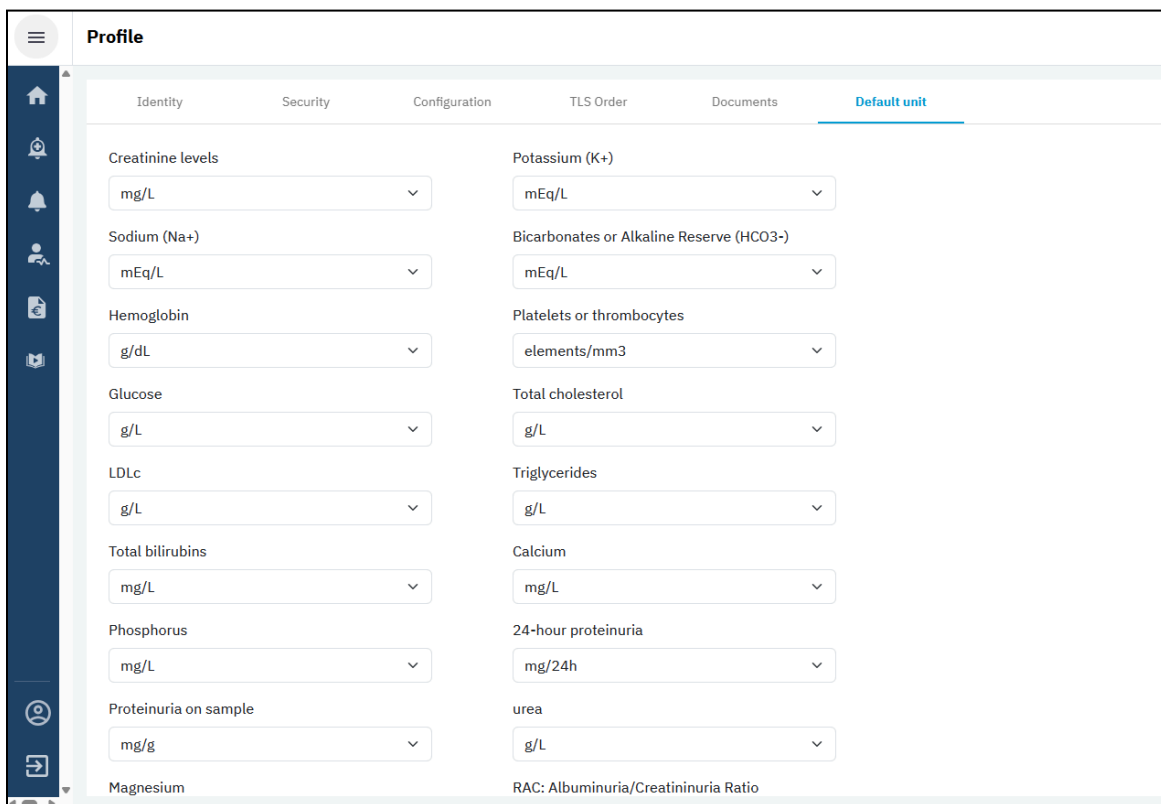
e. Access to the help section

The 'Help' tab gives you access to the Newcard help section, which is divided into two categories (cardiology and nephrology), where you will find the information you need to help you use the platform.

Each tab contains a Q & A, a video tutorial and a downloadable document to guide you through the platform's various features.

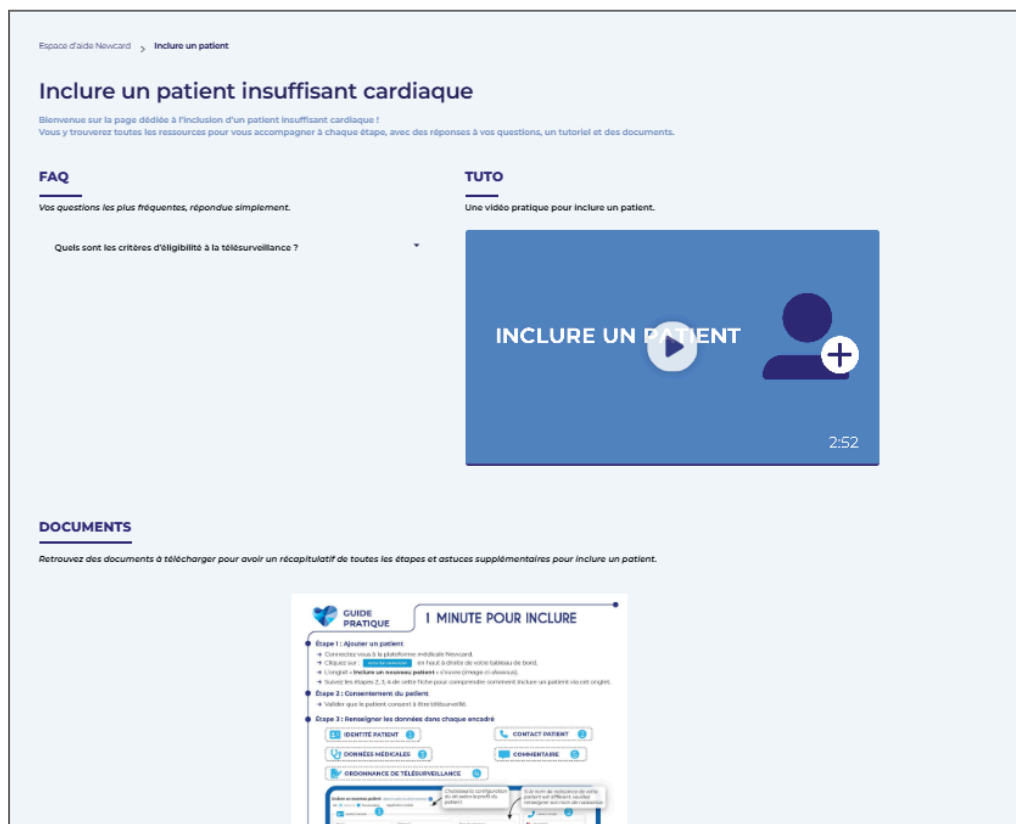5. Patient inclusion

The implementation of remote monitoring for a patient is subject to compliance with the eligibility criteria set out in the standard terms and conditions.

To add a patient to the platform, start by clicking on this button at the top right of your dashboard.

ADD A PATIENT

The dashboard shows a user profile for 'Bonjour, Arnaud Jean' with a search bar and an 'ADD A PATIENT' button. The main content area is divided into several sections:

- Medical alerts:** A list of three alerts for patients BARRE Fanny and BERNARD Mathurin, each with a date of 30/04/2026 and a description of the alert.
- Patients to be installed:** A section indicating that no patients are currently to be installed.
- Patients to be renewed:** A table listing patients Mathurin Bernard, Camillien Bourgeois, Florent David, and Georges Roy, each with a 'RENEW' button.
- Patient statistics:** A summary of patient statistics including 'Patients inclus' (58), 'Patients followed' (23), 'Waiting for training' (0), and 'Waiting for billing' (4690€).

Specifics of cardiology inclusion

A window will open on the right-hand side of your screen. Five fields must be completed with the required information to finalise the patient's inclusion.

The diagram illustrates the five fields required for patient inclusion, each represented by a dashed box with an icon and a numbered circle:

- PATIENT IDENTITY** (1): Represented by a person icon.
- PATIENT CONTACT** (2): Represented by a telephone icon.
- MEDICAL INFORMATION** (3): Represented by a stethoscope icon.
- REMOTE MONITORING ASSIGNMENT** (4): Represented by a document icon with a pulse line.
- COMMENTS** (5): Represented by a speech bubble icon.

Add a new patient in the context of the common law ⓘ

Kit ☒ Scale ☒ Tensiometer Mobile application Choose the kit configuration based on the patient profile.

PATIENT IDENTITY

Last name* First name* Birth name ⓘ

MALE FEMALE UNSPECIFIED Date of birth* Social security number (SSN)*

Identity verified by:* or ☐ Identity not verified Specify if you have verified the identity of your patient.

☐ Consent given by the patient*

CONTACT INFORMATION

Mobile phone* or Landline*

Email

Address, ZIP code, city

Spoken language English

MEDICAL INFORMATION

Current rhythm* Pacemaker* SOCLE package - TVB - \$28 / month ⓘ

Initial weight (kg)* Systolic function (%)

REMOTE MONITORING ASSIGNMENT

CH Valenciennes DOGNIN Elisa In a center, you can select the prescribing physician (name must be authorized for the device).

Mandatory prescription to invoice.
Drag and drop the prescription here or click to upload.

COMMENTS

Useful information for the patient's care pathway.
Example: hospitalized patient, details about the patient's surroundings, contact of the liberal nurse.

ADD PATIENT

! Please note: Do not forget to confirm the patient's consent to remote monitoring and to place the prescription in the designated space (required for billing purposes).

Once you have entered this information, click

ADD PATIENT

to complete the inclusion.

A confirmation message will appear.

You have two options:

- **'Complete inclusion':** the patient is immediately added to the Newcard platform and appears in the list of patients awaiting training. Once included, NEWCARD will contact them and provide them with the necessary equipment.
- **'Complete the information':** you are redirected to the patient's detailed profile

To help Newcard, complete the information.

Otherwise, Newcard will take over.

COMPLETE THE INFORMATION

COMPLETE INCLUSION

Complete details

Editing the patient

1. INCLUSION

2. IDENTITY

3. CONTACT

4. MEDICAL DATA

5. BILLING DATA

6. PATIENT-RELATED USERS

TYPE OF INCLUSION

☐ Common Law

- Hospitalization within the last 12 months for an episode of heart failure;
or
- Currently in NYHA class 2 or higher with elevated natriuretic peptide levels (BNP > 100 pg/mL or NT pro BNP > 1000 pg/mL)

☒ Consent given by the patient

REMOTE MONITORING PRESCRIPTION *

Prescriber*

ARNAUD Jean

Remote monitoring prescription to date 30/04/2026

Add your telemonitoring prescription template here.
Drag and drop a file or click.

FLAT-RATE FEE FOR BILLING THE PROCEDURE

TVB - 28€ / month

DOCUMENTS TO DOWNLOAD

 IC Inclusion Order

 IC Consent Form

 car

CONNECTED DEVICES

☒ Balance

☒ Blood pressure monitor

Application mobile

CANCEL

EDIT

If you choose to add further details, a new window will open, allowing you to provide more detailed information about your patient's admission. If you choose to add further details, you will see the same five tabs:

- Inclusion
- Identity
- Contact
- Medical data
- Patient-related users

Once you have finished, click **'Edit'** to save your changes. Your patient is now included and easily accessible in the list of all your patients.

In the full version, you can download the consent form to give to the patient and have them sign.

DOCUMENTS TO DOWNLOAD

IR Prescription Inclusion

Informed consent form...

Specificities of inclusion in nephrology

A window will open on your screen, containing 7 tabs. Don't worry: inclusion is quick and easy; simply follow the steps.

Tabs:

- Inclusion
- Identity
- Contact
- Medical data
- Alert settings
- Treatment
- Patient-related users

Tab 1 – Inclusion

Include a new patient

1. INCLUSION

2. IDENTITY

3. CONTACT

4. MEDICAL DATA

5. MEDICAL DATA II

6. TREATMENTS

7. PATIENT-RELATED USERS

TYPE OF INCLUSION

☒ STAGE IV
Inclusion under common law:
Patients with stage 4 chronic kidney disease
Glomerular filtration rate between 15 and 30 ml/min/1.73m2
[Exclusion criteria](#)

☐ STAGE V
Inclusion under common law:
Patients with stage 5 chronic kidney disease
Glomerular filtration rate less than 15 ml/min/1.73m2

☐ GRAFT
Inclusion under common law:
Kidney transplant patients who received their transplants more than 3 months ago

☐ Patient consent*

REMOTE MONITORING ORDER *

Prescriber of remote patient monitoring*

Arnaud Adrehilde

Remote monitoring order to be dated 30/04/2026

Prescription required for invoicing.
Drag and drop the prescription here or click.

DOCUMENTS TO DOWNLOAD

IR Inclusion Order

IR Consent Form

KIT

Mobile application

Comment

NEXT STEP >

- Select the inclusion criteria: patients with stage 4 or 5 kidney failure, or transplant recipients.
- Confirm the patient's consent to remote monitoring.
- Upload the prescription (required for billing).
- Specify the prescriber if you are in a healthcare facility.
- Choose the kit configuration: mobile app or tablet.
- Add a comment if necessary.

Once you have entered this information, click on

NEXT STEP >

Please note: Don't forget to obtain consent and submit the prescription.

Tab 2 – Identity

Include a new patient

- 1. INCLUSION
- 2. IDENTITY
- 3. CONTACT
- 4. MEDICAL DATA
- 5. MEDICAL DATA II
- 6. TREATMENTS
- 7. PATIENT-RELATED USERS

IDENTIFY

RETRIEVE THE INS BY IDENTITY TRAITS

RETRIEVE YOUR INS NUMBER USING YOUR HEALTH INSURANCE CARD.

IDENTITY DATA

INS registration number Temporary Identity IPP Registration Number

Sex* MAN WOMEN UNDETERMINED

Name used* Dupont First name used* Jeans

Birth name* Dupont First name(s) at birth* Charles Maurice first name at birth Charles

Date of birth* Place of birth Lille

Social security number (NIR)* Code of the affiliated organization

☐ Not affiliated with French social security
In this case, contact Newcard

Identity verified using* OR ☐ Identity not verified OR ☐ Questionable identity i

< PREVIOUS STEP NEXT STEP >

- If necessary, retrieve the identity details via the INS online service.
- Fill in: gender, surname, first name, maiden name, date of birth, National Insurance number.

- Indicate whether the identity has been verified and by what document (ID card, passport, residence permit, or not verified).

Click on

NEXT STEP >

Tab 3 – Patient contact details

- Mobile and landline numbers (if available)
- Email address and full postal address
- Contact person in case of absence (with link)

Include a new patient

- 1. INCLUSION
- 2. IDENTITY
- 3. CONTACT
- 4. MEDICAL DATA
- 5. MEDICAL DATA II
- 6. TREATMENTS
- 7. PATIENT-RELATED USERS

IDENTIFY

Portable*

Fixed*

Other phone

OR

Address*

Mail

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Spoken language

French

PERSON TO CONTACT IN CASE OF ABSENCE

Name

Phone

Mail

Nature of the link

< PREVIOUS STEP

NEXT STEP >

Click on

NEXT STEP >

Tabs 4 & 5 – Medical data

- **Tab 4:** GFR, height, weight (mandatory), serum creatinine, proteinuria or RAC

Include a new patient

1. INCLUSION

2. IDENTITY

3. CONTACT

4. MEDICAL DATA

5. MEDICAL DATA II

6. TREATMENTS

7. PATIENT-RELATED USERS

DFG(CKD-EPI)

ml/min/1.73 m²

Or

DFG(MDRD)

ml/min/1.73 m²

Size* (cm)

Initial weight (kg)*

PRIMARY DIAGNOSIS OF CAUSATIVE NEPHROPATHY

☐ Diabetic nephropathy

☐ Polycystic kidney disease

☐ Primary glomerulopathy

☐ Others...

☐ Vascular nephropathy

☐ Chronic tubulointerstitial nephropathy

BASELINE LEVELS AT INCLUSION

Creatinine*

mg/L

Hematuria

YES

UNKNOWN

NO

Proteinuria over 24 hours*

mg/24h

Or

Proteinuria on sample*

mg/g

RAC: Albuminuria / Creatininuria ratio*

mg/mmol

< PREVIOUS STEP

NEXT STEP >

- **Tab 5:** Blood pressure and serum potassium levels, any comorbidities

Include a new patient

1. INCLUSION

2. IDENTITY

3. CONTACT

4. MEDICAL DATA

5. MEDICAL DATA II

6. TREATMENTS

7. PATIENT-RELATED USERS

BLOOD PRESSURE RANGES IN MMHG.*

Systolic

180

170

160

150

140

130

120

OTHER

Diastolic

110

100

90

80

70

OTHER

The above limits allow the triggering of a level 1 alert following three consecutive exceedances of systolic and/or diastolic pressure.

POTASSIUM LEVELS IN MEQ/L*

Hypokalemia

2.0

2.5

3.0

3.5

Hyperkalemia

5.0

5.5

6.0

6.5

COMORBIDITY(IES)

Diabetes

YES

UNKNOWN

NO

Cardiovascular disease

YES

UNKNOWN

NO

HTA

YES

UNKNOWN

NO

Hepatic steatosis associated with a metabolic disorder

YES

UNKNOWN

NO

Non-alcoholic fatty liver disease

YES

UNKNOWN

NO

Obesity

YES

UNKNOWN

NO

Sarcopenia

YES

UNKNOWN

NO

Cancer

YES

UNKNOWN

NO

Nephrotoxic treatment

YES

UNKNOWN

NO

< PREVIOUS STEP

NEXT STEP >

Click on

NEXT STEP >

after each tab.

Tab 6 – Treatment

- Enter the patient's treatment details if necessary.

Include a new patient

1. INCLUSION
2. IDENTITY
3. CONTACT
4. MEDICAL DATA
5. MEDICAL DATA II
6. TREATMENTS
7. PATIENT-RELATED USERS

Statins

YES UNKNOWN NO

Antidiabetics

YES UNKNOWN NO

Antiplatelet agents...

YES UNKNOWN NO

Anticoagulants...

YES UNKNOWN NO

Additional Calcium +/- Vitamin D...

YES UNKNOWN NO

IPP

YES UNKNOWN NO

MIMPARA

YES UNKNOWN NO

Iron

YES UNKNOWN NO

Potassium chelating resin...

YES UNKNOWN NO

EPO

YES UNKNOWN NO

Beta-blocker

YES UNKNOWN NO

Another bradycardic

YES UNKNOWN NO

IEC / sartans

YES UNKNOWN NO

Diuretics

YES UNKNOWN NO

Anti-aldosterone

YES UNKNOWN NO

Sacubitril

YES UNKNOWN NO

Gliflozin

YES UNKNOWN NO

Calcium channel blockers

YES UNKNOWN NO

Other antihypertensives

YES UNKNOWN NO

< PREVIOUS STEP

NEXT STEP >

Tab 7 – Patient-related users

- The centre and the prescriber are already entered.
- Add other professionals if necessary: registered nurse, another nephrologist, administrative staff, etc.

Include a new patient

1. INCLUSION
2. IDENTITY
3. CONTACT
4. MEDICAL DATA
5. MEDICAL DATA II
6. TREATMENTS
7. PATIENT-RELATED USERS

USERS LINKED TO THE PATIENT

Healthcare facility*

Roux et Baron

Prescriber of remote patient monitoring*

Arnaud Adrehilde

Healthcare staff or administrative staff

Search for a user

< PREVIOUS STEP

INCLUDE

Once you have finished, click on

INCLUDE

Your patient is now included and can be found in the list of all your patients.

Patient identity verification (offered at inclusion, not mandatory)

What is the INS?

The National Health Identity (INS) is a unique, permanent identifier shared by all healthcare professionals, and has been mandatory since 1 January 2021.

Why verify the patient's identity?

- To minimise patient identification errors
- To facilitate the exchange and sharing of health data and medical information
- To contribute to the quality of care and patient safety
- To avoid duplicate records

How can you check it?

1. You can verify your patient's identity at the time of inclusion.
2. In the section for editing patient details, you can retrieve the patient's identity via the INS online service. This allows you to obtain the patient's INS registration number and ensure the security of their information.

Every patient has a unique INS registration number.

The 5 INS identity details (which you will find in the inclusion section)

1. Birth name
2. Date of birth
3. Gender
4. List of birth first names
5. Place of birth code

**Importance of Patient Identity Verification**

To ensure patient safety and prevent errors:

→ **Verify the patient's identity at least once using an official identity document**
(National Identity Card, passport, permanent residence permit).

During inclusion

During the edition of patient details

→ **A qualified identity (verified + retrieved)** helps secure and ensure the reliability of medical and administrative data exchanges with third-party services.

6. Patient renewal

Please note: patient renewal requires authorisation (refer to chapter 11: Profile and User Management). If you do not have the necessary permissions, the renewal function will not be available.

To renew a patient's subscription on the platform, you have two options:

- **From the dashboard**

Once on the dashboard, the list of patients due for renewal is displayed immediately. Simply click on the blue button

RENEW

The screenshot shows a dashboard for a user named Arnaud Jean. It features a sidebar with navigation icons and a main content area. The main content area is divided into several sections:

- Medical alerts:** A list of three alerts for patients BARRE Fanny and BERNARD Mathurin, each with a date (30/04/2026) and a description of the alert (e.g., 'Very high voltage', 'Moderate increase in blood pressure', 'Rapid weight gain').
- Patients to be installed:** A section indicating that there are no patients to install.
- Patients to be renewed:** A table listing four patients: Mathurin Bernard, Camillien Bourgeois, Florent David, and Georges Roy. Each patient has a blue 'RENEW' button next to their name.
- Patient statistics:** A section showing four statistics: 59 Patients included, 23 Patients followed, 1 Waiting for training, and €4690 Waiting for billing.

- **From the 'Patients' tab**

You can also access this list via the **'Patients'** tab (in the left-hand sidebar), then under the **'Due for renewal'** section. You will find all the relevant patients there; simply click on

RENEW

For your information: Patients are listed in order of date, with the first on the list being the most urgent to renew.

Patients							Add a patient	
All	Waiting for training	Follow-up	To be renewed	In the process of being installed	Disabled			
59	1	23	6	0	35			
LAST NAME, FIRST NAME	PHONE	STATUS	SPECIALIST	DATE	DATE OF INCLUSION	ACTION(S)		
BARRE Fanny Active	+33659520811	Monitored	ARNAUD Jean	07/05/2026	08/11/2025	RENEW		
BERNARD Mathurin Active	+33659520811	Monitored	ARNAUD Jean	09/05/2026	10/11/2025	RENEW		
BOURGEOIS Camillien Active	+33659520811	Monitored	ARNAUD Jean	06/05/2026	07/11/2025	RENEW		
DAVID Florent Active	+33659520811	Absence	ARNAUD Jean	27/04/2026	29/10/2025	RENEW		
RENAULT Rachel Active	+33659520811	Monitored	ARNAUD Jean	28/04/2026	30/10/2025	RENEW		
ROY Georges Active	+33659520811	Monitored	ARNAUD Jean	20/05/2026	21/11/2025	RENEW		
Displaying the element 1 to 6 on 6							1	

Specifics regarding renewal – heart failure

When renewing, a page will open on the right-hand side of your screen.

Patients		Simplified renewal within the framework of common law																																				
All 59 Waiting for training 1 Follow-up 23 To be renewed 6 In the process of being installed 0 Disabled 35 <table border="1"> <thead> <tr> <th>LAST NAME, FIRST NAME</th> <th>PHONE</th> <th>STATUS</th> <th>SPECIALIST</th> <th>DATE</th> </tr> </thead> <tbody> <tr> <td>BARRE Fanny Active</td> <td>+33659520811</td> <td>Monitored</td> <td>ARNAUD Jean</td> <td>07/05/2026</td> </tr> <tr> <td>BERNARD Mathurin Active</td> <td>+33659520811</td> <td>Monitored</td> <td>ARNAUD Jean</td> <td>09/05/2026</td> </tr> <tr> <td>BOURGEOIS Camillien Active</td> <td>+33659520811</td> <td>Monitored</td> <td>ARNAUD Jean</td> <td>06/05/2026</td> </tr> <tr> <td>DAVID Florent Active</td> <td>+33659520811</td> <td>Absence</td> <td>ARNAUD Jean</td> <td>27/04/2026</td> </tr> <tr> <td>RENAULT Rachel Active</td> <td>+33659520811</td> <td>Monitored</td> <td>ARNAUD Jean</td> <td>28/04/2026</td> </tr> <tr> <td>ROY Georges Active</td> <td>+33659520811</td> <td>Monitored</td> <td>ARNAUD Jean</td> <td>20/05/2026</td> </tr> </tbody> </table>		LAST NAME, FIRST NAME	PHONE	STATUS	SPECIALIST	DATE	BARRE Fanny Active	+33659520811	Monitored	ARNAUD Jean	07/05/2026	BERNARD Mathurin Active	+33659520811	Monitored	ARNAUD Jean	09/05/2026	BOURGEOIS Camillien Active	+33659520811	Monitored	ARNAUD Jean	06/05/2026	DAVID Florent Active	+33659520811	Absence	ARNAUD Jean	27/04/2026	RENAULT Rachel Active	+33659520811	Monitored	ARNAUD Jean	28/04/2026	ROY Georges Active	+33659520811	Monitored	ARNAUD Jean	20/05/2026	SINUS FA SV DISORDER Flat rate coefficient TVB 2 - €56/month Initial weight (kg)* 50,71 Date of birth 13/09/1967 Social security number 267097425308730 Health Facility* Gaillard SEM Referring physician* ARNAUD Jean ☒ Consent given by the patient Remote monitoring prescription to date 05/07/2026 Prescription required to bill. Drag and drop the prescription here or click. How COMPLETE AND RENEW RENEW	
LAST NAME, FIRST NAME	PHONE	STATUS	SPECIALIST	DATE																																		
BARRE Fanny Active	+33659520811	Monitored	ARNAUD Jean	07/05/2026																																		
BERNARD Mathurin Active	+33659520811	Monitored	ARNAUD Jean	09/05/2026																																		
BOURGEOIS Camillien Active	+33659520811	Monitored	ARNAUD Jean	06/05/2026																																		
DAVID Florent Active	+33659520811	Absence	ARNAUD Jean	27/04/2026																																		
RENAULT Rachel Active	+33659520811	Monitored	ARNAUD Jean	28/04/2026																																		
ROY Georges Active	+33659520811	Monitored	ARNAUD Jean	20/05/2026																																		

1. Checking the information

- The fields are pre-filled based on the previous entry
- Check that the information is up to date
- Amend it if necessary

2. Remote monitoring package

⚠ Upon renewal, the package automatically reverts to the **basic package (€28)**.

➔ Remember to review it in light of your patient's circumstances

3. Prescription

- Submit a new prescription for this term
- The prescription date must match the **date specified** (end of the term + 1 day)

4. Finalising the renewal

There are two available options:

- **'Renew'** allows you to confirm the renewal immediately using the current details
- **'Complete and renew'** gives you access to the medical data from previous entries, which you can update before confirming

The screenshot displays the 'Patient renewal' form. On the left, a sidebar lists four tabs: 1. INCLUSION, 2. MEDICAL DATA, 3. PATIENT-RELATED USERS, and 4. TREATMENTS. The main form area is titled 'Patient renewal' and contains several sections:

- NYHA / BNP**: Includes a dropdown for NYHA (1, 2, 3, 4) and a checkbox for 'Increased vigilance'.
- BNP (pg/mL)**: A text input field with the value '3500'.
- MEDICAL DATA**: Includes fields for 'Size* (cm)' (162), 'Initial weight (kg)*' (50,71), and 'Systolic function (%)' (60).
- Diabetes**: A row of checkboxes for YES, UNKNOWN, and NO.
- HTA**: A row of checkboxes for YES, UNKNOWN, and NO.
- Respiratory disease**: A row of checkboxes for YES, UNKNOWN, and NO.
- Chronic anemia**: A row of checkboxes for YES, UNKNOWN, and NO.
- Betablocker**: A row of checkboxes for YES, UNKNOWN, and NO.
- Other bradycardia**: A row of checkboxes for YES, UNKNOWN, and NO.
- IEC / sartans**: A row of checkboxes for YES, UNKNOWN, and NO.
- Diuretics**: A row of checkboxes for YES, UNKNOWN, and NO.
- Anti-aldosterone**: A row of checkboxes for YES, UNKNOWN, and NO.
- Sacubitril**: A row of checkboxes for YES, UNKNOWN, and NO.
- Gliflozin**: A row of checkboxes for YES, UNKNOWN, and NO.
- History of supra rhythm disorders v**: A row of checkboxes for YES, UNKNOWN, and NO.
- Current pace**: A row of checkboxes for SINUS, FA, and SV DISORDER.
- Renal insufficiency (clearance in mL/min)**: A row of checkboxes for NO, < 30, 30-60, and > 60.
- Stimulator**: A row of checkboxes for NO, DEFIBRILLATOR, PACEMAKER, and LIFEVEST.

At the bottom of the form, there are two buttons: '< PREVIOUS STEP' and 'NEXT STEP >'.

Once this has been done, your patient will be re-enrolled for a new term.

Specifics of renewal – kidney failure/transplant recipients

When renewing, a page will open on your screen.

It consists of seven tabs

- **Inclusion,**
- **Algorithmic data,**
- **Medical data,**
- **Treatments,**
- and **patient-related users.**

Patient renewal

- 1. INCLUSION
- 2. MEDICAL DATA
- 3. MEDICAL DATA II
- 4. TREATMENTS
- 5. PATIENT-RELATED USERS

TYPE OF INCLUSION

☒ STAGE IV
 Inclusion under common law:
 Patients with stage 4 chronic kidney disease
Glomerular filtration rate between 15 and 30 ml/min/1.73m2
[Exclusion criteria](#)

☐ STAGE V
 Inclusion under common law:
 Patients with stage 5 chronic kidney disease
Glomerular filtration rate less than 15 ml/min/1.73m2

☐ GRAFT
 Inclusion under common law:
 Kidney transplant patients who received their transplants more than 3 months ago

☒ Patient consent*

REMOTE MONITORING ORDER *

Prescriber of remote patient monitoring*

☒ I authorize the use of my remote monitoring order template

The order will be dated **14/05/2026** , contact Newcard if this date is not convenient for you.

DOCUMENTS TO DOWNLOAD

☒ IR Inclusion Order

☒ IR Consent Form

Comment

NEXT STEP >

Please take the time to check that the information entered in each tab is complete and up to date. The pre-filled fields, taken from the previous entry, can be edited if necessary.

Don't forget to submit the prescription for the new term. It must be dated the day after the end of the current term (end date of the term + 1 day).

Once everything has been checked and completed, click on

RE-INCLUDE


Once this has been done, your patient will be included again for the next six months.

7. Patient list

The Patients list makes it easy to navigate between patient records. It allows you to quickly find a specific patient, or view a list of patients due for a follow-up appointment, for example.

The **'Patients'** tab is visible if your profile permissions include **'View heart failure patients'** and/or **'View kidney failure patients'**.

The list of patients under your care is displayed by default:

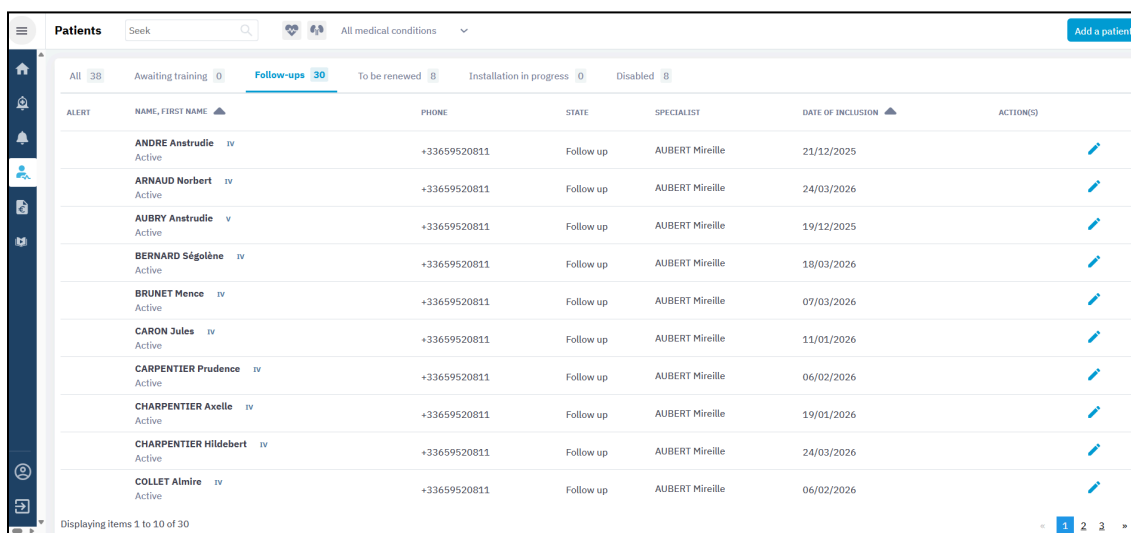Patients						
Seek						
All medical conditions						
Add a patient						
All 38	Awaiting training 0	Follow-ups 30	To be renewed 8	Installation in progress 0	Disabled 8	
ALERT	NAME, FIRST NAME	PHONE	STATE	SPECIALIST	DATE OF INCLUSION	ACTIONS
Active	ANDRE Anstrudie IV	+33659520811	Follow up	AUBERT Mireille	21/12/2025	
Active	ARNAUD Norbert IV	+33659520811	Follow up	AUBERT Mireille	24/03/2026	
Active	AUBRY Anstrudie V	+33659520811	Follow up	AUBERT Mireille	19/12/2025	
Active	BERNARD Ségolène IV	+33659520811	Follow up	AUBERT Mireille	18/03/2026	
Active	BRUNET Mence IV	+33659520811	Follow up	AUBERT Mireille	07/03/2026	
Active	CARON Jules IV	+33659520811	Follow up	AUBERT Mireille	11/01/2026	
Active	CARPENTIER Prudence IV	+33659520811	Follow up	AUBERT Mireille	06/02/2026	
Active	CHARPENTIER Axelle IV	+33659520811	Follow up	AUBERT Mireille	19/01/2026	
Active	CHARPENTIER Hildebert IV	+33659520811	Follow up	AUBERT Mireille	24/03/2026	
Active	COLLET Almire IV	+33659520811	Follow up	AUBERT Mireille	06/02/2026	
Displaying items 1 to 10 of 30						

If your organisation treats patients with both conditions and both permissions are activated, you can filter for patients with heart failure using the icon  or for patients with kidney failure using the icon .

Other filters are also available, such as: search by patient's surname or first name, filter by medical condition, or filter by stage for patients with kidney failure (stage IV, V or transplant recipients).

Almost all lists are sorted alphabetically by default, but it is entirely possible to sort them by date of inclusion by clicking on the arrow to the right of the inclusion date

DATE OF INCLUSION 

The **'All'** patient list includes patients in all different statuses (Setup in progress, Awaiting training, Under review, Deactivated).


The **'Awaiting training'** patient list includes all patients who have just been added and have not yet taken any measurements.

The **'Monitored' / 'Follow-ups'** patient list includes all patients who have taken measurements.

The **'To be renewed'** patient list includes all patients requiring renewal: the prescription's expiry date is within 30 days or less, or the prescription's expiry date has already passed. Here, the list is sorted chronologically: the first entry corresponds to the renewal due earliest. A "Renew" button is directly available in this list.

The **'Installation in progress'** patient list includes all patients whose monitoring was previously stopped and who are currently being re-equipped.

The **'Desactivated' / 'Disabled'** patient list includes all patients for whom monitoring has been stopped (exclusion from criteria, death, request to stop from the patient, etc.)

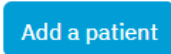If you are authorised to add patients, the button  is also available on this tab.

8. Patient form

The patient form is the main tool for monitoring your patients remotely. On this form, you will find the data regularly transmitted by the patient, their alert history, periods of hospitalisation or absence, and the period of remote monitoring. It is also on the patient record that you can indicate that an alert has been dealt with, and note any changes to the patient's medical data or treatment.

The patient form varies depending on the patient's condition, particularly in terms of data access. Furthermore, the actions available depend on the permissions associated with the user's profile (please refer to 11. Profile and User Management).

a. Administrative information

You will find all the essential information for patient care in the column on the right-hand side of your screen.

Follow-up period 

ORDER IN PROGRESS
21/12/2025

EXPIRATION DATE OF THE ORDER
18/06/2026

PRESCRIPTION
-

TERMS AND CONDITIONS ACCEPTED
02/12/2025

Documents

BIOLOGICAL ASSESSMENTS 

REPORTS 

TREATMENTS 

ANALYSIS PHOTOS

OTHERS 

Identify 

STATUS
Follow Up Active **IV**

SEX *
Women

You will find:

- The follow-up period, including the start date and the end date,
- A copy of the current remote monitoring prescription,
- The date on which the patient accepted the Terms and Conditions,
- Any documents you have attached to your patient's file,
- The patient summary with their contact details and medical condition,
- The healthcare team with a summary of the healthcare professionals involved in their remote monitoring,
- The history of remote monitoring prescriptions.

This button  allows you to edit your patient's key details, such as :

- Their enrolment details,
- Their personal details,
- Their contact details.

DATE OF BIRTH * 17/03/1973 - 53 years old REF. INST. AA-45f6 NEW INST. REF. 22e-45f6 NIR 273033049490922 CODE OF THE AFFILIATED ORGANIZATION - AVERAGE TIME BETWEEN TWO DOSES 6.8 days Contact PORTABLE +33659520811  FIXED - E-MAIL Anstrudie_andre36@Patient.Ncd ADDRESS 268 Place Dauphine 16217 Bourges	Person to contact in case of absence NAME Xavière Laurent E-MAIL Xavire_laurent@Patient.Ncd PHONE +33659520811 NATURE OF THE LINK Spouse Management  HEALTH FACILITY REMY ET HENRY  TLS PRESCRIBER AUBERT Mireille  INSTALLER INSTALLATEUR Officiel  OBSERVER JACQUET Désiré  THERAPEUTIC COMPANION LEMOINE Amande TLS Orders Prescription 24/06/2025 Prescription 26/12/2024
--	--

b. Active alert

You can find details of the alerts at the top of the relevant patient's record. You can indicate on the platform the actions you have taken to resolve the alert.

2

Alert 2

Very high voltage

Moderate increase in blood pressure

30/04/2026

TAKING CHARGE


→ To do this, please refer to Section 9: Alert Management. You will then find these details in your patient's medical history.

c. Access to remote telemonitoring data


If you are not a healthcare professional:

- You do not have access to the various graphs containing medical data, but you do have access to this graph, which shows whether or not the patient has taken their measurements:

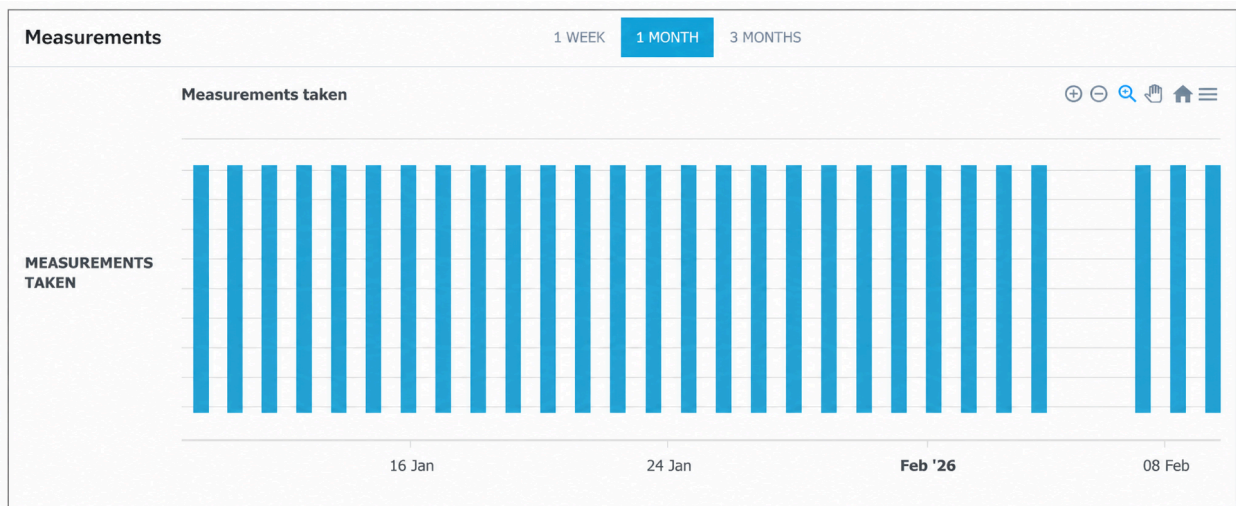
If you are a healthcare professional:

- Patients with heart failure have a weight graph, a blood pressure graph, a heart rate graph and a symptoms graph.
- Patients with kidney failure (stage IV, stage V and transplant recipients) have a serum creatinine graph, a blood pressure graph, a weight graph and their answers to the two questions asked of the user.

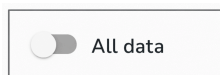Using this feature, you can see the time the measurement was taken, as well as any measurements that have been excluded because they were deemed outliers.

On the graphs, for each reading, you can see whether it was taken manually by the patient or entered manually by Newcard. If nothing is indicated, it means the reading was taken using a connected device (blood pressure monitor or scales).

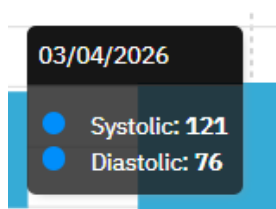For patients with heart failure

On the patient record, you can view the following graphs simultaneously:

- Weight,
- Blood pressure,
- Heart rate
- Symptom tracking:
 - Each time measurements are taken, the patient can select
 - “I feel the same as usual”
 - or “I feel worse than yesterday”.

If the patient feels worse than yesterday, a list of symptoms appears:

- I’m coughing more
- I’m more short of breath
- I’m using an extra pillow when I sleep
- I’m more tired
- My ankles are more swollen
- Nothing to report

If the patient feels as usual, the assessment is complete.

So, if the indicator is **blue**, this means that the patient has replied ‘**I feel as usual**’ or ‘**Nothing to report**’, or has not selected that symptom from the list. If it is **red**, the patient has selected that particular symptom.

You can adjust the data display period by holding down the left mouse button on one of the charts. Then, drag the cursor to change the period.

All the charts are linked, so the period you select will be applied to each of them.


You can also select a predefined period at the top of the patient record, above the first chart.

For kidney transplant patients

On the patient record, you can view the following graphs simultaneously:

- Creatinine levels,
- Immunosuppressant levels: the immunosuppressant level graphs correspond to the immunosuppressant treatment you specified earlier.

- Blood pressure,
- Weight
- Heart rate,
- Immunosuppressant treatment adherence monitoring: each time measurements are taken, the patient answers the questions: “Have you forgotten to take your immunosuppressant treatment in the last week?”, “Have you delayed taking your immunosuppressant treatment by more than an hour in the last week?”

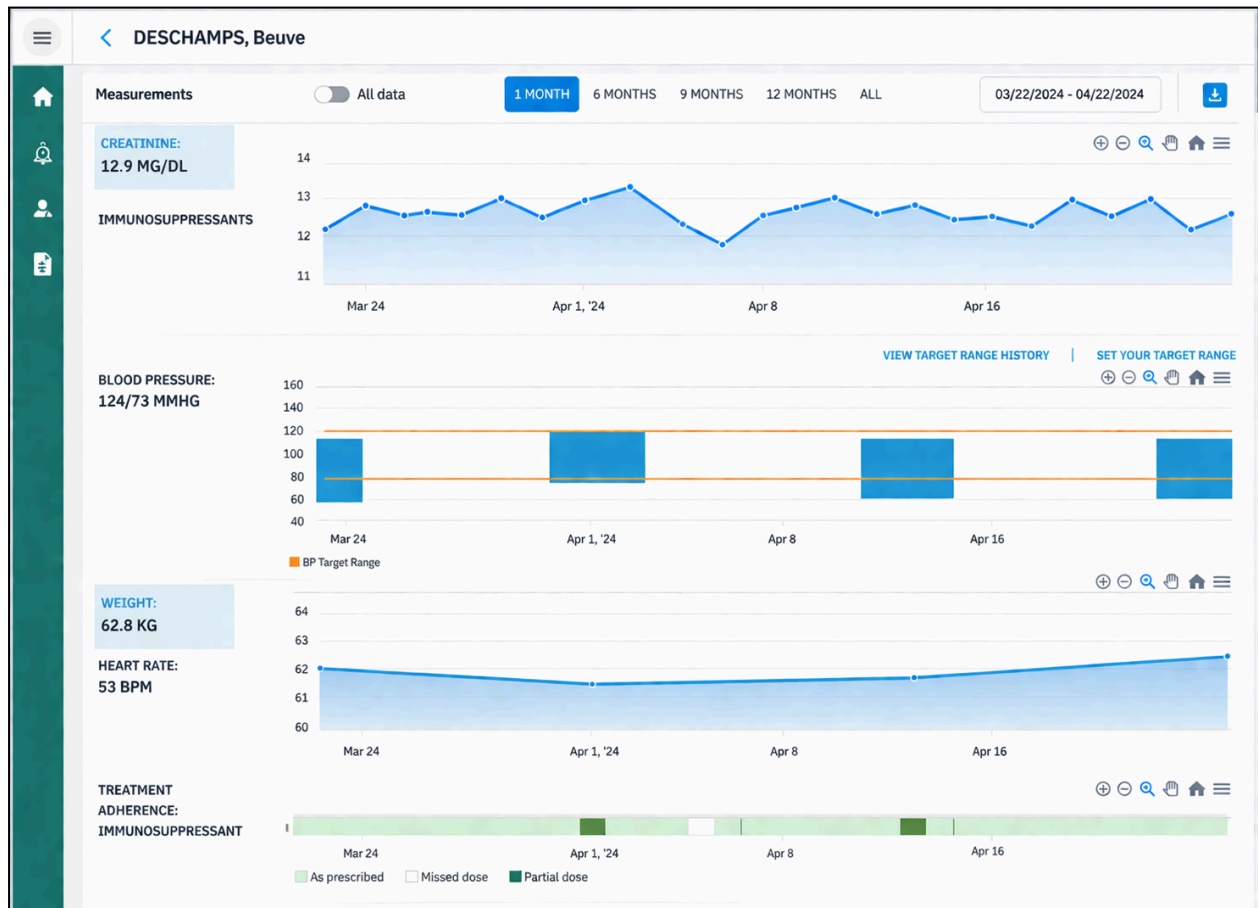If the indicator is light green, ■ this means that the patient answered ‘No’ to both questions. If it is dark green, ■ the patient answered ‘Yes’ to one of the two questions.

If it is a darker shade of green, ■ the patient has answered ‘Yes’ to both questions.

For the weight, blood pressure and heart rate graphs, you can view the date/time and measurement value in the tooltip

For patients in stages IV and V

On the patient record, you can view the following graphs simultaneously:

- Serum creatinine,
- Blood pressure,
- Weight,
- Heart rate,
- Health status compliance:
- Each time measurements are taken,
- The patient answers the following questions:

→ “How are you breathing today?”

→ and “Have you noticed any deterioration in your general condition since your last dose?”

If the indicator is blue, ■ this means that the patient answered ‘As usual’ to the first question and ‘No’ to the second question.

If it is red, ■ the patient answered ‘More short of breath than usual’ to the first question and ‘Yes’ to the second question.

For the weight, blood pressure and heart rate graphs, you can view the date/time and measurement value in a tooltip.

You can adjust the data display period by holding down the left mouse button on one of the graphs. Then, drag the cursor to change the period.

All the graphs are linked, so the period you select will change on each of them.

You can also select a predefined period at the top of the patient record, above the first graph.

d. Visualisation of medical and biological data

Display of medical data

If the data has not been entered, it will not appear in this section.

Medical data 				
02/17/2026				
ORDER DATE	NYHA CLASS	NT-proBNP	HEIGHT	SYSTOLIC FUNCTION
02/17/2026	2	3460 pg/mL	186 cm	63 %
HYPERTENSION (HTN)	DIABETES	CHRONIC ANEMIA	RESPIRATORY FAILURE	RENAL INSUFFICIENCY
No	No	No	No	eGFR 30 – 60
HISTORY OF VENTRICULAR ARRHYTHMIAS (SUPRAV)	CURRENT RHYTHM	DEVICE	BETA-BLOCKER	OTHER BRADYCARDIC DRUG
No	Sinus rhythm	Pacemaker	Yes	No
ACEI / ARB	DIURETICS	ALDOSTERONE ANTAGONIST	SACUBITRIL	SGLT2 INHIBITOR
Yes	Yes	No	No	Yes
SHOW MORE 				

Viewing biological data

Biological data is displayed at the bottom of the patient record in a table organised by test date. Use the arrows to scroll back through the history of biological test results.

Biological data

+

ANALYSIS DATE*	< 01/27/202610/31/2025 >
REPORT DATE	01/27/202610/31/2025
SOURCE	APICRYPTINCLUSION
CREATININE mg/dL	22.722.96
eGFR (CKD-EPI) mL/min/1.73 m²	26
eGFR (MDRD) mL/min/1.73 m²	29
PROTEINURIA (24 H) mg/24h	200
POTASSIUM (K+) mEq/L	5.39
SODIUM (Na+) mEq/L	138

SHOW MORE

□

Only the data that has been entered is displayed.

Biological data can come from a number of sources:

- **Apicrypt:** data automatically uploaded via the Apicrypt messaging service. The laboratory sent us the results and we have incorporated them.
- **Newcard:** data entered by a Newcard nurse.
- **Patient:** data entered into the mobile app by the patient
- **User:** data entered by a user at the healthcare facility
- **Enrolment:** data entered by the prescriber when the patient was enrolled.

To **add biological data**, click either the 'Add biological data' button in the action menu or the button here. 

To **edit them**, you can use the pencil icon. Don't forget to click Save at the bottom of the page to save your changes.

You can change the default units in your profile:

Profile

Identity

Security

Configuration

TLS Order

Documents

Default unit

Creatinine levels

mg/L

Sodium (Na⁺)

mEq/L

Hemoglobin

g/dL

Glucose

g/L

LDLc

g/L

Total bilirubin

mg/L

Phosphorus

mg/L

Proteinuria on sample

mg/g

Magnesium

mg/L

Potassium (K⁺)

mEq/L

Bicarbonates or Alkaline Reserve (HCO₃⁻)

mEq/L

Platelets or thrombocytes

elements/mm3

Total cholesterol

g/L

Triglycerides

g/L

Calcium

mg/L

24-hour proteinuria

mg/24h

urea

g/L

RAC: Albuminuria/Creatininuria Ratio

mg/mmol

REGISTER

e. Action on the graphs

At the top right of each chart, you will find a list of possible actions:

+

−

+

👉

🏠

≡

+

−

Positive or negative zoom on the charts to view the detail of a period.

+

Selection of a specific period that you wish to enlarge on all charts.

👉

Change of the data display period.

🏠

Reset the data display period to the default: 1 month.

≡

Download each chart individually in 3 different visual formats:

- SVG format, which corresponds to a web image
- PNG format, which corresponds to a classic image
- CSV format, which corresponds to Excel format (raw data export)

f. Patient history


On the patient record, you will be able to add medical comments in the **'Medical Events'** tab. You can note instructions for our nurses, reports following consultations, changes to treatment, etc.

In turn, our team can also leave comments for you, such as the postponement of a telephone consultation with a patient or the decision not to trigger a medical alert because the patient was contacted and showed no signs of distress, etc.

This tab also contains the history of medical alerts, technical support sessions, and technical information:


The medical history includes:

- Level 1, 2 and 3 alerts that have been triggered,
- Treatment records and comments that have been added,
- Changes in the patient's medical status, including the associated dates: hospital admissions and discharges.

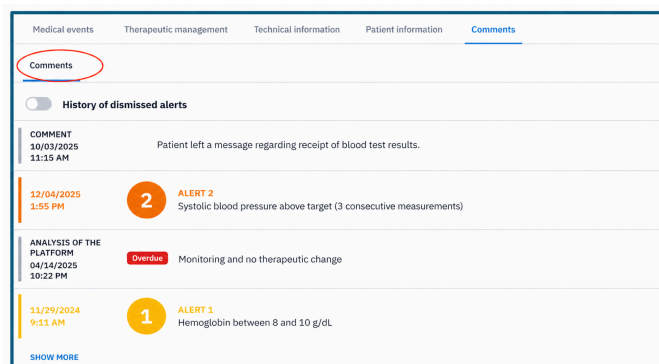The patient's therapeutic support record includes:

- The dates and topics of the sessions,
- Comments provided by the nurse in charge of therapeutic support.

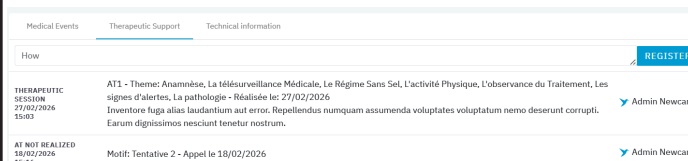The technical history allows you to:

- Find out about any difficulties the customer may have encountered when using their equipment,
- Track the steps taken by Newcard to resolve the issue.

Medical Events	Therapeutic Support	Technical information
How		REGISTER
09/04/2026 15:05	 ALERT 1 Lack of measures	MORE INFO
CREATION 10/11/2025 10:56	The patient was created	Jean Arnaud

You can add a comment to the patient's medical history or treatment history at any time.

This provides additional information to the healthcare professionals involved in the remote monitoring process. Once you have written your comment, click [REGISTER](#) to add it.


Medical events	Therapeutic management	Technical information	Patient information	Comments
Comments				
☐ History of dismissed alerts				
COMMENT 10/03/2025 11:15 AM Patient left a message regarding receipt of blood test results.				
12/04/2025 1:55 PM  ALERT 2 Systolic blood pressure above target (3 consecutive measurements)				
ANALYSIS OF THE PLATFORM 04/14/2025 10:22 PM  Overdue Monitoring and no therapeutic change				
11/29/2024 9:11 AM  ALERT 1 Hemoglobin between 8 and 10 g/dL				
SHOW MORE				

g. Generating a PDF export of remote telemonitoring data

You can generate a PDF export of the remote monitoring data by clicking this button in the patient record  .

This window will open, and you can select:


- the report period
- the sections to be included in the report

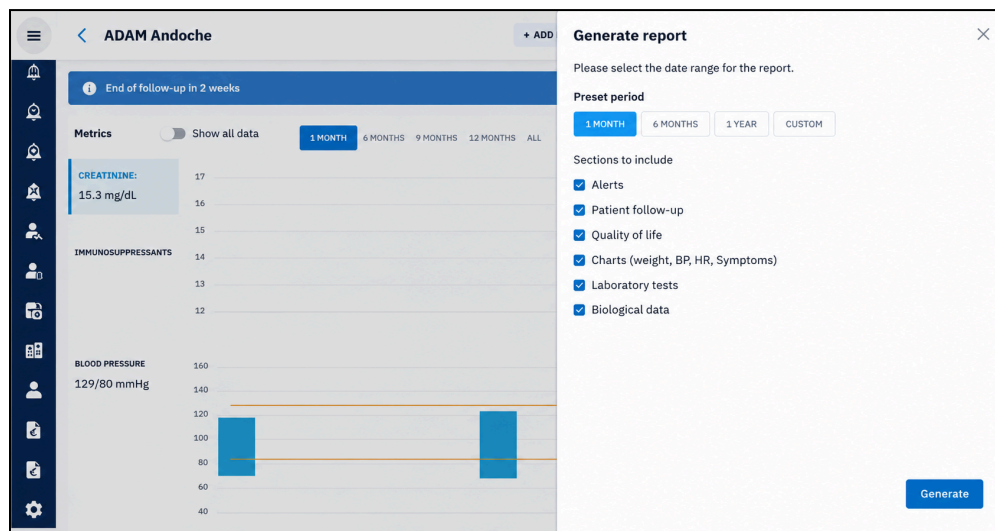Clicking the button  will download the PDF export.

It contains the patient's full medical history: their alerts, clinical data, laboratory results, medical records, adherence, and, for patients with kidney failure, their quality of life.

h. Actions on the patient's profile

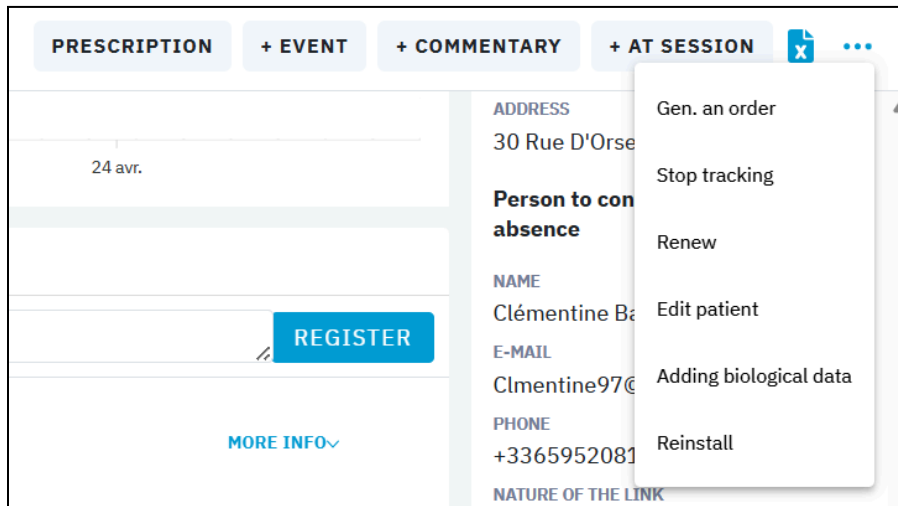At the top of the patient record, you will find several clickable action buttons that provide the following functions:

- The '**ORDO TLS**' button allows you to add a remote monitoring prescription linked to the patient's treatment period
- The '**Event**' button allows you to indicate whether the patient has been hospitalised or is absent for another reason
- The '**Comment**' button allows you to go directly to the bottom of the patient record to add a comment directly to Medical Events, Therapeutic Support or Technical Information, depending on your profile permissions
- The '**TS SESSION**' button allows you to add a new Therapeutic Support session.
- The **X** button allows you to export all your patient's data as an Excel document,
- The button  opens a drop-down menu containing several options.:
 - Generate a prescription for remote monitoring, medication or laboratory tests
 - Stop remote monitoring
 - Renew the patient's remote monitoring,
 - Edit the patient's details to amend their administrative or medical information,
 - Add laboratory test results manually,

Each button is subject to a permission.

→ Please refer to Chapter 11: Profile and User Management

i. Sending text messages

You can send messages to your patients who have a mobile phone directly via the platform. These are one-way text messages; patients cannot reply to you and your personal number will not be displayed. On the patient's record, in the right-hand column, click on the text message icon.

A window will appear asking you to enter your message in the box. The system will prepare the message by including the patient's name and your name as the signature.

You can see the number of characters allowed in this text message at the bottom right of the text box.

Then click on

SEND

The patient receives a message from you without your phone number being exposed.


This feature requires the '**Send text messages**' permission.

9. Alert management

If an alert is triggered, you will be notified via a non-personalised email and, if you wish, via text message as well. The email will contain a link allowing you to access the relevant patient's record directly.

When you log in to the platform, you are taken directly to the dashboard, which displays:

- The severity level of the alert (1, 2 or 3 for patients with kidney failure or kidney transplants, and 1 or 2 for patients with heart failure),
- The patient's name,
- The date the alert was triggered,
- The type of data that triggered the alert (weight, blood pressure, heart rate, symptoms or laboratory results, depending on the patient's condition).

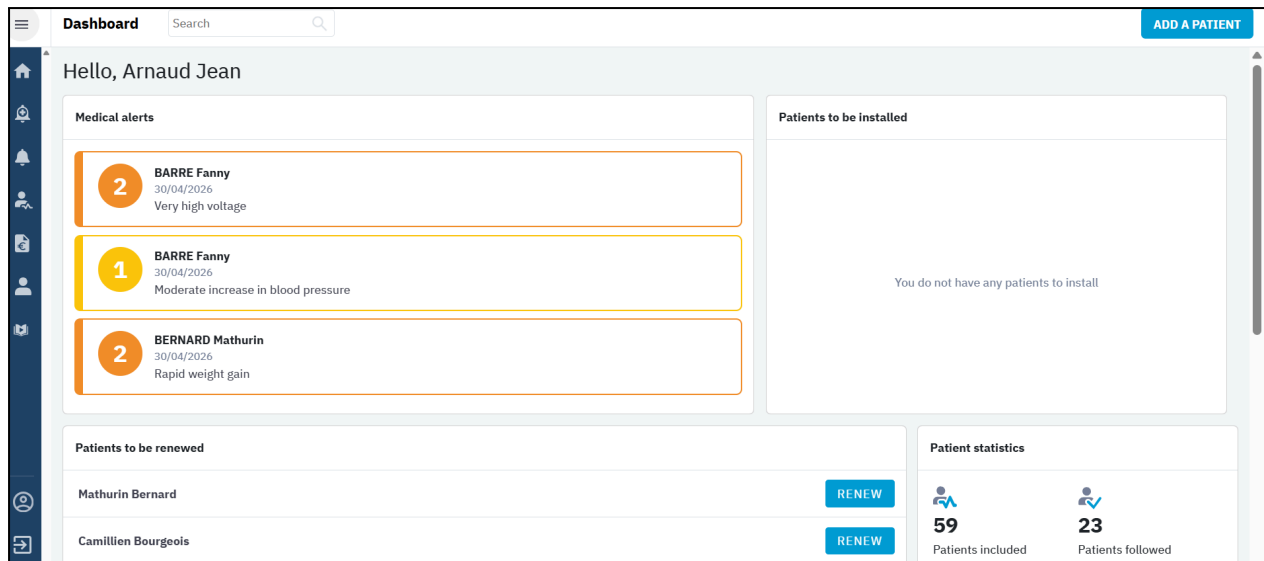
You can also view all alerts directly in the **'Medical Alerts'** tab on the left-hand sidebar.

There are three alert statuses:

- **Triggered:** the alert has been sent but has not yet been acknowledged,
- **Viewed:** you have logged into the platform and opened the relevant patient's record, but you have not acknowledged the alert,
- **In progress:** you have acknowledged an alert, but it has not yet been resolved.

To view a patient's data when an alert is active, click on their name.

The alert is displayed at the very top of their record. Below that, you will find the patient's data charts.

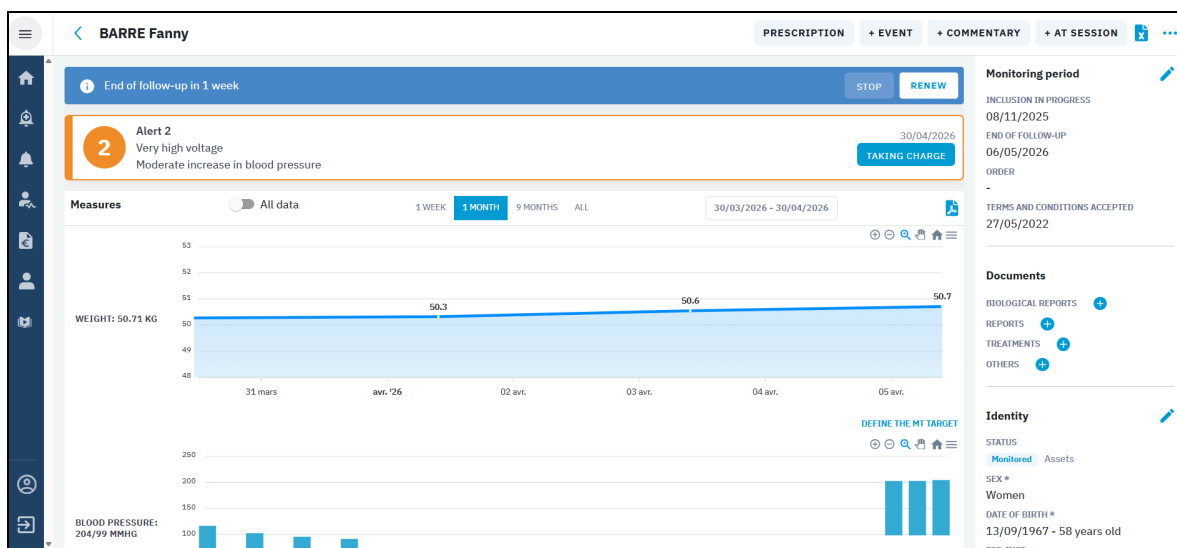Once you have reviewed the alert and the patient's details, you can indicate that you are taking on the case on the platform by clicking on.

A window will open, allowing you to specify **TAKING CHARGE** how this alert should be handled.

The screenshot shows a mobile application interface for managing alerts. The title bar at the top says '< Taking charge of the alert'. On the left is a dark blue sidebar with icons for home, notifications, and other functions. The main content area has a light blue background. It starts with 'Alert 2 Triggered'. Below that, 'PATIENT Fanny Barre' is listed. Then 'TRIGGERED BY 2 Tension très haute'. The 'TAKING CHARGE' section has 'Action taken Please select an action*' with three buttons: 'ANALYSIS ON THE PLATFORM', 'PATIENT CALL', and 'CONSULTATION'. Below this is 'Resolution of the alert' with a 'Select' dropdown menu. A 'Comment' text area follows. Then 'Pertinence de l'alerte' with five star icons. At the bottom is a 'HISTORY' section with 'Alert in progress' and 'Alert over' buttons.

Specifying the patient care actions on the platform makes it easier and more efficient for our team to filter future alerts. On the platform, you must specify several details:

- Action taken:
 - Analysis on the platform,
 - Call to the patient,
 - Consultation.
- Type of alert resolution, depending on the patient's condition:
 - Spontaneous return to normal,
 - Change in treatment
 - Hospitalisation
 - Patient unreachable
 - etc.
- A comment on the handling of this alert
- Assessment of the relevance of the alert (1 star: not very relevant, 5 stars: very relevant).

If the alert has not been resolved, click on **'Active alert'**: you can specify the action to be taken later. However, until an alert has been validated, no other alerts can be triggered.

If the alert has been resolved, click on **'Alert closed'**.

Once you have closed the alert, it will appear in the patient's medical history, as well as in the **'Medical alerts'** tab, under the **'Closed'** sub-tab.

A user's ability to view and manage alerts is subject to the **'View medical alerts'** and **'Manage medical alerts'** permissions.

These permissions are available only to care provider profiles (please refer to Section 11: Profile and User Management).

If you only have access to **'View medical alerts'**, you can see the alert banner and the alert history, but you do not have access to the **'Take charge'** button, which allows you to manage alerts.

10. Notification management

In the Notifications tab, you will find three types of information:

- the addition of a comment regarding a patient;
- the addition of laboratory results;
- the addition of a document.

The aim is to notify you as soon as the Newcard team enters new information relating to a patient, and vice versa.

Please note: once a notification is marked as 'resolved', it disappears. Therefore, if several users are associated with the same organisation, the notification will no longer be visible to others once it has been resolved.

This process also applies to notifications relating to the addition of documents and laboratory results to a patient's record.

Notifications				
Search				
All 3 Comments 1 Lab results 1 Documents 1				
ALERT	NAME	TYPE	NOTIFICATION DATE	ACTIONS
	ANDRE Judicaël	 Lab Results	04/07/2026	PROCESS
	ARNAUD Constance	 Documents	04/07/2026	PROCESS
	ADAM Alexandre	 Comments	04/07/2026	PROCESS
Showing 1 to 3 of 3 entries				« 1 »

11. Profile and user management

User and profile management allows each user within your organisation to have permissions tailored to their role. In particular, this enables you to create non-medical accounts securely. Finally, you can manage these access rights independently. With multi-user management, you can:

- Restrict access to patient data to authorised professionals,
- Provide the right information to the right people
- Tailor the system to your needs
- Deliver a consistent experience for all users

This means that each institution can customize its profiles and permissions to suit its own organisational structure.

USER MANAGEMENT: MANAGE USERS AND THEIR PERMISSIONS

Access the new **Users and Profiles** tab to manage your teams.

Users and Profiles

User and Profile Management

Users Profiles

1 ADD A CAREGIVER USER 2 ADD A NON-CAREGIVER USER

LAST NAME, FIRST NAME	PROFILE	EMAIL	PHONE	ACTIONS
BRENTEAU ALBERT Active	Disabled	a.brenteau@newcard.io	/	
BENOITON GEORGES Active	Disabled	g.benoiton@newcard.io	/	
DELANGUE CARINE Active	Disabled	c.delangue@test1@newcard.io	/	
DOGNIN ELISA Active	Referring cardiologist	e.dognin@test@newcard.io	+33678453245	
DOMI TOM Active	Billing manager	t.domi@test2@newcard.io	/	
DUPONT MONIQUE Active	Billing manager	test.factu@y@mail.com	/	
FRETTE OLIVIER Active	Cardiologist	f.olivier@mail.com	/	

RESET A PASSWORD

For users who have forgotten their password, simply click the **“Reset password”** option.
A new automatic email will be sent.

User and Profile Management

Users Profiles


1 ADD A CAREGIVER USER 2 ADD A NON-CAREGIVER USER

LAST NAME, FIRST NAME	PROFILE	EMAIL	PHONE	ACTIONS
BRENTEAU ALBERT Active	Disabled	a.brenteau@newcard.io	/	
BENOITON GEORGES Active	Disabled	g.benoiton@newcard.io	/	
DELANGUE CARINE Active	Disabled	c.delangue@test1@newcard.io	/	
DOGNIN ELISA Active	Referring cardiologist	e.dognin@test@newcard.io	+33678453245	
DOMI TOM Active	Billing manager	t.domi@test2@newcard.io	/	
DUPONT MONIQUE Active	Billing manager	test.factu@y@mail.com	/	
FRETTE OLIVIER Active	Cardiologist	f.olivier@mail.com	/	

1 ADD A CAREGIVER USER

- 1 Search for the healthcare professional in the RPPS directory
(by last name, first name or RPPS number)
- 2 Select a profile from the list of caregiver profiles

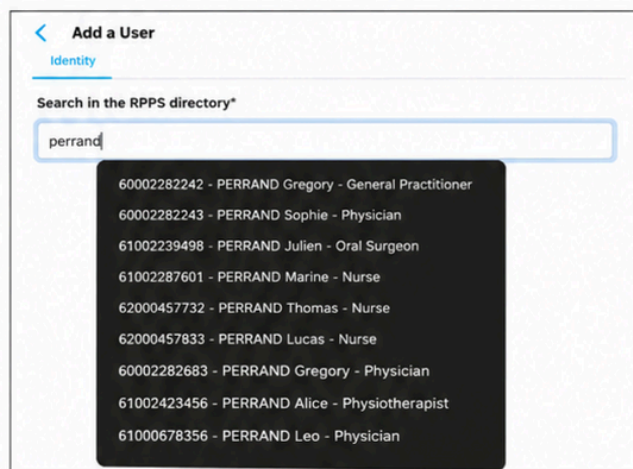
Add a User
Identity

Search in the RPPS directory*

perrand

- 60002282242 - PERRAND Gregory - General Practitioner
- 60002282243 - PERRAND Sophie - Physician
- 61002239498 - PERRAND Julien - Oral Surgeon
- 61002287601 - PERRAND Marine - Nurse
- 62000457732 - PERRAND Thomas - Nurse
- 62000457833 - PERRAND Lucas - Nurse
- 60002282683 - PERRAND Gregory - Physician
- 61002423456 - PERRAND Alice - Physiotherapist
- 61000678356 - PERRAND Leo - Physician

- 3 Complete the requested information (email, phone number, etc.)

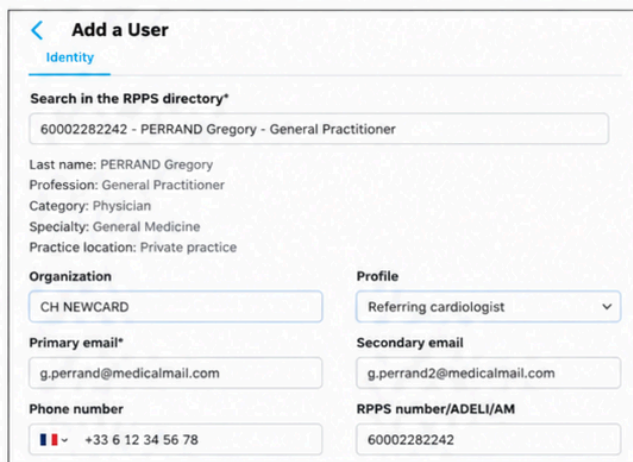Add a User
Identity

Search in the RPPS directory*

60002282242 - PERRAND Gregory - General Practitioner

Last name: PERRAND Gregory
Profession: General Practitioner
Category: Physician
Specialty: General Medicine
Practice location: Private practice

Organization
CH NEWCARD

Profile
Referring cardiologist

Primary email*
g.perrand@medicalmail.com

Secondary email
g.perrand2@medicalmail.com

Phone number
+33 6 12 34 56 78

RPPS number/ADELI/AM
60002282242

4 SAVE

- The user will automatically receive their password by email at the address provided during account creation.

Please note: If none of the existing profiles meets your needs, you can create a new one ([see page 6](#)).

Profile permissions apply to all associated users:

- To assign specific permissions to a user, you must create a dedicated new profile.

2

ADD A NON-CAREGIVER USER

- 1 Choose a profile from the list of non-caregiver profiles
(e.g.: *administrative, billing*)
- 2 Complete the requested information

The screenshot shows the 'Add a User' interface. On the left is a sidebar with navigation links: Dashboard, Medical alerts, Notifications, Patients, Invoicing, Users and Profiles, and Help. Below these is a section for 'CH NEWCARD' with user information: 'Elsa Dognin, Referring cardiologist' and a 'Log out' button. The main area is titled 'Add a User' and has a 'Identity' sub-section. It contains the following fields: 'Organization / Facility' (CH Newcard), 'Last name*' (Doe), 'First name*' (John), 'Primary email*' (john.doe@medicalmail.com), 'Phone number' (France (+33)), and 'Profile' (Billing manager). A red arrow points to the 'SAVE' button at the bottom left.

3 SAVE

- The user will automatically receive their password by email at the address provided during account creation.

Please note: If none of the existing profiles meets your needs, you can create a new one ([see page 6](#)).

Profile permissions apply to all associated users:

- To assign specific permissions to a user, you must create a dedicated new profile.

Deactivate a user

Assign the **"Inactive"** profile to remove all access rights.

[<](#) **REFERRING ENTITY Newcard** **CAREGIVER** Active

Identity

Configuration

Documents

Default unit

Search in the RPPS directory* To match a user's identity, please search by their RPPS directory.

Last name*

réfèrent

First name*

Newcard

Organization / Facility

CH Newcard

Primary email*

ref.cardio@newcard.io

Phone number

 +33 3 27 43 04 61

ADELI/AM number

Profile

Inactive

Cardiologist (referring)

Cardiologist

Health ID in cardiology

Telemonitoring health ID in cardiology

Inactive


Cardiologist (referring)

Dietitian

RPPS number

Edit a profile (caregiver or non-caregiver)

- 1 Select the profile assigned to the user
- 2 Check/uncheck permissions according to their role

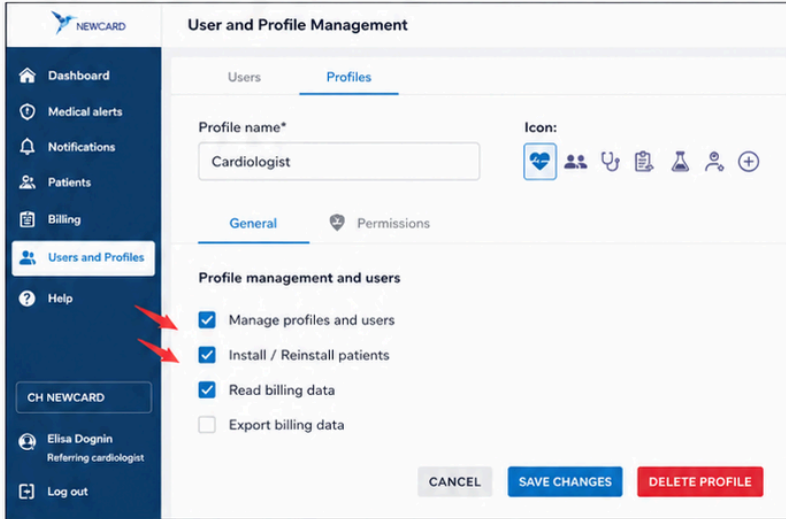Warning:
any modification to a profile's permissions will automatically apply to all users associated with this profile.


 To assign specific permissions to a single user, you must create a dedicated new profile.

- 3 Save the changes

Tip : regularly review permissions to stay aligned with service needs.

Create a new profile (caregiver or non-caregiver)

- 1 Enter a profile name
- 2 Choose an icon
- 3 Select permissions

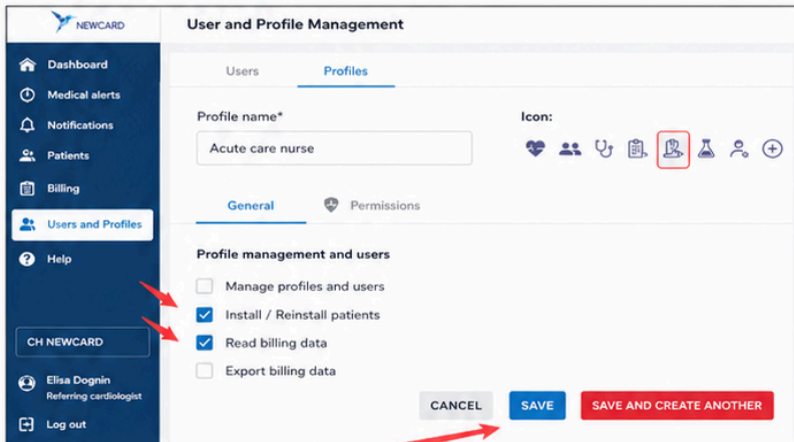- 4 **SAVE**

→ Useful for adapting the platform to new roles or services.

Configure notifications and alerts

Each user can:

- Choose the types of alerts they receive
- Define the patients concerned and the delivery channel

Tip: configure alerts when adding a patient so you don't miss anything.

	Primary email	Secondary email	SMS
Medical alerts	☒	☐	☐
Patient hospitalization	☐	☐	☐

Show only patients for whom I am the prescriber	☐

Filter the patient view

In centers with multiple prescribers, all patients can be visible.

This filter allows you to display only the patients for whom you are the prescriber for clearer navigation.

	Primary email	Secondary email	SMS
Medical alerts	☒	☐	☐
Patient hospitalization	☐	☐	☐

Show only patients for whom I am the prescriber	☒

Only users who have the **'Profile and User Management'** permission in their profile can access the **'Profiles and Users'** tab. We recommend that the organisation's administrator be the one with this permission.

For an individual private practitioner, this tab is not necessary; they will have all rights by default and will not need to configure other people's rights.

List of rights and associated features

Types of rights	Name of the right	Rights of healthcare professionals	Rights of Non-Clinical Professionals	Definition
General	Profile and user management	x	x	The 'Users and Profiles' tab is active. You can create and edit users and profiles within your institution
General	Settle in / Settle back in	x	x	Option to register the patient using this account by entering their email address and password in the patient app; display of the list of patients to be registered; and option to click the 'Re-register' button on the patient record (provided you have access to the patient record with the necessary permissions to enable cardiac or renal monitoring)
General	Billing tab: reading	x	x	The 'Invoicing' tab is active. You can view invoices but cannot make any changes (issue invoices, change the amount, renew, return to the 'To be invoiced' tab, etc.). The invoicing statistics are displayed on the dashboard.
General	Billing tab: invoice	x	x	The 'Invoicing' tab is active. You can view invoices and make changes (issue an invoice, change the amount, renew, return to the 'To be invoiced' tab, etc.). Invoicing statistics are available on the dashboard.
Cardiology	Enable the display of patients with heart failure	x	x	The 'Patients' tab displays only the cardiac patients at the facility. Access is limited to viewing the patient record

Cardiology	May prescribe remote monitoring	x		The user account name appears in the drop-down list of remote monitoring prescribers when adding or editing a patient. It is possible to use the pre-signed prescription in the name of that prescriber. If the user is authorised to prescribe remote monitoring, they automatically have all other rights relating to the condition.
Cardiology	The inclusion of the patient	x		The 'Add a patient' button is available on the dashboard and in the Patients list
Cardiology	Readmit the patient	x		The list of patients due for renewal is available on the dashboard. The 'Renew' option is available on the dashboard; in the patient lists, you can see the 'Due for renewal' list and the renewal banner; the 'RENEW' button is available in the patient lists; and the 'Renew' option is available in the context menu (three dots).
Cardiology	Stop the patient	x		Access to the 'Stop tracking' button in the context menu of the patient record
Cardiology	Report an incident (hospitalisation, absence, etc.)	x		Access to the '+ Event' button, which allows you to record an absence, a hospital admission or a discharge from hospital
Cardiology	Generate and submit a prescription (telemonitoring, organic, therapeutic, etc.)	x		Access to the 'Generate prescription' button in the context menu + upload a prescription in the semesters screen and in the patient report
Cardiology	View medical alerts	x		View the Medical Alerts tab, see the list of medical alerts in this list, on the patient record or on the dashboard, but not be able to manage them
Cardiology	Handling medical alerts	x		When viewing the Medical Alerts tab, whether you are viewing the list of medical alerts here, on the patient record or on the dashboard, the 'Take charge' button is available
Cardiology	Text message	x		Access to the small "Text message" bubble, which allows you to send a text

				message from the platform to a mobile phone
Cardiology	View documents and prescriptions	x	x	Provides access to view documents (Documents section) and orders (Orders section, including access to the semesters tab) + in the notifications area (on the dashboard and in the Notifications tab), the Documents section is visible
Cardiology	Upload a document to the patient file	x	x	Option to add a document, open the document viewer and upload a document (Documents section) and timetables (Timetables section, access to the Terms tab) + in the notifications area (on the dashboard and in the Notifications tab) you can see the Documents section
Cardiology	Edit the patient's administrative information	x	x	Access to the 'Edit patient' button in the context menu of the patient record; access to the pencil icon to edit the patient in the Identity, Contact and Care Team sections. For administrative editing, only the Identity, Contact and Patients linked tabs are available for editing
Cardiology	Edit the patient's medical information	x		Access to the 'Edit patient' button in the context menu of the patient record, and access to the pencil icon to edit the patient in the 'Follow-up period' and 'Medical data' sections. For medical editing, only the 'Inclusion' and 'Medical data' tabs are available for editing. Access to the button to add laboratory results
Cardiology	View medical events and add comments	x		Access to comments in the Notifications section on the dashboard and in the Notifications tab. Access to medical events in the patient record, with the option to add a comment directly in the section or via the '+Comment' button
Cardiology	View technical information and add comments	x	x	Access to technical information on the patient record, with the option to add a comment directly in the section or via the '+Comment' button

Cardiology	Visualisation of therapeutic support sessions	x		Access to Therapeutic Support sessions. You can view the sessions but cannot add a TS session using the '+TS Session' button (only viewing of the TS section is available). However, you can write a comment
Cardiology	Add therapeutic support sessions	x		Access to Therapeutic Support sessions. You can read the sessions and add a comment directly in the section or using the +Comment button. You can also add a TS session using the +TS Session button
Nephrology	Enable display of patients with renal failure	x		The ' Patients ' tab shows only the nephrology patients at the facility. Access is limited to viewing the patient record
Nephrology	May prescribe remote monitoring	x		The user account name appears in the drop-down list of remote monitoring prescribers when adding or editing a patient. It is possible to use the pre-signed prescription in the name of that prescriber. If the user is authorised to prescribe remote monitoring, they automatically have all other rights relating to the condition.
Nephrology	The inclusion of the patient	x		The ' Add a patient ' button is available on the dashboard and in the Patients list
Nephrology	Readmit the patient	x		The list of patients due for renewal is available on the dashboard. The 'Renew' option is available on the dashboard; in the patient lists, you can see the 'Due for renewal' list and the renewal banner; the 'RENEW' button is available in the patient lists; and the 'Renew' option is available in the context menu (three dots).
Nephrology	Stop the patient	x		Access to the 'Stop following' button in the context menu of the patient record
Nephrology	Report an event (Hospitalized, absent,...)	x		Access to the "+ Event" button, which allows you to enter an absence, a hospitalization, or a return from hospitalization

Nephrology	Generate and file a prescription (remote monitoring, bio, therapeutic...)	x		Access to the button ' Generate a prescription ' in the context menu + upload a prescription on the semesters screen and in the patient edition inclusion
Nephrology	View medical alerts	x		View the Medical Alerts tab, see the list of medical alerts in this list, on the patient record or on the dashboard, but not be able to manage them
Nephrology	Handling medical alerts	x		When viewing the Medical Alerts tab, whether you are viewing the list of medical alerts here, on the patient record or on the dashboard, the ' Take charge ' button is available
Nephrology	Text message	x	x	Access to the small " text message " bubble allowing you to send an tet from the platform to a mobile phone
Nephrology	View documents and prescriptions	x	x	Gives the right to read documents (document section) and orders (orders section, access to the semesters tab) + in the notification area (on the dashboard and in the Notifications tab), you can see the Documents area
Nephrology	Upload a document to the patient file	x	x	Possibility to add a document, open the document viewer, and upload a document (document section) and orders (order section, access to the semesters tab) + in the notification area (on the dashboard and in the Notifications tab), you can see the Documents area.
Nephrology	Edit the patient's administrative information	x	x	Access to the ' Edit patient ' button in the context menu of the patient file, access to the pencil button to edit the patient on the Identity, Contact and Supervision section. For the administrative edition, only the Identity, Contact, and Users tabs related to patients are available on the edition.
Nephrology	Edit the patient's medical information	x		Access to the 'Edit patient' button in the context menu of the patient record, access to the pencil button to edit the patient in the follow-up period and medical data sections. For medical

				publishing, only the Inclusion and Medical Data tabs are available for editing. Access to the button to add biological data
Nephrology	View medical events and add comments	x		Access to Comments in the Notifications area on the dashboard and in the Notifications tab. Access to medical events from the patient file with the possibility of adding a comment directly in the section or with +Comment
Nephrology	View technical information technology and add comments	x	x	Access to the technical information of the patient file with the possibility to add a comment directly in the section or with +Comment
Nephrology	Visualisation of therapeutic support sessions	x		Access to Therapeutic Support sessions. You can view the sessions but cannot add a TS session using the '+TS Session' button (only viewing of the TS section is available). However, you can write a comment
Nephrology	Add therapeutic support sessions	x		Access to Therapeutic Support sessions. You can read the sessions and add a comment directly in the section or using the +Comment button. You can also add a TS session using the +TS Session button

12. Invoice


In the Newcard platform, there are two billing systems:

- automated billing: sending electronic health records - FSE - directly from the platform - mainly used by private prescribers.
- the billing follow-up, which is a list of invoices to be issued and made when the invoices are uploaded to the business software in parallel or by paper care sheet – mainly used for centers.

In both cases, you can access this feature via the **‘Billing’** tab if you have the right to **‘Billing tab: reading’** and **‘Billing tab: invoicing’** in your profile.

If you only have the right to read, you can view invoices without being able to generate them or indicate them as done.

If you have both rights, then with automated invoicing you can create invoices, cancel invoices and upload them. With billing tracking, you can indicate that the invoices have been made.

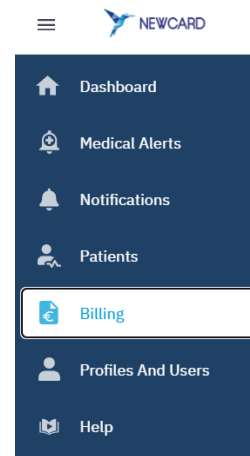a. Automated billing (private practitioners)

1. Installation and configuration process for automated billing

a) Glossary

Term	Definition
jFSE	name of the billing engine
FSE	Feuilles de Soins Électronique = Electronic Treatment Forms
PS	Professionnel de Santé = Healthcare Professional

b) Prerequisite: CPS card reader

Before you begin, please ensure that:

- The CPS card reader is properly connected to the computer
- The CPS card is inserted into the reader


c) Installing the card reader software

Please note: Do not carry out this procedure if you are using any of the following software:
Lifeline - Crossway - MLM

Downloading the module

1. Log in to the Newcard platform.
2. If an orange banner appears at the top of the screen, click on it.
3. Otherwise, go to:
 - Profil
 - Settings
 - Click on the download icon (downward arrow)

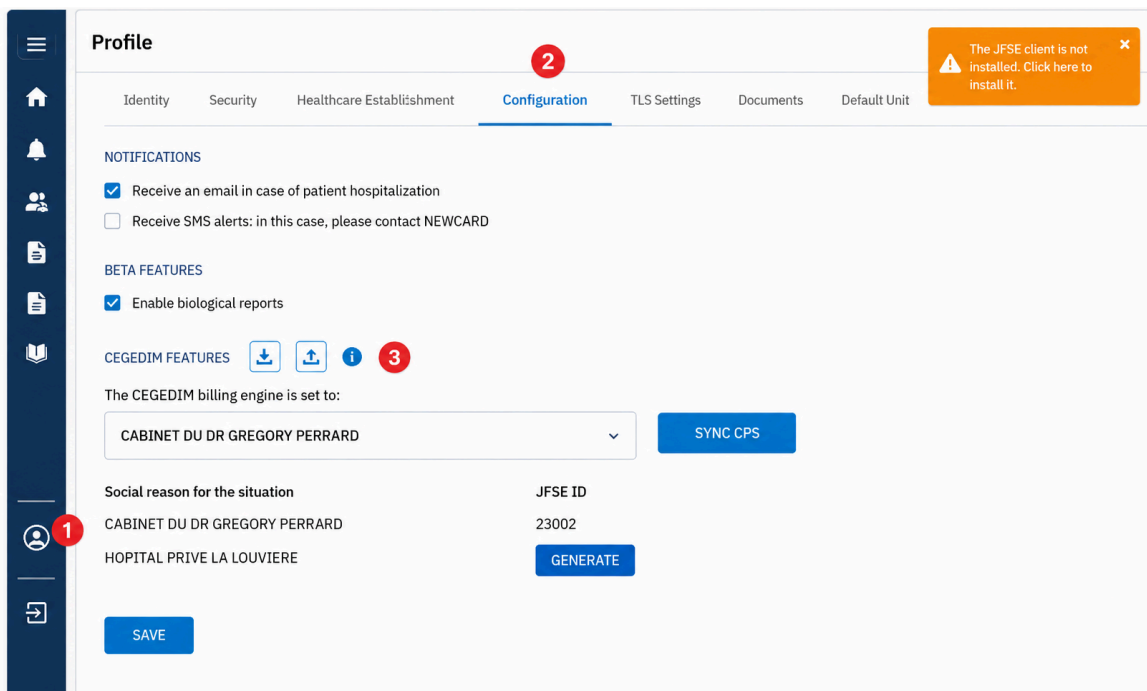Profile

Identity Security Healthcare Establishment **Configuration** TLS Settings Documents Default Unit

NOTIFICATIONS

- ☒ Receive an email in case of patient hospitalization
- ☐ Receive SMS alerts: in this case, please contact NEWCARD

BETA FEATURES

- ☒ Enable biological reports

CEGEDIM FEATURES [Download] [Upload] [Info] **3**

The CEGEDIM billing engine is set to:

CABINET DU DR GREGORY PERRARD [Dropdown] [SYNC CPS]

Social reason for the situation JFSE ID

CABINET DU DR GREGORY PERRARD 23002


HOPITAL PRIVE LA LOUVIERE [GENERATE]

[SAVE]

The JFSE client is not installed. Click here to install it.

A download window will open:

- **Windows** : click on the Windows icon to download the executable file.
- **Apple** : click on the Apple icon to download the executable file.

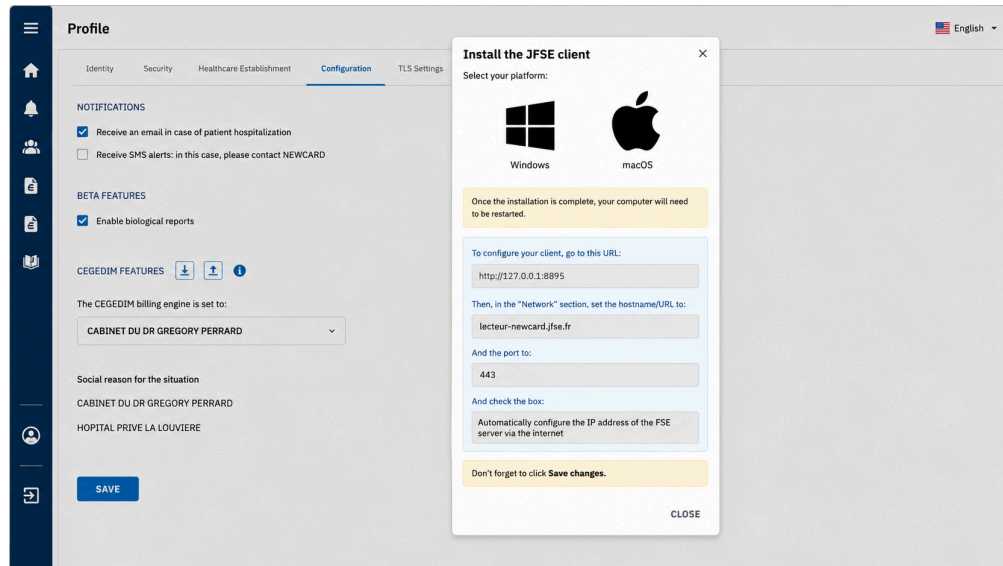
Installation

1. Open the downloaded file (usually located in the Downloads folder).
2. Start the installation.
3. A black window will appear showing the various installation steps.
4. Allow the installation to complete until the message “Press any key to continue” appears.
5. Press the space bar (or any other key).
6. The computer will restart automatically.

d) Connection settings

Go to settings

1. In the notification bar (at the bottom right of the screen), look for the white icon labelled “**SESAM Vitale Reader**”,
2. Right-click on this icon,
3. Select “**Settings**”.

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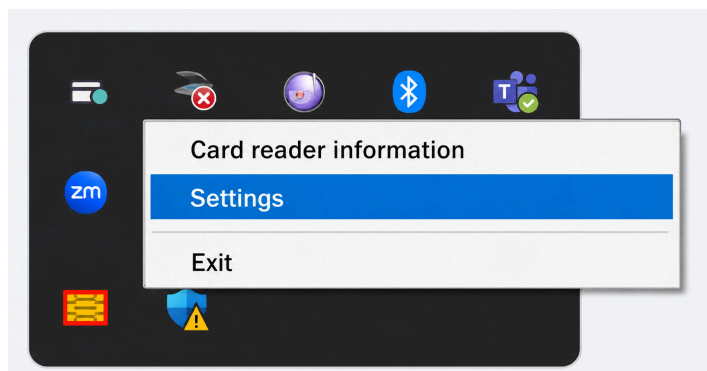1. In the window that opens, select the **Networks** tab from the menu on the left.
2. Enter the following settings:

Parameter	Value
IP address	lecteur-newcard.jfse.fr
Server port	443
TOMCAT port	8895

The screenshot shows the JFSE Networks configuration window. On the left is a sidebar with the JFSE logo, 'SESAM-VITALE READER', a 'Connected to server' status, and a menu with 'Readers', 'Networks' (selected), 'Maintenance', 'Mapping', and 'Libraries'. Below the menu are 'VERSIONS' and 'ENVIRONMENT' details. The main area is titled 'Networks' and contains three panels: 'JFSE Server', 'Local Servers', and 'Proxy'. The 'JFSE Server' panel has fields for IP Address / Hostname (lecteur-newcard.jfse.fr), Port (443), Communication ID (de8674b4-e1eb-4209-9b2c-8f5e9cadaf75), and Compatibility mode (2). It also has three toggle switches: 'Encrypt exchanged data' (checked), 'Automatically configure FSE server IP address via internet' (unchecked), and 'Display animations during reader actions' (checked). The 'Local Servers' panel has three sections: 'Local web server' (checked, port 8083), 'XML-RPC service' (unchecked, port 8080), and 'TOMCAT web server' (checked, port 8895). The 'Proxy' panel has fields for Host (Host) and Port (0). A 'Save changes' button is in the top right corner.

Parameter	Value
IP Address / Hostname	lecteur-newcard.jfse.fr
Port	443
Communication ID	de8674b4-e1eb-4209-9b2c-8f5e9cadaf75
Compatibility mode	2
Encrypt exchanged data	checked
Automatically configure FSE server IP address via internet	unchecked
Display animations during reader actions	checked

Local Servers	Port
Local web server	8083
XML-RPC service	8080
TOMCAT web server	8895


Proxy	Port
Host	0

4. Click **Save**.

e) Checking that everything is working properly


After setup:

1. Check that the CPS card is properly inserted into the reader.
2. Check the icons in the notification bar.

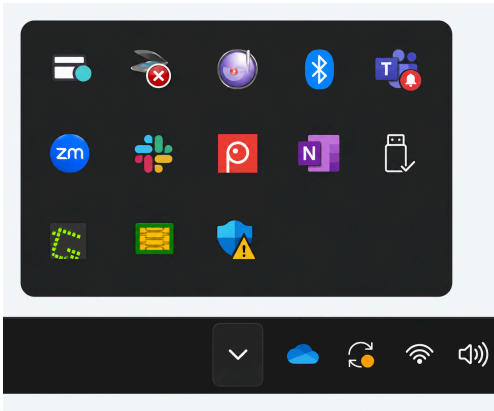The configuration is correct when:

- The drive icon displays a **blue dot**.
- The jFSE icon displays a **green outline**.

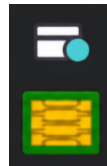f) Setting up your account on the platform

CPS card synchronisation

1. Log in to Newcard.
2. In the left-hand menu, click on your profile.
3. Open the **Settings** tab.
4. Click on **Synchronise CPS**.

Profile

IdentitySecurityHealthcare InstitutionConfigurationTLS SettingsDocuments

NOTIFICATIONS

☐ Receive an email in case of patient hospitalization

☐ Receive SMS alerts: in this case, please contact NEWCARD

CEGEDIM FEATURES [↓](#) [i](#)

The CEGEDIM billing engine is set based on the situation:

Deactivated


SYNC CPS

Social reason for the situationJFSE ID

SAVE

Enter your CPS card number


1. Enter your CPS card number.
2. Confirm.

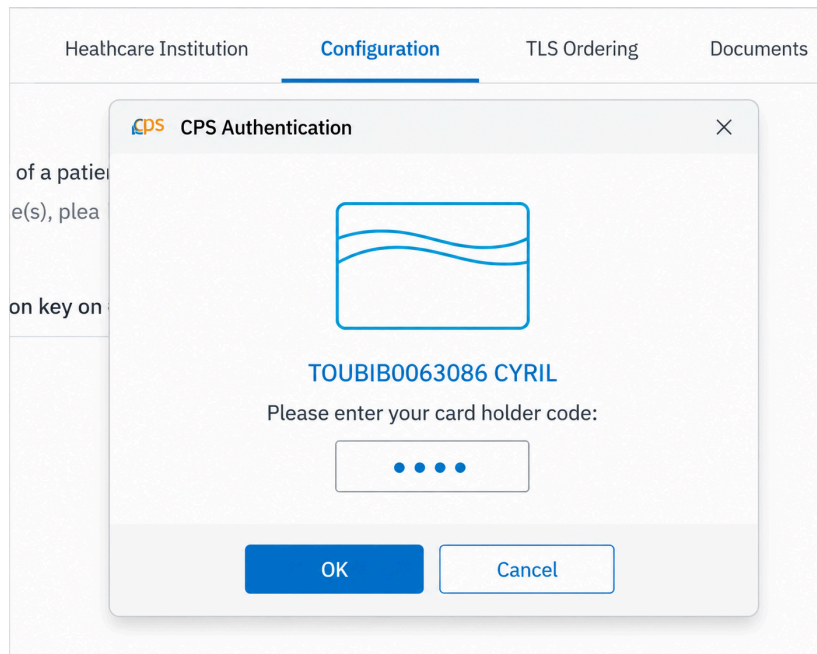The screenshot shows a web application interface with tabs: Healthcare Institution, Configuration (selected), TLS Ordering, and Documents. A modal dialog titled "CPS Authentication" is open. It features a card icon, the text "TOUBIB0063086 CYRIL", and a prompt "Please enter your card holder code:". Below the prompt is a text input field with four dots. At the bottom are "OK" and "Cancel" buttons.

Generating the jFSE ID

1. In the **jFSE ID** column, click **Generate** for the desired billing situation.

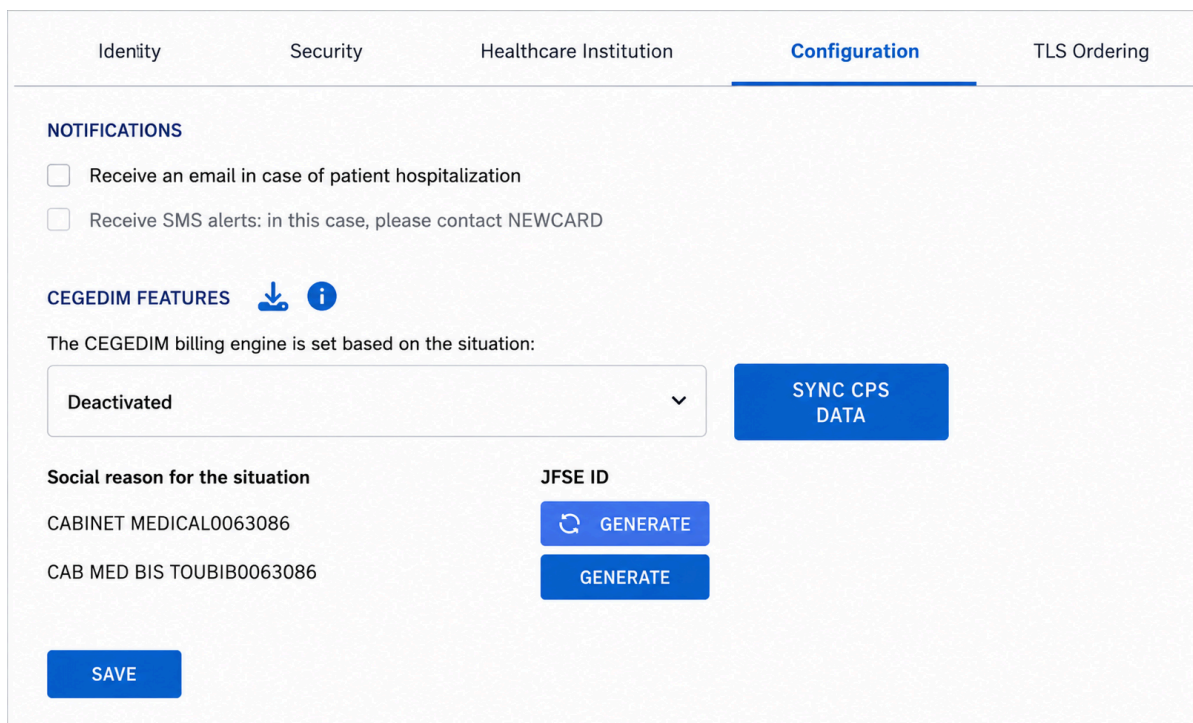The screenshot shows the "Configuration" tab in a web application. It includes sections for "NOTIFICATIONS" (with checkboxes for email and SMS alerts), "CEGEDIM FEATURES" (with a download icon and an info icon), and a dropdown menu for "The CEGEDIM billing engine is set based on the situation:" currently set to "Deactivated". A "SYNC CPS DATA" button is next to the dropdown. Below, under "Social reason for the situation", are two rows: "CABINET MEDICAL0063086" and "CAB MED BIS TOUBIB0063086". To the right of these is a "JFSE ID" column with two "GENERATE" buttons. A "SAVE" button is at the bottom left.

- From the drop-down list, select the billing status used by the healthcare professional.

Profile

Identity Security Healthcare Institution **Configuration** TLS Ordering

NOTIFICATIONS

☐ Receive an email in case of patient hospitalization

☐ Receive SMS alerts: in this case, please contact NEWCARD

CEGEDIM FEATURES  

The CEGEDIM billing engine is set based on the situation:

Deactivated

Deactivated
CABINET MEDICAL0063086
CAB MED BIS TOUBIB0063086

SAVE

SYNC CPS DATA

JFSE ID

113
114

- Click **Save**.

g) Test of creating and sending an invoice

Creating an invoice

- Go to the **Invoicing** tab.
- Click on **Create** an invoice.

When the first invoice is generated, several configuration questions are displayed.

 **ESSAI ADELE**
35 years old - 10/22/1989
1 42 01 62 965 165 02 RÉGIME 01, CAISSE 349, CENTER 9881
DECLARATION OF EXEMPTION FROM COST SHARING (CARTE VITALE): NO

CYRIL TOUBIB0063086
National ID: 001030865
GENERAL PRACTITIONER

Questions

Please answer the following questions to continue

Are the treatments related to a long-term condition (LTC)?

☒ Yes ☐ No

Choice of coverage plan

☐ No coverage
☒ OPTAM
☐ OPTAM-ACO

Select the mode of operation for the healthcare insurance scheme

☒ Standard SESAM-Vitale procedure
☐ Payment coordinated by the healthcare insurance organization

Validate

Recommended answers:

- **Non-contracted doctor:** : do not select a contract.
- **Contracted doctor:** always select:
 - **OPTAM**
 - **SESAM-Vitale**

Electronic submission

Once the invoice has been approved:


1. Go to the **'To be submitted'** tab.
2. Submit the invoice electronically.

Automated invoicing is now up and running.

2. Automated invoicing

Following an installation appointment for the CPS card reader software, you can use automated billing - only for private practitioners).

For this, make sure that the card reader is launched and that your CPS card is in the reader.

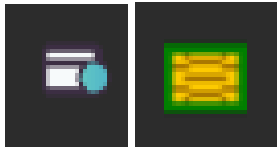It should be noted that the increase in the package only applies to patients with heart failure at present.

Step 1 - Identification of the flat-rate

When adding or renewing a patient, specify the amount of the flat-rate.

- ▶ **€28 - Social flat-rate** : This amount is automatically populated.
- ▶ Three flat-rates are available depending on the patient's health status.
For more information, please refer to the platform.

Social flat-rate - TVB - €28 / month ⓘ


Flat-rate coefficient

Social flat-rate - TVB - €28 / month
Level 1 - TVB 2 - €56 / month
Level 2 - TVB 2.5 - €70 / month

Step 2 - Invoicing

1) Access the invoicing section

- ▶ From the dashboard, go to the Invoicing section using the left-hand sidebar.
- ▶ You will find the list of invoices to be processed for each patient, with details for the period, status and balance.

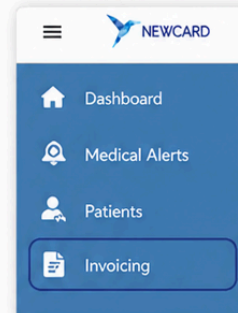2) Choose the billing mode: monthly or in series

- ▶ **Monthly**: Bill each month individually.
- ▶ **In series**: Bill several months in a single operation.
Use the button: Display only complete semesters to bill for 6 months.

3) Review your invoices

Allows you to modify the flat-rate →

Allows you to cancel an invoice if the patient is non-compliant →

4) Create your invoices

Click this button to create an invoice:

CREATE INVOICE

Is the care related to long-term condition (LTC) management?

☒ Yes ☐ No

All fields are pre-filled based on the file information.
It will remain in draft status until validated on the right.

→

Validate

A second validation is required →

Continue validation

All these steps must be completed for each patient.

Step 3 - Electronic transmission

Click the top right button “**Electronic transmission**” →

To invoice **Electronic transmission** Invoiced

Click the **Electronic transmission** button at the bottom right


Electronic transmission is automatic.
You must enter your CPS card PIN code.

Step 4 - Track the status of your invoices

Click the **DASHBOARD** button
at the top right to access your invoices.

Click the **PDF** button to access
the “Select all” checkbox for the displayed batches.

Click the right button to access
the CSV export of all displayed invoices.

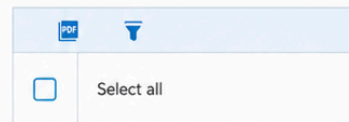ADDITIONAL INFORMATION

Report a missing prescription or the need to renew a patient.
By hovering over the icon, an explanation is displayed.

Indicates that the patient is deceased.

Export invoices in Excel format with period selection:
NEW FEATURE.

Provider = Newcard
Operator = Physician

Registration required
with the ARS
to invoice

For more information,
visit the tutorial page.

An Excel export of all invoices is available with the button located  at the top right of the page.

After specifying the export period, an excel file is generated. This Excel includes all the invoices on screen over a desired period: there is a tab with the invoices to be invoiced, a tab with the invoices to be electronically transmitted, and a tab with the invoices made.

b. Billing tracking list

When billing is not automated, it means that you want to use your billing software to bill for remote monitoring operations.

**PRACTICAL GUIDE**

HOW TO BILL FOR MEDICAL TELEMONTORING

via the Newcard platform

Step 1 - Identification of the flat-rate

When adding or renewing a patient, specify the amount of the flat-rate.

- **€28 - Social flat-rate** : This amount is automatically populated.
- Three flat-rates are available depending on the patient's health status, for more information, please refer to the platform.

Social flat-rate - TVB - €28 / month

Flat-rate coefficient

Social flat-rate - TVB - €28 / month

Level 1 - TVB 2 - €56 / month


Level 2 - TVB 2.5 - €70 / month

Step 2 - Invoicing

1) Access the invoicing section

- From the dashboard, go to the Invoicing section using the left-hand sidebar.
- You will find the list of invoices to be processed for each patient, with details for the period, status and balance.

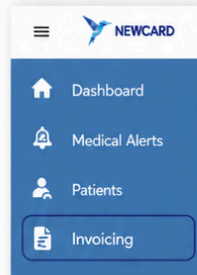2) Choose the billing mode: monthly or in series

- **Monthly**: Bill each month individually.
- **In series**: Bill several months in a single operation.
Use the button: Display only complete semesters to bill for 6 months.

Attention !

- Make sure the patient meets the observation conditions before proceeding with invoicing.
- In case of non-compliance, you can delete the invoice line using the trash icon (this will be recorded in the Invoicing tab).

3) Invoice your telemonitoring services

1 line on the tracking screen = 1 30-day period = 1 line on the "FSE" or 1 line on the paper "CERFA" form

FSE*

In your billing software
Invoicing of up to 6 periods maximum

Invoice date	Invoice code	Amount
07/09/2023	TVB	€28
07/10/2023	TVB	€28
06/11/2023	TVB	€28
06/12/2023	TVB	€28
05/01/2024	TVB	€28
04/02/2024	TVB	€28

! This does not represent all care sheet lines used ("electronic care sheet").

OR


CERFA

Invoicing of up to 4 periods maximum

CARE PROVIDED			
Dates		Code	Amount
From	To		
07/09/2023	TVB		€28
07/10/2023	TVB		€28
06/11/2023	TVB		€28
06/12/2023	TVB		€28
TOTAL			€112

162

4) Check "Invoiced" in the Newcard platform

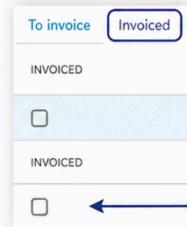In this tab, you can find the billing periods that have already been invoiced.

Check the INVOICED box after the invoice has been processed.

Attention!

- Do not check "Invoiced" for telemonitoring if it is for billing tracking in the platform. To bill, you must use your billing software.

ADDITIONAL INFORMATION

Modifications: allows you to modify the flat-rate amount.

Signals an issue related to a missing prescription or the need to renew a patient. Hover over the icon to view the explanation.

Indicates that the patient is deceased.

Export invoices in Excel format with period selection!

Provider = Newcard
Operator = Physician

Registration required
with the ARS to invoice

FINES Number (Newcard)
Acquirer number* :
59 26 57 118

Prescriber code =
your ADEL code (CDOM)**

**This code is assigned by
your local department.

For more information,
visit the tutorial page.

Status 'TO BE BILLED' / 'TO TELETRANSMIT'

When a month of remote monitoring has elapsed, you will find the name, date of birth, NIR, patient status, prescription date, prescribing physician, and remote monitoring period in the **'To be billed' / 'To teletransmit'** tab.

You can filter on all invoices for a patient by entering their first and last name in the search bar, or on all invoices for a doctor by entering their first and last name.

For deactivated patients, the invoice is stopped at the end of the current month, while for deceased patients the invoice is stopped on the date of death (a red pictogram warns you that the patient has died).

If there is no order filed on the platform and validated by a Newcard administrator, then an orange pictogram appears. In this case, you cannot check "Billed" on this tracking screen.

Billing Search ☐ Show completed billing periods						
To bill		Billed		Total outstanding: €232,966		
INVOICE	LAST NAME, FIRST NAME	DATE OF BIRTH	NATIONAL ID	PAYER CODE	REFERRING CH. NEWCARD	
☐	X-TEST MAEVA	12/06/1956	14444 5555555555	/	Referring Ch. Newcard	
INVOICE	ORDER DATE	INVOICE DATE	OBSERVATION	STATUS	CODE - AMOUNT	
☐	12/13/2023	02/10/2024	40%	Deactivated 06/09/2024	TVB - €28	
☐	12/13/2023	03/11/2024	100%	Deactivated 06/09/2024	TVB - €28	
☐	12/13/2023	04/10/2024	40%	Deactivated 06/09/2024	TVB - €28	
☐	12/13/2023	05/10/2024	80%	Deactivated 06/09/2024	TVB - €28	
☐	12/13/2023	06/09/2024	80%	Deactivated 06/09/2024	TVB - €28	
☐	Z-TEST ALEXIS	11/08/1956	fake SSnum	/	Referring Ch. Newcard	

You can proceed with billing for the month of remote monitoring as usual:

- Electronically with the TVB key, amount of €28 (excluding surcharge*);
- By paper care sheet.

Remote monitoring can be billed at the end of each month or in series for a maximum of 6 months. To invoice in series, you can use the button **'Show only complete semesters of billing'** to invoice only full semesters. These are the semesters for which the 6 months have elapsed or the semesters that are stopped because the patient is deactivated or deceased.

You can change the amount of the package via the pencil button and you can cancel an invoice, if the patient was not observant, with the trash button.

Increase* of the package for patients with heart failure.

The modulations provided for on the basis of the criteria mentioned in 1° and 2° of Article L. 162-54 and in II of Article R. 162-95 of the Social Security Code, applicable to the rates of the Level 2 operator package, give rise to the following premium packages:

- Premium package level 1: €56;
- Premium package level 2: 70 €.

This premium package can be invoiced at the end of each one-month period during which monitoring is carried out. Billing for the same patient for the increased packages provided for in this Article may not exceed a continuous period of: 6 months, renewable once for the telemedicine activity on heart failure listed in generic form, on the list provided for in Article L. 162-52 of the Social Security Code;

In accordance with Article R. 162-100 of the Social Security Code, limiting the prescription of a remote monitoring activity to one year, the billing of the operator package and any adjustments thereto, in the absence of a new prescription, may not exceed a period of one year without interruption.

Eligibility criteria for the premium package level 1 - €56/month

- Patients aged 80 years or over with at least one comorbidity as provided for in 2° of III of this article;
- Patients under 80 years of age with at least two comorbidities listed in 2° of III of this article;
- Patients hospitalized in the last 30 days for cardiac decompensation.

Eligibility criteria for the premium package level 2 - €70/month

- Patients aged 80 years or older with at least two comorbidities as provided for in 2° of III of this article;
- Patients under 80 years of age with at least three comorbidities listed in 2° of III of this article;
- Patients with specific pathologies or pathways:
- Victim of cardiogenic shock in the last 6 months;
- Placed on the waiting list for a heart transplant;
- Cardiorenal syndrome and waiting for dialysis;
- Heart failure that is intractable and life threatening for less than one year;
- Restrictive cardiomyopathy, out of hospital with heart failure corresponding to at least stage II NYHA;
- Cardiac amyloidosis, out of hospitalization with heart failure corresponding to at least stage II NYHA.

List of comorbidities:

- Cancer under systemic treatment or radiotherapy treatment
- Chronic Kidney Failure ((only if stage IV or V - eGFR < 45 mL/min/1.73 m²)
- Severe anemia or severe iron deficiency
- Severe malnutrition associated with sideroblastic anemia or severe protein-energy malnutrition

'BILLED' status

Once the deed is on the side, you can indicate it on the platform by checking the '**BILLED**' box. You can find all your invoices in the '**Invoiced**' tab.

The invoices are sorted by patient. To view all the invoices for a patient, click on the small arrow that allows you to scroll down all the lines.

Billing							
Search							
To bill Billed							
BILLED	LAST NAME, FIRST NAME ▲	NUMBER OF INVOICES	DATE OF BIRTH	NATIONAL ID	PAYER CODE	PRESCRIBING DOCTOR	
> ☐	TESTROUTIER TESTCHARLINE	23	10/10/1987	1xxxx 33xxxx	Not specified	Prescribing doctor	
> ☐	Z - TEST CLARA	14	08/11/1995	5xxxx 33xxxx	Not specified	Prescribing doctor	
> ☐	Z - TEST FAKE SSNUM	3	08/11/1958	5xxxx 33xxxx	Not specified	Prescribing doctor	
> ☐	Z - TEST JULIO	1	01/01/1958	5xxxx 33xxxx	Not specified	Prescribing doctor	
▼ ☒	Z TEST MAEVA	1	06/12/1956	5xxxx 33xxxx	Not specified	Prescribing doctor	
BILLED	ORDER DATE	INVOICE DATE	INVOICE ID	DATE OF TRANMMISION	OBSERVATION	CODE - AMOUNT	RETURN TO INVOICING TAB
☒	12/13/2023	01/11/2024	PAPER	02/01/2024	40%	TVB - €28	
> ☐	ZTEST CLARA	2	08/11/1956	5xxxx 33xxxx	Not specified	Prescribing doctor	

After verifying that you have been paid by the Health Insurance for the remote monitoring service, you can indicate this next to the corresponding invoice.

In case of an error, you can return the invoice to the 'To be billed' tab using this button.

An Excel export of all invoices is available with the button located at the top right of the page.

After specifying the export period, an excel file is generated. This Excel includes all the invoices on the screen over a desired period: one tab with the invoices to be invoiced and one tab with the invoices made.

14. Contact us

If you have any questions, please contact NEWCARD, a simplified joint-stock company, at **03.20.14.71.14**

The head office is located at 236 Avenue Clément Ader 59118 Wambrechies.

NEW-NTU-PLTF-PS-00

*This user manual NEW-NTU-PLTF-PS-ENG-00, created on **June 12, 2026**, may be subject to change depending on updates to the platform. If this is the case, we will be sure to inform you of any changes.*

The patient must be informed that remote monitoring does not constitute emergency care and that recorded data or alerts are not read and interpreted instantly, but only during the days and hours set by the remote monitoring operator. The patient is therefore informed that in case of emergency, he must contact an emergency call number. - <https://www.legifrance.gouv.fr/jorf/id/JORFTEXT000047720200>

The 1 Minute for my heart and 1 Minute for my kidneys solutions are class I medical devices intended for remote monitoring of heart failure and renal insufficiency, in particular for patients with stage IV or V CKD, as well as to kidney transplant patients who have been transplanted more than three months ago. Manufactured by Newcard, these regulated health products carry, under this regulation, the CE marking. Newcard also provides patients with equipment in addition to the remote monitoring software, in order to facilitate data collection and analysis. These devices are not emergency systems and must be used according to the instructions in the user manual. They are 100% covered by the Health Insurance.

Head office: Wambrechies - 821 766 961 - Registered with the Trade and Companies Register of Lille Métropole B.

Creation date: June 12, 2026 - Image: Newcard - Do not litter in public areas - NEW-NTU-PLTF-PS-00.